



AUG 04, 2023

OPEN ACCESS



DOI:
dx.doi.org/10.17504/protocols.io.261ge3qpw147/v1

Protocol Citation: Narayana Yadavalli, Shawn M. Ferguson 2023. iPSC differentiation into Microglia. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.261ge3qpw147/v1>

MANUSCRIPT CITATION:

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
 We use this protocol and it's working

🌐 iPSC differentiation into Microglia

Narayana Yadavalli^{1,2,3,4,5}, Shawn M. Ferguson^{1,2,3,4,6,5}

¹Department of Cell Biology, Yale University School of Medicine, New Haven, Connecticut 06510, USA;

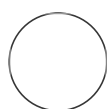
²Neuroscience, Yale University School of Medicine, New Haven, Connecticut 06510, USA;

³Program in Cellular Neuroscience, Neurodegeneration and Repair;

⁴Wu Tsai Institute Yale University School of Medicine, New Haven, Connecticut 06510, USA;

⁵Aligning Science Across Parkinson's (ASAP) Collaborative Research Network, Chevy Chase, MD, 20815, USA;

⁶Kavli Institute for Neuroscience, Yale University School of Medicine, New Haven, Connecticut 06510, USA



berrak.ugur

ABSTRACT

This protocol describes iPSC differentiation into microglia.

ATTACHMENTS

[724-1814.docx](#)

GUIDELINES

This protocol is adapted from the below article.

McQuade A, Coburn M, Tu CH, Hasselmann J, Davtyan H, Blurton-Jones M. Development and validation of a simplified method to generate human microglia from pluripotent stem cells. *Mol Neurodegener.* 2018;13(1):67.

The differentiation protocol involves two steps.

1. iPSC differentiation into CD34 positive hematopoietic progenitors
2. Hematopoietic progenitors' differentiation into mature Microglia

MATERIALS

Reagents required



STEMdiff™ Hematopoietic Kit 1 Kit STEMCELL Technologies Inc. Catalog #5310

Created: May 24, 2023

Last Modified: Aug 04, 2023

PROTOCOL integer ID: 82378

Keywords: iPSC differentiation, microglia



Matrigel Growth Factor Reduced (GFR) Basement Membrane Matrix, LDEV-free Corning Catalog #356230



Gibco™ DMEM/F-12 HEPES Fisher Scientific Catalog #11-330-032



MEAA (MEM Non-Essential Amino Acids) Gibco - Thermo Fischer Catalog #11140050

Glutamax (Gibco,)



Gibco™ N-2 Supplement (100X) Thermo Fisher Scientific Catalog #17502048



Gibco™ B-27™ Supplement (50X) serum free Fisher Scientific Catalog #17-504-044



Recombinant Human M-CSF peprotech Catalog #300-25



Recombinant Human TGF-β1 (HEK293 derived) peprotech Catalog #100-21



Recombinant Human IL-34 peprotech Catalog #200-34



Recombinant Human Fractalkine (CX3CL1) peprotech Catalog #300-31

CD200 (Novo protein #C311)



1-Thioglycerol Merck MilliporeSigma (Sigma-Aldrich) Catalog #M6145



Insulin-Transferrin-Selenium (ITS -G) (100X) Thermo Fisher Catalog #41400045



Insulin Merck MilliporeSigma (Sigma-Aldrich) Catalog #I2643-25MG

Base media:

- DMEM/ F12

- 2X insulin-transferrin-selenite
- 2X B27, 0.5X N2
- 1X glutamax
- 1X non-essential amino acids
- 400µM monothioglycerol
- 5 µg/mL insulin

3 cytokine media:

Base media	
mCSF	25 ng/ml
IL34	100 ng/ml
TGFbeta	50 ng/ml

5 Cytokine media:

Base media	
mCSF	25 ng/ml
IL34	100 ng/ml
TGFbeta	50 ng/ml
CX3CL1	100 ng/ml
CD200	100 ng/ml

Steps for iPSC differentiation into CD34 positive hematopo...

1

Note

Note: This differentiation protocol requires healthy iPSCs approximately 50-60% confluency with nice tight colonies.

Day-2: Coat 3 wells of a 6 well pates with Matrigel. (1 hour coating also works)

2




Day-1:

2.1

Bring iPSC maintenance to plate cell culture hood, remove the media and rinse once with PBS.



2.2

Then add  1 mL of  0.5 micromolar (μM) EDTA to the well and leave in the incubator for about  00:05:00. By this time, you will see colonies lifting from the plates. If not leave plate in the incubator for few more minutes.


5m



Note

Note: Do not tap or pipet to lift off the cells from the plate. This protocol requires a clumps off iPSC colonies for good yield of hematopoietic progenitors.

2.3

Once the 50 to 60% of colonies come off the plate, bring plate into the hood and neutralize the reaction by adding  2 mL of E8+Ri media.




2.4

Now gently swirl the plate to mix EDTA solution containing cells and E8+Ri media.

Note

Never pipet, this mechanical force will disrupt the clumps.

2.5

After swirling the plate pipet  1 mL of the cells with 5 ml pipet into 15 ml falcon tube.



2.6

Spin down the cells at  1 rcf, 00:03:00.

3m




Note



This slow spin is required for avoiding single cells smaller clumps.

2.7


After this spin, take aspirate the supernatant and gently tap the cell pellet.

2.8 Then add  3 mL of E8+Ri media and tap gently once again to mix the clumps.



2.9 For counting the clump number take  5 μ L of this mix into a 96 well plate, add  100 μ L of E8 media. Gently tap the plate and count the colonies under the microscope. Repeat this in 3 wells of 96 well plate and average the clump number.



2.10 Now bring the Matrigel coated plate, aspirate the Matrigel and add  1.5 mL of E8+Ri media. (3 well were coated on day0) Now plate 10 clumps in 1st well, 20 clumps in 2nd well and 30 clumps in 3rd well of the 6 well plate.* Shake the plate up and down and to sides, return the plate to the incubator.



Note

*Since counting of clumps is arbitrary, it is never perfect. This exercise is essential till you get very good experience in plating the clumps.

3 Day 0

3.1 Count the colonies in in each well, ideal colony number is between 10-30. It doesn't have to be perfect number. Wells with as low as 5 colonies and max up to 30 can also be used. Anything above 40 should be avoided.

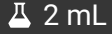

Note

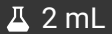

Note:

1. Each clump should have 20-40 cells. If you have smaller colonies change media to E8 only and let the colonies grow for 1-2 days.
2. Above 30 clumps number differentiation may not work because mesodermal cells require space to migrate and differentiate into hematopoietic progenitors.

3.2 Choose on well with desired number of colonies and you can discontinue maintaining

remaining wells.

3.3 After achieving desired colony number prepare Media A ( 2 mL base media +  10 μ L supplement A).


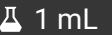

3.4 Aspirate E8+Ri media and add  2 mL Media A and leave plate in incubator for  48:00:00 .



2d


4 Day 2




Add  1 mL A ( 1 mL base media +  5 μ L supplement A).

5 Day 3




5.1 Prepare Media B ( 2 mL base media +  10 μ L supplement B).

5.2 Aspirate media A and add  2 mL Media B.







6 Day 5, 7, 9, 10
supplement the cell with  1 mL of Media B.

7 Day 12: Collection


- 7.1 By day 12 you will see lot of floating hematopoietic progenitor cells.
- 7.2 Collect hematopoietic progenitor cells by gently swirling the plate with a 5 ml pipet into 15 ml falcon tube.
- 7.3 Spin down the cells  3 rcf, 00:03:00 3m
- 
- 7.4 Remove the supernatant and resuspend the cell pellet in Macrophage differentiation media. (RPMI+ 20% FBS+  100 ng/mL M-CSF).
- 7.5 Count the cell by using hemocytometer and plate 100,000 in one well of 6 well plate.

Steps for Microglia differentiation from hematopoietic prog.. 6m

- 8 **Day 12:**
Plate 100,000 hematopoietic progenitor on Matrigel coated 6 well plate in 3 cytokine media.
- 9 On days 14,16,18,20 and 22 supplement with  1 mL of 3 cytokine media.
- 10 **Day 24**


10.1 Collect  6 mL of cells + media into 15 ml falcon by leaving  1 mL conditioned media in the plate. Spin down at  3 rcf, 00:03:00 and remove the supernatant.

3m

10.2 Resuspend the pellet in  2 mL of fresh 3 cytokine media and plate back into the same well containing conditioned media.


11 Day 26,28,30,32,34,36 supplement with  1 mL of 3 cytokine media.


12 Day 37

12.1 Collect  6 mL of cells + media into 15 ml falcon by leaving 1ml conditioned media in the plate.

3m

Spin down at  3 rcf, 00:03:00 and remove the supernatant.

12.2 Resuspend the pellet in  2 mL of fresh 5 cytokine media and plate back into the same well containing conditioned media.

13 Day 39: supplement cells with  1 mL of 5 cytokine media.

14 Day 41: Collect cells for experiment.

