



# GFP-VSV Infection V.1

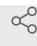
Christopher Rousso<sup>1</sup>, Alison Macdonald<sup>1</sup>

<sup>1</sup>UOttawa


Version 1 ▼

Sep 13, 2022

1 *Works for me*

 Share

This protocol is published without a DOI.

 Kelsey Grimes

## ABSTRACT

GFP-VSV Infection

## PROTOCOL CITATION

Christopher Rousso, Alison Macdonald 2022. GFP-VSV Infection. **protocols.io**  
<https://protocols.io/view/gfp-vsv-infection-cgcitsue>



## LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

## CREATED

Sep 08, 2022

## LAST MODIFIED

Sep 13, 2022

## PROTOCOL INTEGER ID

69738

## MATERIALS TEXT

- Aliquot of GFP-VSV virus (from -80°C)
- Serological pipets
- 1.5 mL Eppendorf tubes
- Serum free DMEM media (high glucose) [depending on cell type]
- DMEM media (high glucose + 10% FBS) [depending on cell type]
- Ice bucket

1 If infecting transfected cells, infect them 24 hours after transfection.

- 2 Remove cells from 37 °C incubator.
- 3 Replace media with fresh media (with FBS) prior to infection with GFP-VSV.
- 4 Thaw VSV-GFP stock aliquot on ice
- 5 Dilute VSV-GFP in serum-free culture media. Perform the necessary calculations to determine the dilution that is required to obtain the desired MOI (multiplicity of infection).  
Note: MOI is the number of viruses/number of cells. For example, a MOI of 1 = 500,000 PFU (plaque forming units) of virus per 500,000 cells.  
Note: We typically start with a dilution range of MOI = 0.01, 0.1, 1 for optimization
- 6 Infect cells by dropping the diluted virus directly into the wells making sure to disperse the virus around the well.
- 7 Gently swirl the plate to help disperse the virus.
- 8 Leave the plate for approximately 24 hours in a 37°C CO<sub>2</sub> incubator. Then proceed with endpoint analysis.

#### Appendix - Sample Calculation

- 9 *Sample calculation:*  
 (MOI = 1; Virus titre =  $8.4 \times 10^8$  virus/mL ; 500,000 cells/well;)  
*Step 1: calculate amount of PFU per well*

Step 2: calculate how much to dilute the viral stock

$$C_1V_1 = C_2V_2$$

Where  $V_1$  = volume/well you have

$C_1$  = virus/mL calculated in step 1

$C_2$  = virus titre (stock concentration)

$V_2$  = x (unknown volume to be calculated)

You will add x mL of your virus stock per well. This volume will be low, and thus must be diluted in media for easier pipetting.