



OCT 04, 2023

OPEN  ACCESS

**Protocol Citation:** Brian Lovett, Kristen Pierce 2023. Primer stock preparation. **protocols.io** <https://protocols.io/view/primer-stock-preparation-c2tyyepw>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
We use this protocol and it's working

**Created:** Oct 02, 2023

**Last Modified:** Oct 04, 2023


**PROTOCOL integer ID:**  
88664

## Primer stock preparation

Brian Lovett<sup>1</sup>, Kristen Pierce<sup>2</sup>

<sup>1</sup>USDA-ARS; <sup>2</sup>West Virginia University

Lovett Lab - USDA ARS EPPRU

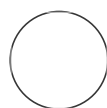

Brian Lovett  
USDA-ARS

### ABSTRACT

General protocol for preparation of lyophilized primer stocks from IDT.

### GUIDELINES

**Recommended:** When you receive your lyophilized primers, check the name of the primer, primer sequence and the concentration listed on the tube. IDT provides a specification sheet you can download/archive for more information.

### MATERIALS

Lyophilized IDT primers, molecular grade water, pipettes and pipette tips.

**Keywords:** lyophilized  
primer reconstitution, primer  
dilution

## Preparation of stock solution

- 1 Reconstitute lyophilized primers in required volume of molecular grade water (i.e., RNase/DNase free) to **100 micromolar ( $\mu\text{M}$ )**.
- 1.1 Practically, the microliters of water required to create a 100  $\mu\text{M}$  solution is 10x the nmol concentration of lyophilized primer (i.e., if **23.5 Nanomoles** is noted on side of tube, then add **235  $\mu\text{L}$**  of water to stock tube).
- 2 After water has been added, vortex for at least **00:00:10** to ensure the lyophilized primers are fully dissolved. 10s
- 3 Once reconstituted, the stock solution can be stored in a labeled box at **-20  $^{\circ}\text{C}$** , and thawed completely before use.

## Preparation of working solution

- 4 Allow stock solution to thaw completely at **Room temperature**.
- 4.1 While waiting for your stock to thaw, prepare clean microcentrifuge tubes by labeling them with the primer name and date of preparation. Label the lid so it is easy to identify tubes in freezer box.
- 5 Vortex the stock solution for at least **00:00:10** to ensure homogenization before adding it to the working solution. 10s

6 To prepare a **10 Molarity (M)** working solution from the **100 micromolar ( $\mu\text{M}$ )**, you will dilute 1:10 (i.e., **1  $\mu\text{L}$**  of stock for every **9  $\mu\text{L}$**  of water).

6.1 For a **100  $\mu\text{L}$**  working solution, pipette **90  $\mu\text{L}$**  of molecular grade water into the newly labeled microcentrifuge tube and then add **10  $\mu\text{L}$**  of the stock solution.

7 Once the concentrated stock is added, vortex the new working solution for at least **00:00:10** and store at **-20  $^{\circ}\text{C}$**  to use as needed.

10s

## Use of primers

8 For a typical **25  $\mu\text{L}$**  PCR reaction, **1  $\mu\text{L}$**  of each primer is used (i.e., **0.4 micromolar ( $\mu\text{M}$ )** final concentration).

8.1 To maintain optimal use of your primers, limit freeze/thaw cycles, as this can lead to degradation over time.