




Jun 30, 2021

Murine CD8 T cell transduction

Kristin Anderson¹¹University of Washington, Fred Hutchinson Cancer Research Center

1 Works for me

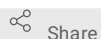

Share

dx.doi.org/10.17504/protocols.io.smrec56

Kristin Anderson

ABSTRACT

CD8 T cells can be transduced to express genes transferred by retroviral transduction. This protocol outlines the steps for generating retrovirus from PLAT-E packaging cells to transduce primary murine CD8 T cells isolated from donor mouse spleens.

In brief: The ectopic retroviral packaging cell line, Platinum (Plat)-E cells (Cell Biolabs, Inc.), were plated (2.2×10^6 cells) in PLAT-E media on 10 mm culture plates. After 24 hours, Plat-E cells were transfected with DNA encoding the desired TCR and/or construct using the Effectene transfection kit (Qiagen). At 48 hours, media was replaced with T cell media. On days 2 and 3 post-transfection (at 72 and 96 hours) virus-containing supernatant was collected. Splenic T cells isolated from P14 Thy1.1+ female mice were stimulated with α -CD3e (clone: 145-2C11) and α -CD28 (clone: 37.51 (both from BD Pharmingen) and IL-2 (50 IU/ml Aldesleukin, UW Pharmacy), and transduced with retroviral supernatant by spinfection in polybrene (5mg/mL, 90 minutes at 1000xg) at 24 and 48 hours after T cell activation.

DOI

dx.doi.org/10.17504/protocols.io.smrec56

PROTOCOL CITATION

Kristin Anderson 2021. Murine CD8 T cell transduction. **protocols.io**<https://dx.doi.org/10.17504/protocols.io.smrec56>

LICENSE

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Aug 14, 2018

LAST MODIFIED

Jun 30, 2021

PROTOCOL INTEGER ID

14737

MATERIALS TEXT

MATERIALS

[DPBS \(no Ca, no](#)

[Mg\) ThermoFisher Catalog #14190144](#)

[Purified Hamster Anti-Mouse CD3e \(Clone 145-2C11\) BD](#)

[Biosciences Catalog #553058](#)

[Purified Hamster Anti-Mouse CD28 \(Clone 37.51\) BD](#)

[Biosciences Catalog #553295](#)

[IL-2 \(Aldesleukin\) Contributed by users](#)

[Effectene Transfection](#)

[Reagent Qiagen Catalog #301427](#)

[PLAT-E \(Platinum-E\) ectropic retroviral packaging cell line Cell Biolabs,](#)

[Inc Catalog #RV-101](#)

[0.05% Trypsin-EDTA \(1x\) Gibco - Thermo](#)

[Fisher Catalog #25300-054](#)

[ACK Lysing Buffer Gibco - Thermo](#)

[Fisher Catalog #A10492-01](#)

[Falcon 5mL polypropylene round-bottom](#)

[tube Corning Catalog #352063](#)

[Falcon 50 mL polypropylene conical](#)

[tube Corning Catalog #352070](#)

[15 mL polypropylene centrifuge](#)

[tube Corning Catalog #430791](#)

[Tissue Culture Dish Techno Plastic Products](#)

[\(tpp\) Catalog #93100](#)

[T25 cell culture flask, canted neck, 0.2um vent](#)

[vap Corning Catalog #430639](#)

[Paraformaldehyde, 16% w/v aq. soln, methanol free liquid Alfa](#)

[Aesar Catalog #43368](#)

[Falcon cell strainer, 70um,](#)

[nylon Corning Catalog #352350](#)

[BD 10 mL syringe, Luer-Lok Tip BD](#)

[Biosciences Catalog #302995](#)

[Nalgene syringe filter, 0.45um \(surfactant free cellulose acetate membrane\) Thermo](#)

[Scientific Catalog #723-9945](#)

SAFETY WARNINGS

Retroviral supernatant requires BSL-2 precautions. Work in a biosafety cabinet and wear proper PPE.

Day 0: Prepare PLAT-E cells for transfection

1. Gently remove PLAT-E media and rinse PLAT-E cells with 1x DPBS.
2. Gently remove PBS and add 0.05% Trypsin. Trypsinize at 37C for 1-2 minutes or until cells begin to lift off the plate with a gentle tap.

3. Dilute and inactivate trypsin with 10ml PLAT-E media. Transfer cells to a conical tube, spin down @ 369 g for 4-5 minutes at 4C and remove supernatant.
4. Resuspend PLAT-E cells in 10mL PLAT-E media. Count cells. Make new 10 cm plates with 2.2e6 cells per plate in 10ml PLAT-E media.
Cells:
Quadrants:
Dilution:
Calculation of concentration:
Calculation of cells to add:
Volume of cells per plate:
Media per plate:
5. Return plates to 37C incubator.
DNA to transduce:
Number of plates prepared:

Day 1: Transfect PLAT-E cells

- 2
 1. Thaw tubes of DNA. And prepare 5mL polypropylene tubes (1/plate).
 2. Add 300 uL EC Buffer.
 3. Add 2ug DNA to EC Buffer.
Record DNA amounts to add for each construct:
 4. Add 16uL of Enhancer and vortex for 1 second.
 5. Incubate 5 min at Room Temp.
 6. Add 60uL of Effectene. Vortex for 10 seconds.
 7. Incubate 10 min at Room Temp.
 8. Aspirate media from PLAT-E plates. Add 7mL of T cell media to each plate.
 9. Add 3 mL of T cell media to each tube of Effectene mixture. Pipette up and down to mix, then add slowly (dropwise) to each PLAT-E plate.
 10. Swirl gently to mix and return to 37C incubator overnight.

Day 2: Change PLAT-E media and prepare T cells for transduction

- 3
 1. Aspirate media from PLAT-E plates and replace with 10.5 mL T cell media.
 2. Place the plates in the 32C incubator overnight.
 3. Harvest T cells by smashing donor spleens through a 40um filter. Rinse and transfer to a 15 or 50 mL tube.
Number of spleens:
Donor mice information:
 4. Spin down T cells down @369 g for 4-5 minutes at 4C.
 5. Aspirate media and lyse spleens with 2ml/spleen ACK Lysis buffer for 2 minutes.
 6. Dilute out the ACK Buffer with equal or greater volume of T cell media. Spin down T cells down @ 369 g for 4-5 minutes at 4C.
 7. Resuspend pellet in T cell media.
Volume of T cells:
 8. Add 1ug/ml of anti-CD3 and anti-CD28 (e.g. 20ul of Ab per 10 mL of splenocytes).
Anti-CD3 added:
Anti-CD28 added:
 9. Add 50 IU/mL of IL-2 (e.g. 10ul of Ab per 10 mL splenocytes).
IL-2 added:
 10. Transfer 10mL of splenocytes per T25 flask and place in 37C incubator overnight.

Day 3: Transduction 1

- 4
 1. Transfer T cells into 50mL tubes (need 1 tube for each independent viral transduction. Number of tubes:
 2. Harvest virus from PLAT-E plates and filter through 0.45um Nalgene Syringe Filter.
 3. Add 10.5 mL T cell media per PLAT-E plate and return to 32C incubator overnight.
 4. Resuspend T cells in filtered viral supernatant. Calculate volume of each tube.
Volumes:
 5. Add IL-2 at 50 IU/mL and Polybrene at 5ug/mL.
Volume of IL-2 added to each tube:

Volume of Polybrene added to each tube:

6. Plate T cell/virus mixture in 12 well plates (Add approximately 1.5-2 mL per well). Balance the plates by weight. Spin @ 1000 \times g for 90 minutes at 30C.
7. Place plates in 37C incubator overnight.

Day 4: Transduction 2

- 5
 1. Collect T cells from wells, combine similar wells into 50 mL tube, and spin down @1350 RPM for 4-5 min at 4C.
 2. Harvest virus from PLAT-E plates and filter through 0.45um Nalgene Syringe Filter.
 3. Resuspend T cells in filtered viral supernatant. Calculate volume of each tube.
Volumes:
 4. Add IL-2 at 50 IU/mL and Polybrene at 5ug/mL.
Volume of IL-2 added to each tube:
Volume of Polybrene added to each tube:
 5. Plate T cell/virus mixture in 12 well plates (~1.5-2.5 mL per well). Balance the plates by weight. Spin @ 1000 \times g for 90 minutes at 30C.
 6. Place plates in 37C incubator overnight.

Day 6: Add media (as needed) and IL-2

- 6
 1. If the wells are starting to change from orange to yellow in color, add 1-2 mL of T cell media per well of cells.
 2. Add 50 IU/ml IL-2 per well.

Day 7 or 8: Screen T cells for transduction efficiency

- 7
 1. Remove a small sample of cells for screening.
 2. Continue to add media (as needed) and IL-2 (50 IU/ml) every 48 hours.