




FEB 10, 2023

# 🌐 Preparation and transformation of chemically super-competent Escherichia coli

Andreas Sagen<sup>1</sup>

<sup>1</sup>University of Oslo

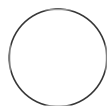

Andreas Sagen

University of Oslo, The National Institute of Occupational H...

## ABSTRACT

Based on a method produced by Inoue, et al.

Inoue, Nojima, H., & Okayama, H. (1990). High efficiency transformation of Escherichia coli with plasmids. *Gene*, 96(1), 23–28. [https://doi.org/10.1016/0378-1119\(90\)90336-P](https://doi.org/10.1016/0378-1119(90)90336-P)

## MATERIALS

LAF

Scale

Ice

Centrifuge

Incubator

## OPEN ACCESS

**Protocol Citation:** Andreas Sagen 2023. Preparation and transformation of chemically super-competent Escherichia coli. **protocols.io** <https://protocols.io/view/preparation-and-transformation-of-chemically-super-cn5uv6w>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** In development  
We are still developing and optimizing this protocol

**Created:** Feb 08, 2023







**Last Modified:** Feb 10, 2023

**PROTOCOL integer ID:**  
76692

**Keywords:** Transformation, Inoue method, E. coli

## Preparation of Inoue transformation buffer

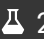1 In a sterile flask, add  200 mL distilled water

- 2 Measure  2.5 mL (  1 Molarity (M) ) PIPES,  2.177 g Manganese(II) chloride tetrahydrate,  3.75 mL (  1 Molarity (M) ) Calcium chloride and  4.660 g Potassium chloride.

Materials:

 PIPES **Sigma-aldrich Catalog #P1851**

 Manganese(II) chloride tetrahydrate **Sigma-aldrich Catalog #M3634**

 Calcium chloride **Sigma-aldrich Catalog #C3881**

 Potassium chloride **Sigma-aldrich Catalog #P4504**

- 3 Add measured reagents and mix for  00:05:00


5m

- 4 Adjust pH to  6.7 with Potassium hydroxide solution



Materials:

 Potassium hydroxide solution **Supelco Catalog #P4494**

- 5 Fill flask with distilled water to  250 mL

- 6 Filter sterilize solution with a filter (0.2  $\mu$ m) and store refrigerated (  4 °C )



## Preparation of chemically super-competent cells

- 7 Prepare a culture of *E. coli* on an LB agar plate. Pick a single colony, and inoculate in  500 mL S. O. C. broth in a  1000 mL flask.

Incubate at  18 °C with shaking at 100 rpm overnight, until OD<sub>600</sub> reaches 0.6

- 8 Aliquot entire culture volume into  50 mL canonical tubes

- 9 Place tubes on ice for  00:10:00 10m
- 10 Centrifuge tubes with  5000 rcf, 4°C, 00:10:00 10m
- 11 Discard supernatant and resuspend with  16 mL CRM transformation buffer
- 12 Incubate cells on ice for  00:10:00 10m
- 13 Centrifuge tubes with  5000 rcf, 4°C, 00:10:00 10m
- 14 Discard supernatant and resuspend with  8 mL Inoue transformation buffer. Pool into two tubes
- 15 Centrifuge tubes with  5000 rcf, 4°C, 00:10:00 . Meanwhile, prepare  50 mL DMSO-Inoue transformation buffer by diluting  3.5 mL DMSO in  46.5 mL Inoue transformation buffer 10m
- Materials:
-  **DMSO MP Biomedicals Catalog #196055**
- 16 Discard supernatant and resuspend with  10 mL DMSO-Inoue per tube
- 17 Incubate cells on ice for  00:30:00 30m

- 18 Aliquot  100  $\mu\text{L}$  cell suspension into sterile  500  $\mu\text{L}$  screw cap reaction tubes (Sarstedt #72.704.200).

**Note**

While creating aliquots, keep original tubes, and aliquots on ice, until snap-freeze take place

- 19 Snap freeze tubes in liquid nitrogen using a floating foam tube rack (Southern labware #HS2166)

- 20 Transfer aliquots storage box, and place in an ultra-low temperature freezer or vapor-phase nitrogen tank

**Note**

Store tubes in 50 mL canonical tubes, or similar containers

## Transformation

1h 32m 30s

- 21 Quickly thaw a single reaction tube with  100  $\mu\text{L}$  hyper-competent cells in Inoue-DMSO









- 22 Mix 1-5  $\mu\text{L}$  plasmid (ligation product)

**Note**

Do not exceed 5% of the volume competent cells

**Note**

Use up-to  25 ng per  50  $\mu\text{L}$  of competent cells

- 23 Incubate cells on ice for  00:30:00 30m
- 24 Heat-shock cells at  42 °C for  00:00:30, followed by  00:02:00 at  4 °C immediatly 2m 30s
- 25 Add  500 µL prewarmed S. O. C. medium and incubate for  37 °C at 200 rpm for  01:00:00 1h
- 26 Add desired amount of suspension on LB plates with ampicillin (100 µg/mL) and incubate overnight