



Sep 06, 2022

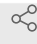
🌐 Preparing mRNA for nucleofection of hPSCs

 In 1 collection

Hanqin Li¹, Oriol Busquets², Steven Poser², Dirk Hockemeyer¹,
Frank Soldner²

¹University of California, Berkeley; ²Albert Einstein College of Medicine

1 *Works for me*

 Share

dx.doi.org/10.17504/protocols.io.b4peqvje



Devin E Snyder

ABSTRACT

This protocol describes the standard procedure for preparing mRNA to be delivered into human pluripotent stem cells (hPSCs) using nucleofection.

General notes

1. Throughout this protocol, the term hPSC is used to collectively refer to both hiPSCs and hESCs. All described procedures have been tested and work equally well for hiPSCs and hESCs.

DOI

dx.doi.org/10.17504/protocols.io.b4peqvje

PROTOCOL CITATION

Hanqin Li, Oriol Busquets, Steven Poser, Dirk Hockemeyer, Frank Soldner 2022.
Preparing mRNA for nucleofection of hPSCs. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.b4peqvje>



FUNDERS ACKNOWLEDGEMENT

Aligning Science Across Parkinson's
Grant ID: ASAP-000486

COLLECTIONS ⓘ



Nucleofection (Amaxa) and electroporation (Biorad) of hPSCs

KEYWORDS

ASAPCRN

LICENSE

_____ This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Feb 03, 2022

LAST MODIFIED

Sep 06, 2022

PROTOCOL INTEGER ID

57798

PARENT PROTOCOLS

Part of collection

[Nucleofection \(Amaxa\) and electroporation \(Biorad\) of hPSCs](#)

MATERIALS TEXT

Item	Vendor	Catalog #
Synthetic pegRNAs	IDT or Synthego	
Synthetic sgRNAs	Synthego	

Note: This protocol makes reference to other protocols. Please check for any materials found in those protocols, which might not be listed here

- 1 Thaw mRNA and synthetic pegRNA/ngRNA [On ice](#)
- 2 In each nucleofection, use 4 µg total of mRNA
- 3 For prime editing **PE2 strategy**, use:
4 µg, IVT PE2 mRNA
1.5 µl, 100 µM Synthetic pegRNA

- 4 For prime editing **PE3 strategy**, use:
 - 4 µg, IVT PE2 mRNA
 - 1 µl, 100 µM Synthetic pegRNA
 - 0.5 µl, 100 µM, Synthetic ngRNA
- 5 Pipet the proper amount of each component into a microcentrifuge tube.