



Version 2 ▾

Sep 11, 2020

covid 19 indirect detection thru rise of 100nm filter fluid resistance V.2

Ladislav Nevery

In Development

dx.doi.org/10.17504/protocols.io.bmajk2cn

Coronavirus Method Development Community

XPRIZE Rapid Covid Testing

1 more workspace



Ladislav Nevery

ABSTRACT

By passing saliva sample thru 150nm filter to remove particles larger than avg 125nm sized covid 19 virus particles and concentrating resulting fluid to very small 1 mm² area of 110nm filter. Pores in this small area should be clogged up fast raising filter air/fluid resistance. Due to high virus density 5.2 log₁₀ ml saliva of infected person should clogg filter much faster than healthy person. Test is setup in such way that we let gravity pass 2l of fluid thru clogged filter and record rate of drops from its bottom as audio on mobile phone placed bellow glass cup. rate of drops thru clogged filter should be measurably slower than clean filter.

THIS PROTOCOL ACCOMPANIES THE FOLLOWING PUBLICATION

dx.doi.org/10.17504/protocols.io.bkznkx5e

ATTACHMENTS

[DropRate.zip](#)

DOI

dx.doi.org/10.17504/protocols.io.bmajk2cn

PROTOCOL CITATION

Ladislav Nevery 2020. covid 19 indirect detection thru rise of 100nm filter fluid resistance. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.bmajk2cn>
 Version created by [Ladislav Nevery](#)



MANUSCRIPT CITATION please remember to cite the following publication along with this protocol

dx.doi.org/10.17504/protocols.io.bkznkx5e

LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Sep 11, 2020

LAST MODIFIED

Sep 11, 2020

PROTOCOL INTEGER ID

42027

MATERIALS

NAME

CATALOG #

VENDOR

NAME	CATALOG #	VENDOR
pragopor 9		
pragopor 10		
120nm nanoparticles		


STEPS MATERIALS

NAME	CATALOG #	VENDOR
3d print filament 8m		
pragopor 10		pragochema
pragopor 9		
2l soda bottle		
glass cup		
mobile phone		
120nm nanoparticles		

DISCLAIMER:

DISCLAIMER – FOR INFORMATIONAL PURPOSES ONLY; USE AT YOUR OWN RISK

The protocol content here is for informational purposes only and does not constitute legal, medical, clinical, or safety advice, or otherwise; content added to protocols.io is not peer reviewed and may not have undergone a formal approval of any kind. Information presented in this protocol should not substitute for independent professional judgment, advice, diagnosis, or treatment. Any action you take or refrain from taking using or relying upon the information presented here is strictly at your own risk. You agree that neither the Company nor any of the authors, contributors, administrators, or anyone else associated with protocols.io, can be held responsible for your use of the information contained in or linked to this protocol or any of our Sites/Apps and Services.

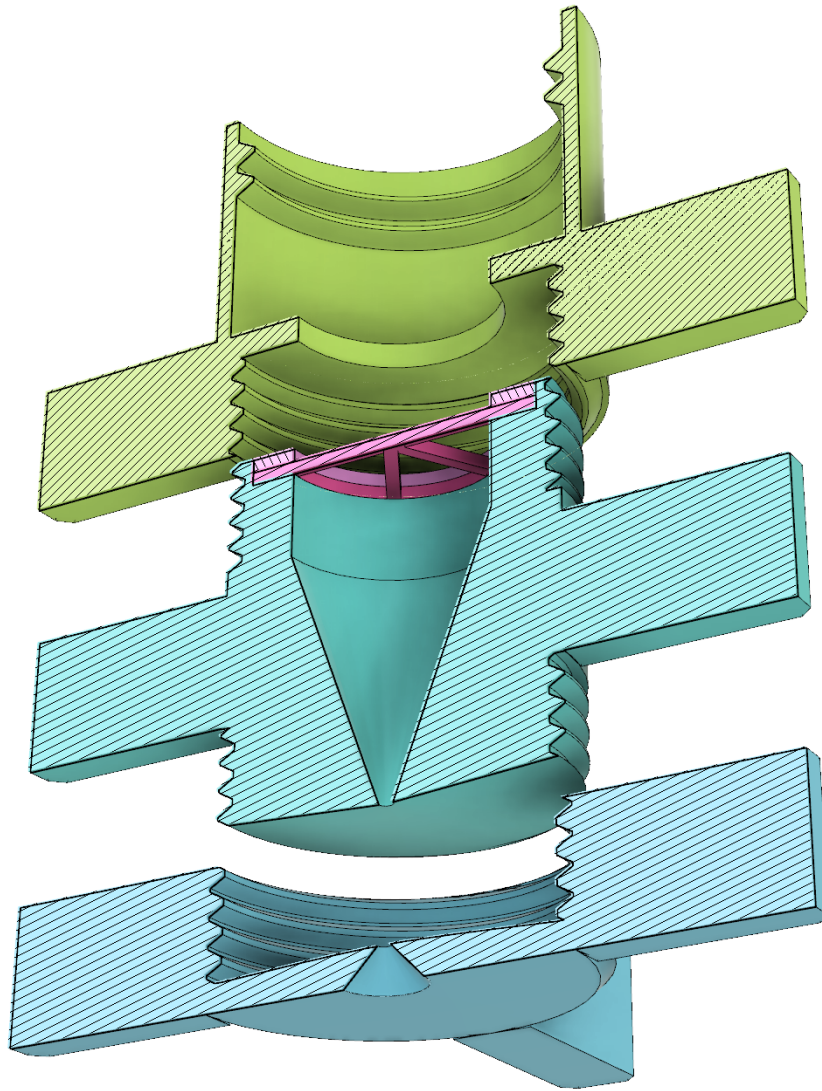

3d print all required parts



3d print filament 8m

2 place 24mm ideally 110nm pc or nylon membrane filter between blue parts. I used what I had



pragopor 10

by pragochema

[View](#)




0.24eur piece

- 3 place 34mm ideally 150nm membrane filter between red parts. I used what I had



- 4 install green saliva holder

- 5 .attach

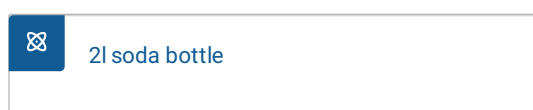

and squeeze until all fluid passes thru both filters


- 6 remove 2l bottle remove red parts including filter. clean

- 7 fill 2l with clean fluid and reinstal it.

- 8 place whole assembly over glass cup



with bottle on top. place mobile phone bellow glass cup and start audio recording to count and record sounds of drops.

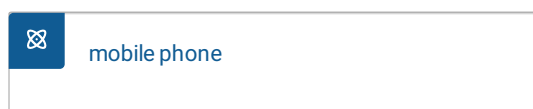

- 9 fill sample holder in step 4 with fluid with



in concentration resembling covid19 as
5.2log10ml and repeat whole process once more ↻

10 compare three recordings.



drop frequency with 120nm nanoparticle clogged filter should be measurably different