



Jun 19, 2022

DNA extraction protocol for historical toe pad samples from birds

Dave Lutgen^{1,2,3}, Reto Burri^{3,4,2}¹Institute of Ecology and Evolution, University of Bern, Switzerland;²Institut of Ecology and Evolution, Friedrich-Schiller-University Jena, Germany;³Swiss Ornithological Institute, Sempach, Switzerland;⁴Institute of Ecology and Evolution, University of Bern, Bern, Switzerland

1

dx.doi.org/10.17504/protocols.io.dm6gpwrplzp/v1 lutgendave

Museums hold collections of specimens from vast taxonomic and geographic ranges that constitute rich resources for research into the origins, evolution, and maintenance of biodiversity on earth. However, specimens are often not prepared for specific research purposes. Especially the isolation of adequate quantities of DNA for genomic research can be challenging. Here, we report a modified protocol of a commercially available kit optimized to reliably extract adequate quantities of DNA from bird toepads and dried skin from sutures of museum specimens for whole-genome resequencing.

DOI

dx.doi.org/10.17504/protocols.io.dm6gpwrplzp/v1

Dave Lutgen, Reto Burri 2022. DNA extraction protocol for historical toe pad samples from birds. **protocols.io**

<https://dx.doi.org/10.17504/protocols.io.dm6gpwrplzp/v1>



German Research Foundation (DFG)

Grant ID: BU3456/3-1

National Research Fund (FNR) Luxembourg

Grant ID: 14575729

Birds, Museomics, Illumina sequencing, historical DNA, Toe pads, whole-genome resequencing

protocol ,

Oct 07, 2020

Jun 18, 2022

42861

Material

- [QiAmp DNA Micro](#)
kit Qiagen Catalog #56304 Step 1
- [Micropestle Sigma](#)
Aldrich Catalog #BAF199230001 Step 2
- [1.5 mL Safe Lock](#)
Tube Eppendorf Catalog #0030120086 In 5 steps

Digestion

- [Proteinase K,](#)
1. **40 µL 2x20µl** **2mL Qiagen Catalog #19131** In 3 steps
[Buffer](#)
- 2. **360 µL** **ATL Qiagen Catalog #19076** Step 6
[RNase](#)
- 3. **4 µL** **A Qiagen Catalog #19101** In 2 steps
[AL lysis](#)
- 4. **400 µL** **buffer Qiagen Catalog #19075** In 4 steps

Extraction

- [Buffer AW1, Wash buffer \(1\),](#)
concentrate Qiagen Catalog #19081 Step 22
- [Buffer AW2, Wash buffer \(2\),](#)
concentrate Qiagen Catalog #19072 Step 23
- [AE](#)
buffer Qiagen Catalog #19077 In 3 steps
- [Ultrapure Distilled, Nuclease Free Water Contributed by users](#) In 2 steps

Quality control

 [Qubit® 3.0 Fluorometer](#) **Thermo Fisher**

Scientific Catalog #Q33216

 [NanoDrop spectrophotometer](#) **Thermo Fisher**

Scientific Catalog #ND-1000

 [TapeStation 4200](#) **Agilent**

Technologies Catalog #G2991BA

 [D1000 Screen Tape](#) **Agilent Technologies**

Preparation

1 

 [QiAmp DNA Micro](#)

This protocol uses the [kit Qiagen Catalog #56304](#) with 2x the amount/volumes of

 [Proteinase K,](#)

2mL Qiagen Catalog #19131

 [Buffer ATL \(tissue lysis](#)

buffer) Qiagen Catalog #19076

 [AL lysis](#)

buffer Qiagen Catalog #19075

 [RNase](#)

A Qiagen Catalog #19101

2

 [Micropestle](#) **Sigma**

Prepare a **Aldrich Catalog #BAF199230001**

per sample.

 [1.5 mL Safe Lock](#)

Prepare a **Tube Eppendorf Catalog #0030120086** to store the micropestle separately for each sample.

 [AL lysis](#)

Preheat **buffer Qiagen Catalog #19075**

at **56 °C**. Adding preheated

 [AL lysis](#)

buffer Qiagen Catalog #19075

prevents the formation of potential

precipitates.

Tissue preparation

- 3 Put toe pad/breast skin tissue (2-3mg) into [1.5 mL Safe Lock](#) [Tube Eppendorf Catalog #0030120086](#) (use scalpel tip to transfer).
- 4 Museum toe pad samples might have been treated with chemicals to better conserve them, therefore, we will rinse the tissue before starting the digestion
 - 4.1 Add [1000 µL](#) of [Ultrapure Distilled, Nuclease Free Water Contributed by users](#) to wash down chemicals on the toe pad. Incubate at [1000 rpm, 56°C, 01:00:00](#). Discard all [Ultrapure Distilled, Nuclease Free Water Contributed by users](#). Make sure no [Ultrapure Distilled, Nuclease Free Water Contributed by users](#) is left.
- 5 Mince the tissue with a pointed scalpel blade. Use a new blade for each sample!

Digestion

- 6 To make sure no tissue is left on the blade, pipet [80 µL](#) of [Buffer](#) [ATL Qiagen Catalog #19076](#) over the scalpel blade, then remove it and ensure all tissue remains in the digestion.
- 7 [Proteinase K](#),
Add [20 µL](#) of [2mL Qiagen Catalog #19131](#) and start crushing the tissue with a micropestle. Keep the pestle in the tube.
- 8 [Buffer ATL \(tissue lysis](#)
Pipet [100 µL](#) of [buffer\) Qiagen Catalog #19076](#) over the micropestle to ensure all tissue remains in the digestion. Transfer the pestle into the prepared

☒ 1.5 mL Safe Lock

Tube Eppendorf Catalog #0030120086

for storage.

- 9 Mix well by vortexing. Incubate tubes at 900 rpm, 56°C, 06:00:00 . Make sure the timer is on infinity.
- 10 After 06:00:00 , centrifuge the tubes and add 20 µL of
☒ Proteinase K,
2mL Qiagen Catalog #19131 to the remaining tissue. Use the pestle to further crush tissue parts. Pipet a second volume of 180 µL
☒ Buffer ATL (tissue lysis
buffer) Qiagen Catalog #19076 over the pestle to ensure all tissue remains in the digestion.
- 11
- Vortex well. Make sure all pieces of tissue are in solution and continue digestion overnight at 56 °C . Leave to digest 900 rpm, 56°C, 40:00:00 . (Note: Tests with extraction after 24 and 40 hours show that DNA amount doubles within the additional 16 hours).

Pre-extraction preparations

- 12 ☒ AL lysis
Preheat buffer Qiagen Catalog #19075 at 56 °C . This step prevents the
☒ AL lysis
formation of precipitates when buffer Qiagen Catalog #19075 is added to the digestion.
- 13 Prepare a 1:1 dilution of ☒ Ultrapure Distilled, Nuclease Free Water Contributed by users
☒ AE
and buffer Qiagen Catalog #19077 . Preheat at 37 °C .
- 14 ☒ 1.5 mL Safe Lock
Prepare and label Tube Eppendorf Catalog #0030120086 for each

extraction (2 for 2 elutions)

Digestion

- 15 Once all tissue is properly digested, in case there are undigested parts (bones, feathers, etc)^{5m}, centrifuge the lysate at max speed for **00:05:00**, and transfer the lysate into a new tube. Else, continue with the next step.

15.1 Note: it is possible that with increased digestion time bigger pieces of tissue are becoming visible. This is likely due to increasing hydration/swelling of tissue pieces with prolonged digestion, and exactly what you want)

- 16 Add **4 µL** **RNase** **A Qiagen Catalog #19101**. Mix by vortexing, and incubate at **Room temperature** for **00:10:00**. 10m

- 17 Vortex for **00:00:15**. Add **400 µL** of preheated **56 °C** **AL lysis buffer Qiagen Catalog #19075**. Mix very quickly. A coagulate will form. 15s
Dissolve this in the heat block at **900 rpm, 56°C, 00:10:00**

Extraction

8m


- 18 Turn temperature on the heat block down to **37 °C**. Preheat 1:1 diluted Buffer AE at **37 °C**.



- 19 Add **400 µL** ice-cold (**-20 °C**) **Ethanol, absolute 99.8% Contributed by users Catalog #10342652**. Mix very well by vortexing. Incubate at room temperature for 10 min.


- 20 Apply **650 µL** of the lysate to the column. Centrifuge at **8000 rpm, 00:01:00**. Discard flow-through, dry the collection tube, and place the column back into the collection tube. 1m



- 21 Apply the rest of the lysate. Centrifuge at **8000 rpm, 00:01:00**. Discard flow-through and 1m


collection tube.





- 22 Place the column into a new collection tube. Add  **500 µL** 1m



 [Buffer AW1, Wash buffer \(1\),](#)
[concentrate Qiagen Catalog #19081](#) . Centrifuge at
 **8000 rpm, 00:01:00** . Discard flow-through and collection tube.

- 23 Place the column into a new collection tube. Add  **500 µL** 1m

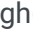



 [Buffer AW2, Wash buffer \(2\),](#)
[concentrate Qiagen Catalog #19072](#) . Centrifuge at
 **8000 rpm, 00:01:00** . Discard flow-through, dry collection tube, and place column back into it.

- 24 Centrifuge at  **14000 rpm, 00:02:00** to dry the membrane. Make sure no ethanol is left 2m
anywhere in or on the column. Discard the collection tube.

- 25 Place the column into a labeled
 [1.5 mL Safe Lock](#)
[Tube Eppendorf Catalog #0030120086](#) . Apply  **20 µL** of
 [AE](#)
preheated ( **37 °C**) 1:1 diluted [buffer Qiagen Catalog #19077](#) to the column.

- 26 Put the column into the thermoshaker  **0 rpm, 37°C, 00:07:00** Then spin at 2m
 **8000 rpm, 00:02:00** .

- 27 

If you expect high(er) amounts of DNA, repeat steps  **go to step #25** and
 **go to step #26** with another  **20 µL** of 1:1 diluted
 [AE](#)
[buffer Qiagen Catalog #19077](#) .



Individual	Collection year	Toe pad	
		ng/mg	Average size (bp)
OENLRSSCH_BMNH-KE-1902-1965.M.12161	1902	528	282
Museum_43	1904	612	279
Museum_58	1913	340	278
OENPICCAP_ZMUC-AF-HAI-1949-29495	1949	621	344
OENPICOP1_ZMUC-AF-MAI-1949-29578	-	808	327
Museum_48	1952	1033	237
OENPHI_YPM-SO-TOG-1954-ORN035210	1954	709	313
OENLRVAU_AMNH-SML-ERI-1954-461151	-	443	239
MYRTHO_YPM-AN-BUN-1958-ORN95640	1958	791	321
Museum_46	1959	1854	359
Museum_44	1960	886	325
Museum_45	-	891	307
Museum_47	-	646	264
Museum_59	-	459	322
Museum_60	-	990	363

Examples of DNA yield per milligram input tissue for samples of varying age.