



NOV 02, 2023

OPEN ACCESS



DOI:
dx.doi.org/10.17504/protocols.io.bp2l69rn5lqe/v1

Protocol Citation: Sean M Cook, Jennifer Jones, Joshua A Welsh 2023. Resource 1: Scatter Detector Setting Incrementation for FCMPASS. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.bp2l69rn5lqe/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

Created: Jan 30, 2023

Last Modified: Nov 02, 2023

Resource 1: Scatter Detector Setting Incrementation for FCMPASS

Forked from [Flow Cytometer Fluorescence Voltration for FCMPASS](#)


In 1 collection

Sean M

Cook¹, Jennifer Jones¹, Joshua A Welsh¹

¹Translational Nanobiology Section, Laboratory of Pathology, Center for Cancer Research, National Cancer Institute, National Institutes of Health

Translational Nanobiology Section

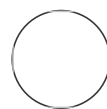

Jennifer Jones

DISCLAIMER

This protocol summarizes key steps for a specific type of method, which is one of a collection of methods and assays used for EV analysis in the NCI Translational Nanobiology Section at the time of submission of this protocol. Appropriate use of this protocol requires careful, cohesive integration with other methods for EV production, isolation, and characterization.

ABSTRACT

Flow cytometry (FCM) is a common extracellular particles (EPs), including viruses and extracellular vesicles (EVs), characterization method. Frameworks such as MIFlowCyt-EV exist to provide reporting guidelines for metadata, controls, and data reporting. However, tools to optimize FCM for EP analysis in a systematic and quantitative way are lacking. Here, we demonstrate a cohesive set of methods and software tools that optimize FCM settings and facilitate cross-platform comparisons for EP studies. We introduce an automated small particle optimization (SPOT) pipeline to optimize FCM fluorescence and light scatter detector settings for EP analysis and leverage quantitative FCM (qFCM) as a tool to further enable FCM optimization of fluorophore panel selection, laser power, pulse statistics, and window extensions. Finally, we demonstrate the value of qFCM to facilitate standardized cross-platform comparisons, irrespective of instrument configuration, settings, and sensitivity in a cross-platform standardization study utilizing a commercially available EV reference material.

PROTOCOL integer ID:
76100


**Funders
Acknowledgement:**


NIH
Grant ID: ZIA BC011502

NIH
Grant ID: ZIA BC011503


NIH
Grant ID: 4UH3TR002881-03

MATERIALS

 FluoSpheres Carnboxylate Thermo Fisher Scientific Catalog #F8803 , referred to as 'FluoSpheres'


 5 mL Round-bottom tube Corning Catalog #352052 , referred to as 'FACS tubes'

 DPBS Thermo Fisher Scientific Catalog #14190144 , referred to as 'DPBS'

 QbSure Cytex Catalog #B7-10005 , referred to as 'QbSure'

Sample preparation

10s

- 1 Vortex FluoSpheres bottle on a high setting for  00:00:05 .

5s

Note

This stock bottle is 2% solids containing 110 nm beads, equal to $\sim 2.73 \times 10^{13}$ / mL


Note







This method in principle is compatible with any fluorescent nanosphere that is bright enough to be fully resolved when using a fluorescent trigger.

- 2 Create 3 mL of 5×10^6 p/mL solution of FluoSpheres in a FACS tube.

Note

An observation from our protocol development is that using low protein binding tubes for this step may negatively impact the efficacy of this protocol as it can result in excess unbound fluorophore from the beads increasing background noise leading to excessive event rate.

- 2.1 Pipette  2998.90 μ L of DPBS into a FACS tube. Label this tube '1E10 Intermediate'.

- 2.2 Pipette  1.10 μL stock FluoSpheres into the tube labelled '1E10 Intermediate'. Reverse pipette the tube to mix.
- 2.3 Pipette  990 μL of DPBS into a FACS tube. Label this tube '1E8 Intermediate'. Pipette  10 μL "1E10 Intermediate" into the tube. Reverse pipette the tube to mix.
- 2.4 Pipette  2850 μL of DPBS into a FACS tube. Label this tube 'Beads. Pipette  150 μL "1E8 Intermediate" into the tube. Reverse pipette the tube to mix. This should result in a concentration of $\sim 5\text{E}6$ p/mL of FluoSpheres, and this tube will be used for acquisition.
- 3 Pipette  500 μL of DPBS into a FACS tube. Label this tube 'DPBS'.

Cytometer Setup

- 4 Ensure cytometer is clean and that -Height and -Area statistics are set to be collected on all parameters and that all parameters are on.
- 5 On the Cytex Aurora, set window extension to 0. On the CytoFLEX platform turn on 'High Acquisition Mode'.


Note

To set the window extension to 0 on the Cytex Aurora, navigate to the 'Lasers' tab under the sample acquisition settings and set the window extension to 0.

To change the acquisition mode on the CytoFLEX S, click on the 'Advanced' menu on the top-of-the-screen ribbon. Next click on 'Event Rate Settings' menu, change the acquisition mode to 'High', and press 'Ok'.

- 6 Create a histogram plot with the FITC height parameter (CytoFLEX: FL1-H | B-525-40, Aurora: B2 | B-524.5-17) on the X-Axis and make sure it is plotted on a log-scale.

- 7 Create a histogram plot with (405 nm) violet SSC-A (CytoFLEX: SSC_1-A, Aurora: SSC-A) on the X-Axis and make sure it is plotted on a log-scale.

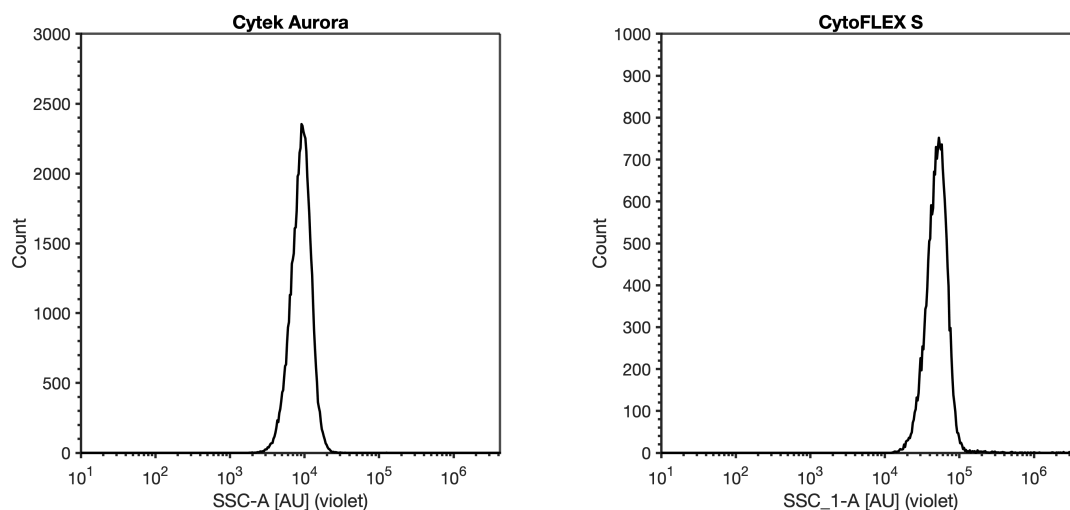

Example of FluoSpheres acquired on Cytek Aurora and CytoFLEX S. Data has been gated to remove the noise population.

- 8 Set the cytometer triggering threshold to the FITC parameter (CytoFLEX: FL1 | B-525-40, Aurora: B2 | B-524.5-17) . **All samples should be acquired with the lowest flow rate, typically ~10-15 $\mu\text{L min}^{-1}$.**

Note

Cytometer Voltage/Gain and threshold settings are subjective due to their dependency on alignment, and the scatter filters in place, amongst other variables. The following are guide values to start with and may need adjustment for optimal acquisition.


Beckman Coulter, CytoFLEX [405 nm OD0 filter, 488 nm OD2 filter]

- Threshold FL1 | B-525-40-H = 550;
- FL1 | B-525-40 Gain = 500

Cytek Bioscience, Aurora [405 nm OD0 filter, 488 nm OD2 filter]

- Threshold B2 | B-524.5-17-H = 500;
- B2 | B-524.5-17 Gain = 1750

- 9 Acquire the 'DPBS' tube while viewing the FITC histogram plot from [go to step #6](#). Adjust the detector gain or trigger threshold until the instrument noise is being **acquired at ~1000 events/sec.**

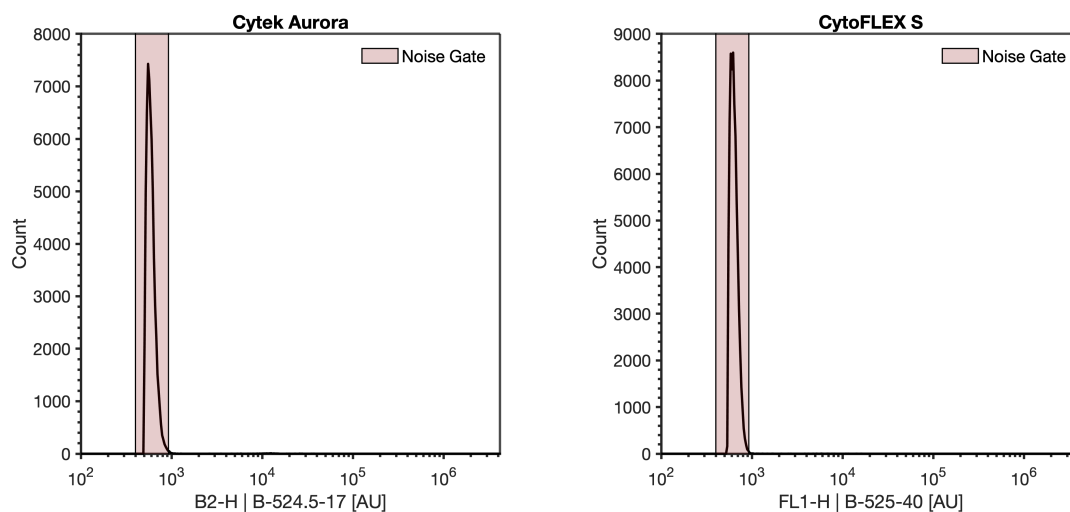

Acquisition of DPBS on the Cytek Aurora and CytoFLEX S, with the Aurora triggering on instrument noise on the B2-H channel at gain 1750 and threshold 500 and the CytoFLEX triggering on instrument noise on the FL1-H channel at gain 500 and threshold 550.


Note

The instrument noise floor is distinct from detected background events in sheath as it has a sharp increase. In a system with debris there may be a tail that elongates out of this sharp peak.

- 9.1 Recording this noise is not necessary as this step is identifying optimal settings.

- 10 Acquire the "Beads" tube from [go to step #2](#). Using the plot from [go to step #6](#), ensure the FluoSpheres are visible on the FITC and violet SSC parameter. Use the FITC trigger

settings identified in [⇒ go to step #9](#)

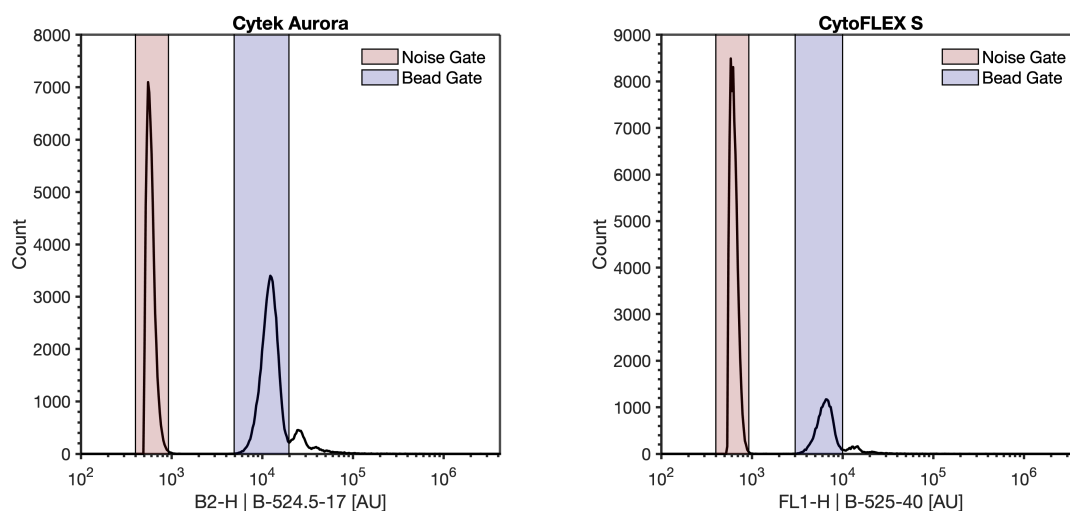

Acquisition of 110 nm FluoSpheres in DPBS on the Cytek Aurora and CytoFLEX S, with the Aurora triggering on instrument noise on the B2-H channel at gain 1750 and threshold 500 and the CytoFLEX triggering on instrument noise on the FL1-H channel at gain 500 and threshold 550.

Note

The total event rate when acquiring the FluoSpheres should not exceed 6000 events/sec. If the event rate is higher than 6000 events/sec, further dilution of the FluoSpheres is necessary prior to acquisition.

- 11 Draw a gate around the FluoSpheres on the FITC parameter. Label this gate 'Bead Gate'
- 12 Adjust the stopping criteria of the instrument to record until at least 2,000 events are acquired on 'Bead Gate' drawn in [⇒ go to step #11](#).

Performing Voltration

- 13 Voltration can now be performed by recording the 'Beads' tube at multiple light scatter detector

gains, leaving the trigger threshold and fluorescent gains consistent. It is recommended that a recording of at least 10 light scatter detector settings is taken. Including more increments within a volttration will result in being more confident of the subsequent optimal detector settings.

 Acquisition Template.xlsx

Note

For flow cytometers with avalanche photodiodes its is recommended that the detector settings have more incrementation at lower gains than higher gains while for instruments with photomultiplier tubes they should be spaced evenly. See template for example of settings for each tube analyses

Example Gain Volttration for CytoFLEX & Aurora

1. 50
2. 100
3. 200
4. 300
5. 400
6. 500
7. 750
8. 1000
9. 1250
10. 1500
11. 2000
12. 2500
13. 3000

13.1



To ensure accurate data analysis, the FluoSpheres must separate from the noise population on the light scatter parameter on at least two of the selected gains. If they do not, repeat this protocol using a 200 nm FluoSphere population.