

VERSION 2
OCT 11, 2023

OPEN ACCESS



DOI:
dx.doi.org/10.17504/protocols.io.5qpvorz99v4o/v2

Protocol Citation: Natalia Martagón, Gurdrun Kliem 2023. Wnt-3a and R-spo1 conditioned media reporter assay. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.5qpvorz99v4o/v2> version created by Natalia Martagón

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working


Wnt-3a and R-spo1 conditioned media reporter assay V.2

Natalia

Martagón¹, Gurdrun Kliem²

¹Charité Universitätsmedizin Berlin; ²Robert Koch Institute

Rezvani Lab Charite

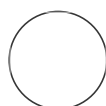

Natalia Martagón

ABSTRACT

Protocol designed to measure the activity of Wnt-3a or R-spondin-1 (Rspo1) conditioned media.

A reporter HEK cell line expressing luciferase under Wnt-3a stimulation is cultured with conditioned media followed by cell lysis and a luciferase reporter assay. Activity is compared to previous media batches or references.

IMAGE ATTRIBUTION

Created with BioRender.com

MATERIALS

Assay-reagent solution:

- [M] 20 millimolar (mM) Tricine (MW 179,2)
- [M] 1.07 millimolar (mM) Mg Carbonate x 5 H₂O (MW 485,7)
- [M] 2.67 millimolar (mM) Mg Sulfate x 7 H₂O (MW 246,5)
- [M] 0.1 millimolar (mM) EDTA (MW 372,2)

Preparation (warm up for preparation):

Dissolve:

3.584 g Tricine

0.519 g Mg-CO₃

in ca. 850 mL water

Add 54 µL of the [M] 0.4 Molarity (m) MgSO₄ pH 8 solution

adjust pH to pH 07.8

bring volume to 1 L Store at RT in glass bottles

Lysis solution

Created: Oct 11, 2023

Last Modified: Oct 11, 2023

PROTOCOL integer ID:
89164

Keywords: Wnt-3a,
conditioned media, activity
assay, luciferase reporter
assay, R-spondin conditioned
media

Funders
Acknowledgement:

Deutsche
Forschungsgemeinschaft
(DFG)
Grant ID: 437638965

luciferase solution stocks preparation

Co Enzym A: .

Sigma C3019-100 mg

Dissolve 100 mg in 1,27 ml in water. Make aliquotes of 100 µL

Luciferin

Dissolve 25 mg D-Luciferin-Na salt in 4,135 ml water. [Stock]=20 mM. Make aliquotes of 500 µL

0.5 M CDTA pH 8:

Dissolve 8,7 g in 50 ml. Buffer the solution first with NaOH-pearls and fine adjustment with liquid NaOH ca. 10 M

0.5 M MgSO₄-Stock, pH 8:

Dissolve 6,15 g in 50 ml in water

1M EDTA-Stock:

Dissolve 29,2 g in 100 ml water

0.5 M Tris-Stocksolution, pH 7,8:

Dissolve 15,14 g Tris-Base(MW 121,14) in 250 ml water. Adjust pH with phosphoric acid

1M DTT-Stock:

Dissolve 1,54 g (MW 154,3) in 10 ml water

Lysis solution stocks

125 mM Tris-Phosphate

10 mM DTT

10 mM CDTA

50% Glycerin

5 % Triton X 100

5x Lysis-Solution preparation

12,5 ml Tris-Phosphat 0,5 M, pH 7,8

500 µL DTT 1 M

1 ml CDTA 0.5 M, pH 8

25 g Glycerin

2,5 g Triton X 100

Bring to a total vol. of 50 ml with water. Make 10 ml aliquotes and store at -20°C

HEK-medium

DMEM F12 HEPES (Gibco Cat. No. 11330),

20% FCS

200 µg/ml G418

material for cell counting

Microplatte for luminescence readers

Greiner Bio-one Cat. no. 655094

BEFORE START INSTRUCTIONS

Be sure to have access to:

- Luminometer with or without injection system.
- 96 well plates appropriate for luminescence signal measurements
- Reagents for culturing reporter cell line. See reference (1)

*reagents might take a while to gather and prepare but once prepared they last for a long time

*time to produce conditioned media around 3 weeks. See reference protocols (2), (3).

Day 1: Seeding of Hek 293 STF cells

- 1 Culture HEK 293 STF CRL-3249TM (from now on referred as HEK-STF) according to the company specifications (1) until confluent and not too many passages old.

Seed 3 wells of HEK-STF cells per sample and following controls:

- negative control: conditioned media from L-cells (not transfected with luciferase construct) (4)
- positive control: previous batch with known activity, HEK-STF cells with recombinant Wnt-3a or agonist stimulation
- control lysates as blank for luminescence: HEK-SFT cells with HEK-medium only


- 1.1 Start with one almost confluent T75 culture bottle of HEK-STF cells

- 1.2 1x wash with DPBS: take out medium, add 5 mL DPBS, turn gently the bottle, take out DPBS



- 1.3 Detach cells with 1-2 mL Trypsin/EDTA ( 37 °C) and transfer to a conical tube with 8-9 mL HEK-medium

5m


- 1.4 Centrifuge at 100-200 g for  00:05:00
count cells (Neubauer chamber or automatic cell counter)

5m



- 1.5 count cells (Newbauer chamber or automatic cell counter)

5m

- 1.6 Seed cells in a 24 well plate. HEK-STF :
24 well Platte (0.05×10^6 cells/well). $1,3 \times 10^6$ cells / 13ml for the whole plate
cover each well with  0.5 mL of HEK medium.


15m




Day 2: cell stimulation


1d

- 2 Aspirate medium and discard.

30m

Add a total of  500 μ L of HEK medium with the following amount of conditioned medium (CM) to test:

- For R-spondin CM:  12.5 % volume Wnt-3a CM +  2.5 % volume Rspo1 CM (3)
- For Wnt3a CM .  50 % volume

- 2.1 Incubate at stadard culture conditions for aprox.  24:00:00

1d



Day3 : luicferasse reaction

20m

- 3** **Luiferase solution**
Prepare on the day luciferase solution. Keep in the fridge until ready to use
When using a plate reader with automatic injection, calculate the extra volume needed for it (ca 3 mL)
fill to final volume with **Assay reagent** solution

A	B	C	D
Reagent	units	[stock]	[working]
DTT	M	1	0,033
CoA	mM	10	0,266
Luciferin	mM	20	0,467
ATP	mM	100	0,633

- 4** Prepare 1x Lysis medium from the stock with destiled water



20m

- 5** **Cell lysis**
Aspirate medium from HEK-STF cells

2m

- 5.1** Add  150 µL /well of 1x Lysis buffer.

2m

- 6** Leave the plate for  00:20:00 at  Room temperature . (Note: pippeting seems to make cell aggregates and bubbles)



20m

- 7** Check under the microscope that most cells are lysed



8 Shake gently for  00:05:00 on plate shaker

5m

9 Transfer  20 μL of each sample to a 96 well plate for luminescence read.

10m

Avoid bubbles or cell clumps. Recommended pipetting scheme as in the 24 well plate of HEK-STF culture-

10 If using a plate reader without automatic injection:

20m

Take to the plate reader with luci. solution, multichanel pipette, and reagent reservoir

Once the plate reader is set, add  100 μL of luciferase solution per well and read.

For accuracy, use a multichanel pipette and add the solution covering not one sample at a time but one of the triplicates of all samples at one time.

11



(signal decreases ca. 20% in the first 10 minutes)

12 Check that the relative luminescence units from negative controls are orders of magnitude lower than test samples. Average values from triplicates.

Values below more than 50-60% lower than a working batch of conditioned media are considered of poor quality.

Examples of equipment and methods: Berthold Tristar-GAS ISRE Luciferase Assay, SpectraMax i3x-SpectraMax Glo Steady-Luc Reporter Assay