



APR 17, 2023

OPEN ACCESS

DOI:
dx.doi.org/10.17504/protocols.io.x54v9dw3qg3e/v1

Protocol Citation: Linda Armbrecht 2023. Ambrecht et al. 2020: An optimized method for the extraction of ancient eukaryote DNA from marine sediments. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.x54v9dw3qg3e/v1>

MANUSCRIPT CITATION:

CITATION


Linda Armbrecht, Salvador Herrando-Pérez, Raphael Eisenhofer, Gustaaf M. Hallegraeff, Christopher J. S. Bolch, Alan Cooper (2020). An optimized method for the extraction of ancient eukaryote DNA from marine sediments. Molecular Ecology Resources.

LINK
[10.1111/1755-0998.13162](https://doi.org/10.1111/1755-0998.13162)

Ambrecht et al. 2020: An optimized method for the extraction of ancient eukaryote DNA from marine sediments

Linda Armbrecht¹

¹School of Biological Sciences, Faculty of Sciences, Australian Centre for Ancient DNA, The University of Adelaide, Adelaide, SA, Australia

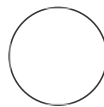

Mark Louie Lopez

ABSTRACT

Four combinations of sedaDNA extraction treatments using marine sediments collected at a water depth of 104 m off Maria Island in Tasmania are compared. These methods contrast frozen versus refrigerated sediment, bead-beating induced cell lysis versus ethylenediaminetetraacetic acid (EDTA) incubation, DNA binding in silica spin columns versus in silica-solution.

All four methods worked to varying degrees; see paper for recommended shotgun library preparation

ATTACHMENTS

[Ambrecht et al. 2020 - MolEcolRes sedaDNA extractions.pdf](#)

GUIDELINES

DNA extractions:

- **Bead-beating + spin column (DNeasy PowerLyzer PowerSoil Kit, Qiagen; “Kit”)**
- **Bead-beating + liquid silica in QG Buffer (“Si4” and “Si20”)**
- **EDTA + MinElute (“EDTA”)**
- **EDTA + bead-beating + liquid silica in QG Buffer (“Combined”, or “Com”)**

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Feb 13, 2023

Last Modified: Apr 17, 2023

PROTOCOL integer ID:
76894



Keywords: ancient DNA, diatoms, dinoflagellates, haptophytes, Maria Island, metagenomics, plankton, seafloor, Tasmania

BEFORE START INSTRUCTIONS


Sediment processing and pretreatment:

Core section processing, sedaDNA extractions and sequencing library preparations took place at ACAD's ultraclean forensic facilities following aDNA decontamination standards (Willerslev & Cooper, 2005). We placed the three sediment core sections into zip-lock bags sterilised with UV light and manually homogenized them for ~5 min. From each section, two 1 cm³ subsamples were transferred into two separate 15 ml centrifuge tubes using a sterile disposable spatula. One subsample was kept at 4°C and the other at -20°C for one month. The samples were prepared for the different extraction methods in a glove box decontaminated (3% bleach) between consecutive subsamples.



Method 1: Bead-beating + spin column (DNeasy PowerLyze... 3m 25s

- 1 This technique was applied to  0.25 g of sediment subsamples stored at  4 °C, following the **manufacturer's protocol** with the some modifications:

- 1.1 **TRANSFER** sediment into individual bead-tubes using a disposable, sterile spatula

- 1.2 **APPLY** bead-beating in three runs of  00:00:20 with  00:00:05 breaks using a Precellys 24 homogenizer

CENTRIFUGE at  10319 undetermined for  00:00:30

- 1.3 **RETAIN** all optional  00:05:00 incubation steps at  4 °C per the kit's protocol

- 1.4 ELUTE DNA in 80 µL of Buffer EB instead of the customary C6 solution and store at -20 °C

Method 2: Bead-beating + liquid silica in QG Buffer ("SI4" a.. 1h 32m

- 2 This lysis process was applied to 0.25 g of the subsamples stored at both 4 °C and -20 °C

- 2.1 FOLLOW the same protocol as described in Section 1 down to step 10 of the manufacturer's instructions (addition of Solution C3 and subsequent centrifugation)

AFTER this step...

- 2.2 TRANSFER the supernatant to 15 mL centrifuge tubes containing a DNA-binding buffer

Note

DNA Binding Buffer

- 100 µL silica-solution (Sigma Aldrich)
- 3 mL modified Buffer QG (2.7 mL Buffer QG, 46 µL H₂O, 39.08 µL Triton X-100, 24.66 Molarity NaCl, and 164.5 Molarity NaOAc (Brotherton et al., 2013)

- 2.3 STIR on a rotary mixer for 00:00:00 at Room temperature 5m


CENTRIFUGE at 4500 rpm for 00:05:00



DISCARD supernatant

RESUSPEND pellet in 900 µL of DNA-binding buffer




- 2.4 RE-CENTRIFUGE at 14000 rpm for 00:01:00 16m


DISCARD the supernatant

WASH the pellet twice in  80 % ethanol

DRY pellet for  00:15:00 at  37 °C



RESUSPEND pellet in  80 µL Buffer EB


2.5 FOLLOWING incubation for  00:10:00 at  50 °C, centrifuge at  14.000 rpm for  00:01:00 11m




STORE the supernatant (free of silica) in a sterile Lo-bind tube (Eppendorf) at  -20 °C

Method 3: EDTA + MinElute (“EDTA”)

3m

3 This technique was applied to  0.25 g sediment subsamples stored at  4 °C following **Slon et al. (2017)** with minor modifications

3.1 ADD  1 mL of ethylenediaminetetraacetic acid (EDTA) to the sediment in a 2 mL screw-cap tube

PLACE samples on a rotary mixer and mix at  25 rpm,  Overnight at  Room temperature

3.2 CENTRIFUGE at  13.000 rpm for  00:03:00 3m



PURIFY the DNA using the MinElute Kit (Qiagen) as **per the manufacturer's instructions**

3.3 BIND DNA using the kit's spin column

ELUTE the DNA in  80 µL of Buffer EB


Method 4: EDTA + bead-beating + liquid silica in QG Buffer ..

3m


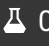

4 INCUBATE  0.25 g of three frozen sediment subsamples in EDTA overnight as in **step 3.1 of the above section "EDTA + MinElute (“EDTA”)**, EXCEPT use only  0.75 mL to keep volumes consistent with a subsequent step (below)

4.1 **CENTRIFUGE** at  13.000 rpm for  00:03:00

3m


KEEP supernatant at  4 °C

PROCESS pellet separately using bead-beating and DNA purification, as in **Method 2** (above)

4.2 **RECOMBINE** the resulting  0.75 mL DNA-solution purified from the pellet (step 10 of DNeasy Kit protocol) with  0.75 mL EDTA supernatant to make  1.5 mL total

ADD  6 mL modified QG buffer with  100 µL liquid silica

PROCEED as described in **Method 2** (above)

4.3 **ELUTE** the DNA in  100 µL Buffer EB