



Aug 17, 2021

# Optogenetic experiments with iLID system

Andrés Guillén-Samander<sup>1,2</sup>, Pietro De Camilli<sup>1,2</sup>

<sup>1</sup>Departments of Neuroscience and of Cell Biology, Howard Hughes Medical Institute, Program in Cellular Neuroscience, Neurodegeneration and Repair, Yale University School of Medicine, New Haven, Connecticut 06510, USA;

<sup>2</sup>Aligning Science Across Parkinson's (ASAP) Collaborative Research Network, Chevy Chase, MD, 20815

1 Works for me



Share

dx.doi.org/10.17504/protocols.io.bvgvn3w6

Andrés Guillén-Samander

## ABSTRACT

This protocol details experiments with the iLID optogenetic system as performed to acutely recruit Miro to mitochondria in <https://doi.org/10.1083/jcb.202010004>.

## ATTACHMENTS

[dn3xbgtzx.pdf](#)

## DOI

dx.doi.org/10.17504/protocols.io.bvgvn3w6

## PROTOCOL CITATION

Andrés Guillén-Samander, Pietro De Camilli 2021. Optogenetic experiments with iLID system. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.bvgvn3w6>

## KEYWORDS

Optogenetic experiments, iLID system, ASAPCRN

## LICENSE

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

## CREATED

Jun 02, 2021

## LAST MODIFIED

Jul 02, 2021

## OWNERSHIP HISTORY

Jun 02, 2021 Urmilas

Jun 08, 2021 Andrés Guillén-Samander

## PROTOCOL INTEGER ID




50421

## Optogenetic experiments with iLID system

3d 12h

- 1 Plate the COS-7 cells for imaging and transfect each MatTek dish (see Cell Culture, Transfection and Imaging Protocol) with the plasmids containing the bait (construct encoding iLID) and the prey (construct encoding SspB peptide).




As an example, in the case of Guillen-Samander A, Leonzino M et al 2021 (<https://doi.org/10.1083/jcb.202010004>), the plasmids and amounts transfected were:

-  **0.25 ug** mCh-Miro1ΔTM-SspB
-  **0.5 ug** Venus-Mito-iLID
-  **0.8 ug** VPS13D<sup>Δ</sup>Halo

**Suggestion:** When dealing with an optogenetic pair of constructs encoding a membrane anchored protein and a cytosolic one, we recommend transfecting them in a ratio of 2:1 for membrane: cytosolic proteins.



3d 12h

After adding the transfection mix, incubate the dishes at  **37 °C** in the dark for  **36:00:00** -  **48:00:00**.

3 Just before imaging, remove the growth medium and replace with Live Cell Imaging solution (Life technologies).

4 Image the cells at  **37 °C** and 5% CO<sub>2</sub>.

Warning: Remember to keep the cells in the dark while handling them and as much as possible while changing the medium. Exposure to light could preactivate the optogenetic system which could affect the interpretation of results.

Step 4 includes a Step case.

**For whole cell activation:**

**For localized activation:**

step case

**For whole cell activation:**



Carry out imaging in an Andor Dragonfly system equipped with a PlanApo objective (63X, 1.4NA, Oil) and a Zyla sCMOS camera.

6 Identify the cells to be imaged by scanning the dish looking for red fluorescence in the 564nm channel.

Remember that looking for green fluorescence would preactivate your sample.

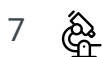


Image the cells at a rate of 0.5Hz. Acquire the 5 frames before the activation, and achieve the activation with a single 200ms pulse of the 488nm laser. Acquire the 90 frames after activation.