



DEC 29, 2023


DoTA-seq V3.1

Forked from [DoTA-seq V3.1](#)

freeman.lan¹

¹UW Madison

DoTA-seq

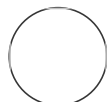

freeman.lan

OPEN ACCESS



DOI:
dx.doi.org/10.17504/protocols.io.81wgbxx3ylpk/v1

Protocol Citation: freeman.lan 2023. DoTA-seq V3.1. [protocols.io](https://dx.doi.org/10.17504/protocols.io.81wgbxx3ylpk/v1)
<https://dx.doi.org/10.17504/protocols.io.81wgbxx3ylpk/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Aug 24, 2023

Last Modified: Dec 29, 2023

ABSTRACT

This protocol describes the process of DoTA-seq generating a single cell sequencing library from a cell suspension. This workflow can be performed in two days, with the PCR step happening overnight. Before beginning this workflow make sure to have:

1. The necessary microfluidics devices prepared and ready to go
2. The multiplex DoTA-seq target primers validated to work together without generating large molecular weight primer dimers.

Please read the publication for further details.

GUIDELINES

Strongly recommend all pre-PCR steps (setting up reagents, washing gels) to be done in a PCR Clean hood. This has two purposes:

1. Reduce PCR contamination of templates which can strongly effect single-cell PCR reactions.
2. Reduce dust contamination of reagents which can clog devices and cause failures.

MATERIALS

ddPCR Supermix for probes (no dUTP) BioRad
Sciences Catalog #1863024

MetaPolyzyme Sigma
Aldrich Catalog #MAC4L-5MG

Lysozyme from chicken egg white Sigma
Aldrich Catalog #L6876



HFE 7500 Perfluorinated
Oil



Perfluorooctanol **Sigma**
Aldrich Catalog #370533



TCEP-HCl Gold
Biotechnology Catalog #TCEP



NN'-Bis(acryloyl)cystamine **Santa Cruz**
Biotechnology Catalog #sc-215506



Ammonium persulfate Catalog #A3678



Acrylamide P212121



TEMED (Tetramethyl-ethylenediamine) **Sigma-**
aldrich Catalog #T9281



Biorad Evagreen Droplet Oil **BioRad**
Sciences Catalog ##1864005



DNA Clean & Concentrator™-5 **Zymo**
Research Catalog #D4003



DNA Clean & Concentrator™-5 **Zymo**
Research Catalog #D4003



Axygen® 0.2 mL Maxymum Recovery® Thin Wall PCR
Tubes **Corning** Catalog #PCR-02-L-C



NEBNext Library Quant Kit for Illumina - 100 rxns **New England**
Biolabs Catalog #E7630S



SYBR Green **Thermo Fisher**
Scientific



Proteinase K solution, 20 mg ml -
1 **Ambion** Catalog #AM2546



Lysozyme from chicken egg white **Sigma**
Aldrich Catalog #L6876



Pre-injection buffer 10mM HEPES pH 7.5 Pluronic
0.1%

Safety information

Unpolymerized Acrylamide is toxic, handle with care and dispose according to regulations

PROTOCOL MATERIALS

⊗ DTT Merck MilliporeSigma (Sigma-Aldrich) Catalog #D0632 Step 30

⊗ TEMED (Tetramethyl-ethylenediamine) Merck MilliporeSigma (Sigma-Aldrich) Catalog #T9281

Materials, Step 10

⊗ Proteinase K solution, 20 mg ml⁻¹
1 Ambion Catalog #AM2546

Materials, Step 24

⊗ Pre-injection buffer 10mM HEPES pH 7.5 Pluronic 0.1% Materials

⊗ TCEP-HCl Gold Biotechnology Catalog #TCEP Materials, Step 38

⊗ HFE 7500 Perfluorinated Oil Materials, Step 16

⊗ Acetone Step 17

⊗ TE Buffer (Tris 10mM EDTA 1mM pH 8) Step 24

⊗ SPRIselect reagent kit Beckman Coulter Catalog #B23317 Step 40

⊗ MetaPolyzyme Merck MilliporeSigma (Sigma-Aldrich) Catalog #MAC4L-5MG

Materials, Step 21

⊗ High Sensitivity D1000 ScreenTape Agilent Technologies Catalog #5067-5584

Step 41

⊗ DNA Clean & Concentrator™-5 Zymo Research Catalog #D4003

In Materials, Materials, Step 39

⊗ PBS EDTA 1mM Tween20 2% w/v Step 19

⊗ TE 2% Tween-20 In 2 steps

⊗ Pre-injection buffer 10mM HEPES pH 7.5 Pluronic 0.1% Step 28

⊗ Biorad Evagreen Droplet Oil Bio-Rad Laboratories Catalog #1864005

Materials

⊗ PBS 0.1% Tween 20 Step 2

✕ Acrylamide P212121 Materials, Step 6

✕ PBS 0.1% Tween20 Step 1

✕ Biorad Evagreen Droplet Oil Bio-Rad Laboratories Catalog ##1864005 Step 10

✕ PBS 0.1% Tween20 Step 23

✕ Perfluorooctanol Merck MilliporeSigma (Sigma-Aldrich) Catalog #370533

In Materials and [2 steps](#)

✕ 10% SDS Bio-Rad Laboratories Catalog #161-0146 Step 24

✕ Axygen® 0.2 mL Maxymum Recovery® Thin Wall PCR Tubes Corning Catalog #PCR-02-L-C

Step 32

✕ NN'-Bis(acryloyl)cystamine Santa Cruz Biotechnology Catalog #sc-215506

Materials

✕ NEBNext Library Quant Kit for Illumina - 100 rxns New England Biolabs Catalog #E7630S

Materials, Step 42

✕ Biorad Evagreen Droplet Oil Bio-Rad Laboratories Catalog ##1864005 Step 7

✕ Biorad Evagreen Droplet Oil Bio-Rad Laboratories Catalog ##1864005 Step 31

✕ ddPCR Supermix for probes (no dUTP) Bio-Rad Laboratories Catalog #1863024

Materials, Step 30

✕ Ammonium persulfate Catalog #A3678 Materials, Step 6


✕ Axygen® 0.2 mL Maxymum Recovery® Thin Wall PCR Tubes Corning Catalog #PCR-02-L-C

Materials

✕ Lysozyme from chicken egg white Merck MilliporeSigma (Sigma-Aldrich) Catalog #L6876

In Materials, Materials, Step 21

✕ Cellbrite Fix 555 Biotium Catalog ##30088 Step 5


 NN'-Bis(acryloyl)cystamine **Santa Cruz**
Biotechnology Catalog #sc-215506

Step 6

 Isopropanol **Step 18**

 Pre-injection buffer 10mM HEPES pH 7.5 Pluronic
0.1%


Step 27

 SYBR Green **Thermo Fisher**
Scientific

Materials, Step 3

 HFE 7500 Perfluorinated
Oil








Step 36

 EDTA (0.5 M), pH 8.0 **Life**
Technologies Catalog #AM9260G

Step 35












Preparing Cells

10m

- 1 Prepare a cell suspension by washing twice in 1mL of  PBS 0.1% Tween20 **Sigma**
Aldrich by spinning down at  5000 x g, 00:01:00
- 2 Resuspend cells in  100 μ L  PBS 0.1% Tween 20 **Sigma**
Aldrich
- 3 Add  1 μ L  SYBR Green **Sigma**
Aldrich 10,000X dye to the cells to stain them
- 4 Count cells using a hemacytometer using the SYBR signal, calculate concentration of the cell suspension.
- 5 (Optional) Stain with  Cellbrite Fix 555 **Sigma**
Aldrich Catalog ##30088 to get a cell membrane/wall stain

Preparing Gel



30m

- 6 Make  200 μL Hydrogel Precursor Solution - Mix together in a tube:
-  100 μL  Acrylamide Sigma Aldrich monomer in water  25 Mass / % volume
-  15 μL  NN'-Bis(acryloyl)cystamine Sigma Aldrich Catalog #sc-215506 in Methanol
-  5 Mass / % volume
-  10 μL  Ammonium persulfate Sigma Aldrich Catalog #A3678  10 Mass / % volume
-  75 μL Cell suspension diluted in PBS (a total of **7e6 cells** to achieve a final concentration of **3.5e7 cells/mL** in the total solution)

Vortex Vigorously to Mix

Generate Gel Droplets

10m

- 7 Prepare and Load the Syringes with the gel sample and  600 μL  Biorad Evagreen Droplet Oil Sigma Aldrich Catalog ##1864005 and connect to the microfluidic devices by following this protocol

Protocol




NAME

Loading Syringes to Inject into Microfluidics Device

CREATED BY

freeman.lan

PREVIEW

- 8 Connect the syringes to the inlets of the **DoTA-seq Step 1** Microfluidics Device and Run the syringe pumps at **500uL/hr** for the gel syringe, and **900uL/hr** for the oil syringe.
- 9 Collect gel droplets for  00:20:00 in a 1mL tube.


20m

Note


Sometimes the initial droplet formation produces polydisperse droplets. In this case, wait 2 min for the bad emulsion to leave the outlet tubing into a waste tube, then begin collecting in the collection tube.



10 Make  200 μL Gel Polymerization Oil - mix together in a tube:

 195 μL

 Biorad Evagreen Droplet Oil Sigma
Aldrich Catalog ##1864005


 5 μL

 TEMED (Tetramethyl-ethylenediamine) Sigma
Aldrich Catalog #T9281

11 Add the Gel polymerization oil to the collected droplets, invert slowly 3 times to mix, and Incubate the tube containing droplets at  37 $^{\circ}\text{C}$ for  00:10:00 to complete polymerization of the gel matrix. 10m

Note




You can now look at the emulsion under the microscope using

 Countess slides Thermo Fisher Scientific Catalog #C10228 to determine the encapsulation ratio of your cells. SYBRGreen and CF555 signal should be concordant and correspond to cells.

Breaking out gels from emulsion

12 Pulse spin the emulsion in a centrifuge to close pack the emulsion and drain the oil to the bottom of the tube.



13 Use a pipette to remove the oil from the bottom of the tube, leaving just the emulsion

14 Add  200 μL  Perfluorooctanol Sigma
Aldrich Catalog #370533 to break the emulsion 1m
Vortex, then Wait  00:01:00 for the emulsion to break.


Pulse spin again and remove the oil in the bottom of the tube with a pipette.

15



16

Add  1000 μ L of  HFE 7500 Perfluorinated Oil Contributed by users to the tube.


10s

Mix by inverting 5 times, Wait  00:00:10 then remove with a pipette.

17

Add  1000 μ L of  Acetone Sigma Aldrich to the tube.

10s

Mix by inverting 5 times, Wait  00:00:10 then remove with a pipette.

The gels should begin to flocculate and dehydrate.

18

Add  1000 μ L of  Isopropanol Sigma Aldrich



10s

Mix by inverting 5 times, Wait  00:00:10 then remove with a pipette.

The gels should dehydrate and become hard.

Note: Do not wait too long as it could cause the gels to irreversibly aggregate into clumps.


19

Resuspend in  1000 μ L  PBS EDTA 1mM Tween20 2% w/v Sigma Aldrich

The gels can be stored at  4 $^{\circ}$ C for several days without changing DoTA-seq results.

Note

You can now look at the gels under the microscope using



 Countess slides Thermo Fisher Scientific Catalog #C10228 to determine the encapsulation ratio of your cells.

For **SYBR Staining**, you should first wash a small aliquot in PBS 2% Tween to remove background SYBR+ oil droplets before visualization on the microscope.

You should see some loss in CF555 signal as the acetone and alcohol wash removes some bacterial membranes.


Lysing Bacteria

- 20 Wash gels 3 times in  1000 μL  1X PBS (Phosphate-buffered saline) Sigma Aldrich (No Tween) by centrifugation at  500 x g, 00:00:30 each time 30s
- 21 Make a Enzymatic Lysis Solution by adding:
-  20 mg  Lysozyme from chicken egg white Sigma Aldrich Catalog #L6876
-  100 μL [M] 1 mg/mL  MetaPolyzyme Sigma Aldrich Catalog #MAC4L-5MG
-  900 μL  1X PBS (Phosphate-buffered saline) Sigma Aldrich (No Tween)
- 22 Resuspend the gels in this lysis solution. Incubate at  37 $^{\circ}\text{C}$ for  02:00:00 2h
- 23 Wash the gels 3 times in  1000 μL  PBS 0.1% Tween20 Sigma Aldrich
- 24 Make a SDS Lysis solution by adding:
-  20 μL [M] 20 mg/mL  Proteinase K solution, 20 mg ml – 1 Ambion Catalog #AM2546
-  100 μL  10% SDS Bio-Rad Laboratories Catalog #161-0146
-  880 μL  TE Buffer (Tris 10mM EDTA 1mM pH 8)
- 25 Resuspend the gels in this SDS Lysis solution, incubate at  55 $^{\circ}\text{C}$ for  01:00:00 1h
- 26 Wash the gels three times in  1000 μL  TE 2% Tween-20
Note: Use **2% Tween 20**, not 0.1% Tween

These gels can be stored at  4 °C in  TE 2% Tween-20 Contributed by users for several days without impacting DoTA-seq results.

Note

You can now look at the gels under the microscope using


 Countess slides Thermo Fisher Scientific Catalog #C10228 to determine the encapsulation ratio and lysis efficiency of your cells.


It is advised to restain with SYBR and CF555 to get best signal. Lysed cells should exhibit SYBR signal but no CF555 Signal.

Barcoding the Cells


7m


27

Wash the gels three times in  1000 µL

 Pre-injection buffer 10mM HEPES pH 7.5 Pluronic 0.1% Sigma Aldrich

28

Resuspend gels in  100 µL

 Pre-injection buffer 10mM HEPES pH 7.5 Pluronic 0.1% Sigma Aldrich

29

Load the gels into a syringe following the protocol described in this excellent visual protocol.

CITATION

Demaree B, Weisgerber D, Lan F, Abate AR (2018). An Ultrahigh-throughput Microfluidic Platform for Single-cell Genome Sequencing.. Journal of visualized experiments : JoVE.

LINK

<https://doi.org/10.3791/57598>

Alternatively, for a simpler version, you can also use a P200 pipette to directly pipette the gel into a syringe backfilled with HFE7500.

30

Generate a PCR Master Mix (This mix gives about ~10,000 cells per library - Scale up as required)

20m

🧪 25 μ L	🔗 ddPCR Supermix for probes (no dUTP) Bio-Rad Laboratories Catalog #1863024
🧪 0.4 μ L	📏 50 micromolar (μ M) P7 Primer mixed with Barrev primer at 5 micromolar (μ M)
🧪 0.4 μ L	📏 50 micromolar (μ M) P5 Primer with appropriate I5 index
🧪 0.2 μ L Variable	📏 10 micromolar (μ M) DoTA-seq multiplex primer mix (10-20nM final concentration per primer)
🧪 0.2 μ L Variable	📏 10 micromolar (μ M) 16S DoTA-seq primers (10-20nM final concentration)
🧪 0.5 μ L Variable	📏 1 picomolar (pM) Freshly diluted from 500pM stock Barcode Oligo
🧪 0.25 μ L	📏 500 millimolar (mM) Single use aliquots
🔗	DTT Merck MilliporeSigma (Sigma-Aldrich) Catalog #D0632

Note

The ratio of 16S to DoTA-seq target primers mix can be varied depending on the relative amplification efficiencies. The best way to determine is to start from equal concentrations, then adjust based on the sequencing results (do most cells contain more 16S reads than target reads?)

Note



IMPORTANT - BARCODE CONCENTRATIONS MAY NEED TO BE MEASURED

Typically, 0.5 μ L of 1pM barcode will give approximately 1 barcode for every 10 droplets. However, it is best to measure the barcode encapsulation rate by making PCR droplets containing the barcodes at the expected dilution and P7 and Barrev primers targeting the barcode for amplification. Visualize the resulting PCR emulsion using SYBRgreen staining under the microscope to obtain the real encapsulation ratio. Typically, the real barcode concentration can be ~5 fold off from the expected concentration based on manufacturer's labelling.

Note

Barcode oligos should always be freshly diluted from 500pM to 1pM before use, as we have found gradual loss of barcodes over time in a 1pM solution.

31 Load the PCR mastermix into the syringe following this protocol

Load  500 μL of  Biorad Evagreen Droplet Oil Sigma Aldrich Catalog ##1864005 into a syringe following this protocol

Protocol



NAME

Loading Syringes to Inject into Microfluidics Device

CREATED BY

freeman.lan

PREVIEW

32 Connect the syringes to the **DoTA-seq Step 2** microfluidics device.

7m

Run the syringe pumps at **200uL/hr** for the gel and PCR mastermix, and **800uL/hr** for the oil syringe. Collect droplets in an



Axygen® 0.2 mL Maxymum Recovery® Thin Wall PCR TubesSigma Aldrich Catalog #PCR-02-L-C

for  00:07:00 for every  25 μL of PCR mastermix or until the PCR mastermix runs out.

33 Use a pipette to remove the oil in the PCR tube, leaving just the emulsion layer (it's okay to have a little bit of oil remaining).

34 Thermocycle the PCR emulsion as follows:

4h

 95 °C 5 min

40 cycles of:

 95 °C 30s

 72 °C 10s

 60 °C 5 min

 72 °C 30s

Final incubation of:

72 °C 10min

12 °C Hold

All ramp times are at 1 °C per second

PCR Cleanup

1m

35 Keep the emulsion on ice to prevent polymerase activity

Add 25 µL EDTA (0.5 M), pH 8.0 Sigma
Aldrich Catalog #AM9260G to the emulsion

Vortex the emulsion to mix

36 Add 25 µL HFE 7500 Perfluorinated Oil Sigma
Aldrich to the emulsion

Add 25 µL Perfluorooctanol Sigma
Aldrich Catalog #370533 to the emulsion

Vortex the emulsion to mix

37 Wait 00:01:00 , then pulse centrifuge to separate the PCR mix from the oil.

1m

Note: If you do not see a clear separation of two clear phases and no emulsion remaining, repeat step 35.

Transfer the top aqueous phase to a new 1mL tube.

38 Add 20 µL 1 Molarity (M) TCEP-HCl Sigma
Aldrich Catalog #TCEP to the tube and
vortex to completely decrosslink any remaining gels.


Note

You should be unable to obtain any "jellyish" substance by centrifugation! If there is any jellyish substance left it is not fully de-crosslinked. **ADD MORE TCEP.**

39 Clean up the PCR reaction using the DNA Clean & Concentrator™-5 Sigma
Aldrich Catalog #D4003
kit.


10m

Elute in  50 µL Elution Buffer.

- 40
- Remove primer dimers and free barcodes using the
- 10m
-  SPRIsselect reagent kit Sigma
Aldrich Catalog #B23317
- with 0.7X volume of beads.


- 41
- Check the resulting library for primer dimers using
- 10m

Equipment	
TapeStation	NAME
Agilent	BRAND
G2991AA	SKU
https://www.agilent.com/en/product/tapestation-automated-electrophoresis/tapestation-instruments/4200-tapestation-system-228263	LINK

with a  High Sensitivity D1000 ScreenTape Sigma
Aldrich Catalog #5067-5584

Other high sensitivity capillary electrophoresis methods will also work.

There should be minimal primer dimers on the trace. Below is an example of an acceptable trace.

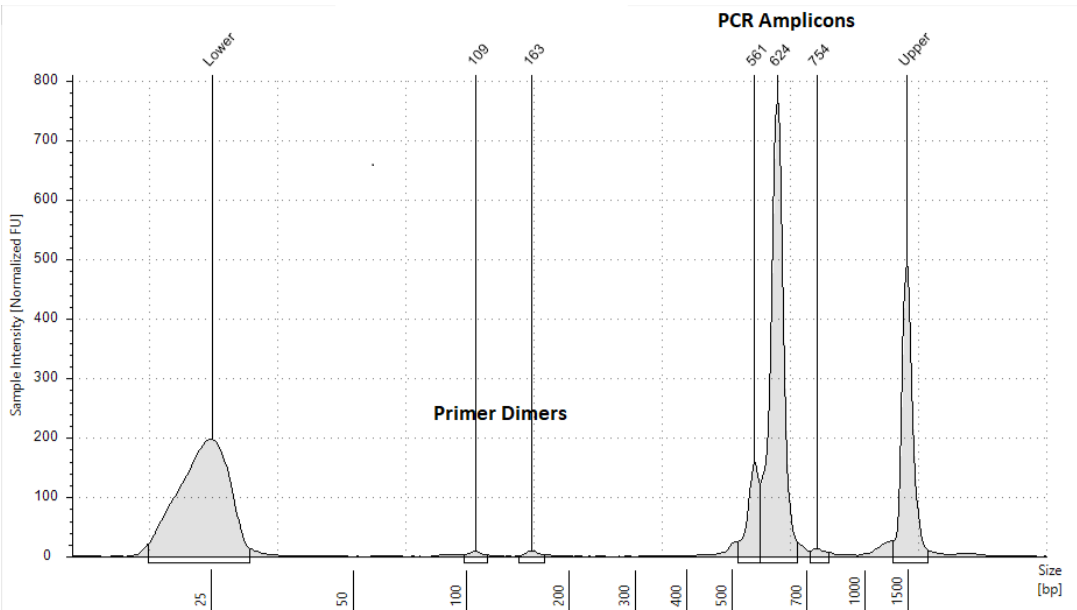

Example of an acceptable TapeStation trace.

42

Quantify the library using a qPCR library quantification kit such as

1h



NEBNext Library Quant Kit for Illumina - 100 rxns Sigma

Aldrich Catalog #E7630S

Note

Note that you must use a PCR based library quantification kit as not all amplicons contain all the adaptors for sequencing and therefore will throw off sequence non-specific forms of quantification!

43

Sequence the library on an Illumina sequencer using Custom Sequencing Primers listed here.



DoTA-seq-Oligo-Sequences.xlsx