




APR 03, 2023

Virus production

 In 2 collections

michela.deleidi¹

¹German Center for Neurodegenerative Diseases (DZNE), Tübingen, 72076 Germany

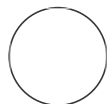

Federico Bertoli

ABSTRACT

Lentiviral production protocol

OPEN  ACCESS

DOI:
dx.doi.org/10.17504/protocols.io.8epv5jpkjl1b/v1

Protocol Citation: michela.deleidi 2023. Virus production.
protocols.io
<https://dx.doi.org/10.17504/protocols.io.8epv5jpkjl1b/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Apr 01, 2023

Last Modified: Apr 03, 2023

PROTOCOL integer ID:
79877

Virus production

- 1 For lentiviral production, HEK cells were transfected with the plasmids expressing the shRNAs together with the lentiviral packaging plasmid psPAX2(RRID:Addgene_12260) and the envelope

plasmid pMD2.G (RRID:Addgene_12259) using TransIT-X2 (Mirus).

- 2 HEK 293 T cells were seeded in a 10-cm dish to reach 80% confluency the day after transfection. The next day, the medium was changed.
- 3 On days 4 and 5 after transfection, the medium was collected and filtered through a 0.45- μ m PVDF membrane.
- 4 To concentrate the virus, the filtered supernatant was centrifuged in a Vivaspin column (Sartorius Stedim Biotech) at 3000 RCF at 4 °C until the volume reached 500–1000 μ l.
- 5 The virus-containing supernatant was collected, and the concentration of p24 particles was determined with Lenti-X GoStix Plus (Takara).
- 6 Cells were infected with equal concentrations of p24 particles for the scramble and respective shRNA.