



OCT 20, 2023

OPEN ACCESS



**DOI:**  
[dx.doi.org/10.17504/protocols.io.eq2lyje1plx9/v1](https://dx.doi.org/10.17504/protocols.io.eq2lyje1plx9/v1)

**Protocol Citation:** Núria Peñuelas 2023. Intestinal transit time. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.eq2lyje1plx9/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
We use this protocol and it's working

**Created:** Oct 19, 2023

**Last Modified:** Oct 20, 2023


## Intestinal transit time

Núria Peñuelas<sup>1</sup>

<sup>1</sup>Vall d'Hebron Research Institute

Vilalab Public

Vilalab Public

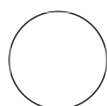

Núria Peñuelas  
Vall d'Hebron Research Institute

### ABSTRACT

Intestinal transit time measurement in mice

- 1 Prepare red dye (Carminered dye, Sigma #C1022) in 0,5% methylcellulose (Sigma #M7027).
- 2 Take the animal and administer the red dye by oral gavage.
- 3 After oral gavage, individualize the animals in separate boxes without bedding.
- 4 Monitor the animal boxes for a total time of 4h.
- 5 Register the time when the animals defecate the first dyed fece (stained red).
- 6 Represent the time until first dyed fece as the transit time of the animals.