



VERSION 1


OCT 17, 2023

🌐 Set Shifting ASAP Operant Behavior_LernerLab V.1

jillian.seiler¹

¹Northwestern University

Lerner Lab

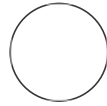

jillian.seiler

ABSTRACT

Protocol for set shifting task used by the Lerner Lab for ASAP work

MATERIALS

Med Associates operant chambers

OPEN  ACCESS



Protocol Citation: jillian.seiler 2023. Set Shifting ASAP Operant Behavior_LernerLab. protocols.io

<https://protocols.io/view/set-shifting-asap-operant-behavior-lernerlab-c3hcyj2w>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Oct 16, 2023

Last Modified: Oct 17, 2023

PROTOCOL integer ID:
89348

Magazine Training

- 1 Food restrict animals overnight prior to beginning magazine training. For the remainder of training animals should be restricted to 85% of their free feeding weight.
- 2 Connect 20 mL syringe with >10 mLs of 1:3 diluted sweetened condensed milk. Add ~15 inches of tubing with tubing adapter to syringe. Load in pump, attach to reward port. Prime pump.
- 3 Place animals in operant box with dummy patch cords attached and allow ~5 minutes to habituate
- 4 Load JN_ASAP_A_MagTraining_30min into MED-PC and start program. The script for the program is attached as a comment. The overview is as follows:
 - 4.1 House light, fan, corner lights, and cue lights above the left and right nosepoke ports are turned on.
 - 4.2 On average, every 60 seconds a noncontingent reward is delivered
 - 4.3 The animal must consume the reward (or make a port entry) in order for the next reward to be delivered
 - 4.4 After 30 minutes, the program ends
- 5 This program occurs once, regardless of performance

- 6 After animal has finished, record weight, return to homecage and clean operant chamber tray with 70% ethanol
- 7 Flush reward port and tubing with ethanol and hang to dry
- 8 Store any remaining sweetened condensed milk in refrigerator

FR1 Both

- 9 Repeat steps 2+3 [⇒ go to step #2](#)
- 10 Load JN_ASAP_B_FR1_BOTH and start program. The script for the program is attached as a comment. The overview is as follows:
 - 10.1 Fan and Near-Infrared Lights are turned on