



AUG 10, 2023

🌐 Measurement of GLP-1 release in cell supernatant from Hutu-80 enteroendocrine cells via ELISA

r.mezabrovschi¹, rachel.bates¹

¹University College London, Aligning Science Across Parkinson's (ASAP) Collaborative Research Network, Chevy Chase, MD, 20815.



Sofia Koletsi

University College London, University of London

DISCLAIMER

This protocol was adapted from Cat. # EGLP-35K to optimise the seeding density of the assay based on the cell line used (Hutu 80 enteroendocrine cells RRID: CVCL_1301), moreover, cell line used and media conditions used to characterise GLP-1 release.

The lowest level of GLP-1 that can be detected by this assay is 2 pM (100ul plasma sample size).

OPEN  ACCESS



DOI:

dx.doi.org/10.17504/protocols.io.j8nlkokm6v5r/v1

Protocol Citation: r.mezabrovschi, rachel.bates 2023. Measurement of GLP-1 release in cell supernatant from Hutu-80 enteroendocrine cells via ELISA . **protocols.io** <https://dx.doi.org/10.17504/protocols.io.j8nlkokm6v5r/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Jul 26, 2023

ABSTRACT

This protocol describes the method of use for the Glucagon-Like Peptide-1 (Active) ELISA Kit
96-Well Plate (Cat. # EGLP-35K) to measure GLP-1 release (pM) from supernatant. The assay should be run in duplicate.

MATERIALS

GLP-1 (Active) ELISA Plate

Coated with anti-GLP-1 Monoclonal Antibody Quantity: 1 plate
Preparation: Ready to use

Adhesive Plate Sealer

Quantity: 1 Sheet
Preparation: Ready to use

10X Wash Buffer Concentrate

10X concentrate of 10 mM PBS Buffer containing Tween 20 and Sodium Azide.
Quantity: 50 mL
Preparation: Dilute 1:10 with deionized water

GLP-1 (7-36) amide ELISA Standards

Last Modified: Aug 10, 2023

PROTOCOL integer ID: 85525

GLP-1 (7-36 amide) in Assay Buffer: 2, 5, 10, 20, 50 and 100 pM Quantity: 1 mL/vial
Preparation: Ready to use

ELISA GLP-1 (Active) Quality Controls 1 and 2

Various peptides including GLP-1 (7-36 amide) in QC Buffer. Quantity: 1 mL/vial
Preparation: Ready to use

GLP-1 (Active) Assay Buffer

0.05M PBS, pH 6.8, containing proprietary protease inhibitors, with Tween 20, 0.08% Sodium Azide and 1% BSA.

Quantity: 25 mL

Preparation: Ready to use

GLP-1 (Active) Detection Conjugate

Anti GLP-1-Alkaline Phosphate Conjugate. Quantity: 21 mL

Preparation: Ready to use

Substrate (Light sensitive, avoid unnecessary exposure to light)

Quantity: 10 mg

Preparation: Hydrate in 1 mL deionized water just before use. Use at 1:200 dilution in substrate diluent (e.g. 100 uL hydrated substrate in 20 mL substrate diluent).

Dilute fresh each time just before use.

Substrate Diluent (Light sensitive, avoid unnecessary exposure to light)

Quantity: 21 mL

Preparation: Ready to use

Stop Solution

Quantity: 6 mL

Preparation: Bring to room temperature before use. Mix thoroughly to ensure no precipitate remains.

Cell Culture Medium

Specific to cell type.

Sterile 24-well plates

Quantity: 2

PBS (1X)

Preparation: Ready to use

Reagents for you experimental media conditions

Varies between experiments.

DPP-1V Inhibitor



Quantity: 10ml (Cat. DPP4-010)

Preparation: Ready to use - to be purchased separately. (Store -20 degrees).

BEFORE START INSTRUCTIONS






All reagents should be warmed to room temperature before proceeding.











Day 1


- 1 Plate cells in sterile 24-well plates at 1.5×10^5 per well and leave  Overnight in a  37 °C incubator with stable CO2 conditions overnight.


Day 2


5m

- 2 Check cells are healthy under a microscope.
- 3 In a laminar flow hood aspirate medium and wash with PBS 3 times, adding  500 µL each time.
- 4 Once the final PBS wash is complete add in  500 µL of the test media made up for your desired experiment (e.g., high glucose media or 2-Deoxy-D-glucose) and place back in the incubator for  02:00:00 . 2h
- 5 Make up the **Wash Buffer** using the **10X Wash Buffer Concentrate** and dilute 1:10 with deionised water.
- 6 Once the plate has finished incubating, move over the  500 µL supernatant into a new sterile 24-well plate and add DPP-IV inhibitor immediately (1:100).to be carried out in a laminar flood hood. The cells may be stored at  -80 °C for further analysis at a later date.

- 7 In each well of the ELISA plate, add  300 μL of diluted **Wash Buffer** (or in two intervals of  150 μL). Incubate at  Room temperature for  00:05:00 . Decant excess buffer and blot with absorbent towels.
- 8 Add  200 μL **Assay Buffer** to NSB (non-specific binding) wells A1, A2.
- 9 Add  100 μL **Assay Buffer** to the remaining wells you wish to load your samples in.
- 10 Add  100 μL standards in ascending order to wells - standards come preprepared.
- 11 Then load  100 μL of QC1 and QC2 to the plate in separate wells.
- 12 Add  100 μL of your samples from the 24-well plate in the remaining wells (already containing  100 μL of the Assay Buffer). Plate layout should resemble the below.

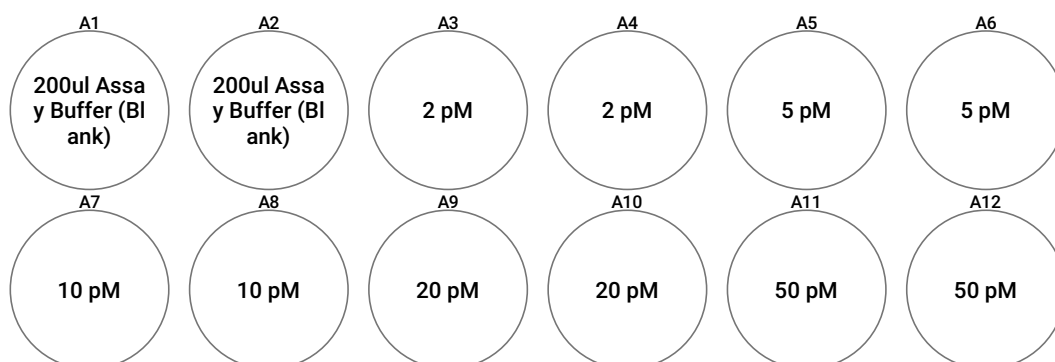

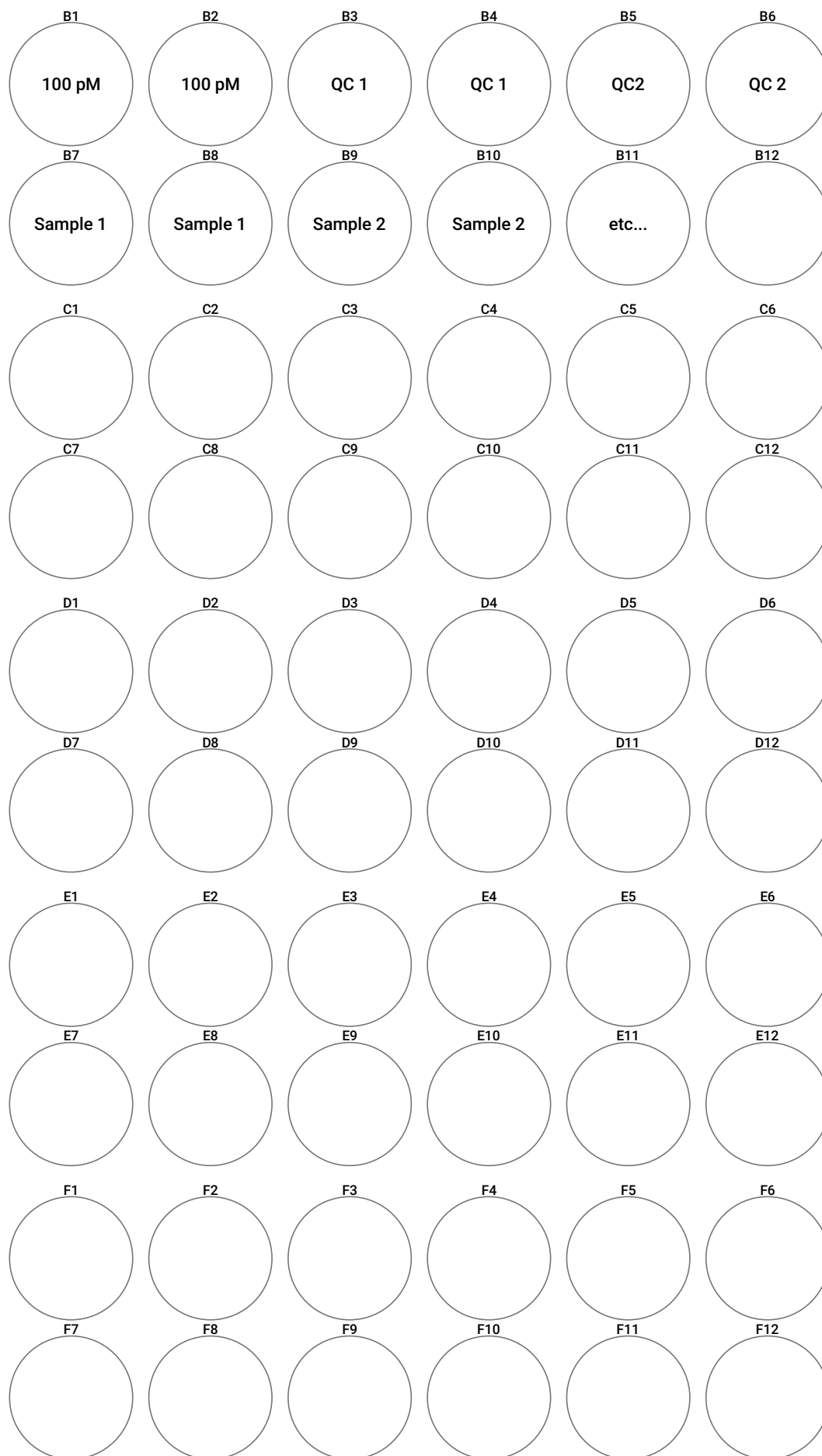

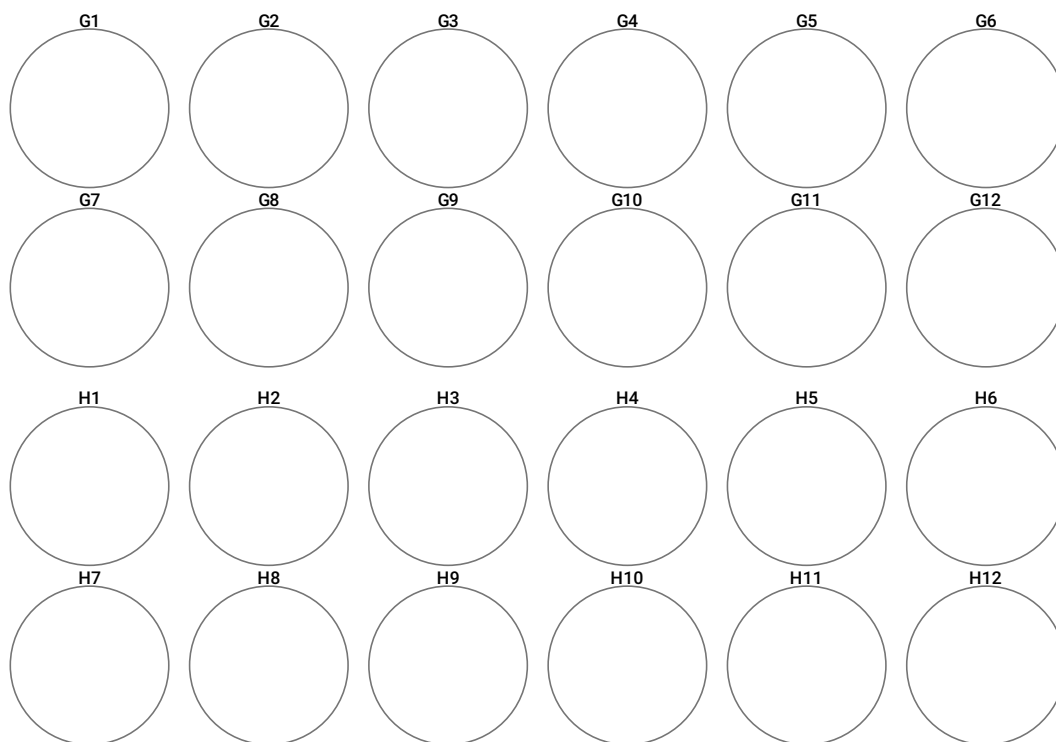

13 For good mixing, lightly agitate the plate.

14 Cover the plate with plate sealer. Incubate Overnight (20 to 24 hours) at 4 °C .

Day 3

2h 35m

15 *All reagents should be warmed to Room temperature before proceeding.*

Decant liquid from plate and tap out excess fluid on absorbent towels.

16 Wash the plate 5 times with 300 µL pre-diluted **Wash Buffer** per well with 00:05:00 incubation at Room temperature in **Wash Buffer** with the fourth wash. Tap out excess buffer on absorbent towels after the fifth wash.

5m

- 17 Immediately add  200 μL **Detection Conjugate** (is ready to use) in each well. Incubate  02:00:00 at  Room temperature then decant. 2h
- 18 Whilst the plate is incubating dilute the **Substrate**.
- 18.1 Hydrate in  1 mL deionized water just before use.
- 18.2 Use at 1:200 dilution in **Substrate Diluent** (e.g.  100 μL hydrated substrate in  20 mL substrate diluent). Dilute fresh each time just before use.
- 19 Wash the wells 3 times with  300 μL diluted **Wash Buffer**. Tap out excess buffer on absorbent towels.
- 20 Add  200 μL diluted **Substrate** in each well.
- 21 Measure fluorescence on a plate reader at an excitation/emission wavelength of 360/460 every  00:05:00 for a minimum of  00:20:00. 25m
- 22 If sufficient fluorochrome has been generated, add  50 μL **Stop Solution** (mix thoroughly to ensure no precipitate remains) to each well in the same order as the **Substrate** was added. Incubate  00:05:00 at  Room temperature in the dark to arrest phosphatase activity. 5m
- 23 Read plate on a fluorescence plate reader with an excitation/emission wavelength of 360/460.

