



VERSION 2

MAY 04, 2023

OPEN ACCESS

DOI:

dx.doi.org/10.17504/protocols.io.3byl4k1jzvo5/v2

Protocol Citation: Bjørn Peare Bartholdy, a.g.henry 2023. Artificial saliva.

protocols.io

<https://dx.doi.org/10.17504/protocols.io.3byl4k1jzvo5/v2> Version created by Bjørn Peare Bartholdy

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Dec 10, 2021

Last Modified: May 04, 2023

PROTOCOL integer ID:
55830


Artificial saliva V.2

In 1 collection

Bjørn Peare Bartholdy¹, a.g.henry¹

¹Leiden University

BYOC

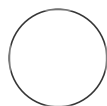

Bjørn Peare Bartholdy

Delft University of Technology, Leiden University

ABSTRACT

Creating an artificial saliva solution for oral biofilm growth, specifically dental calculus.

This is a modified version of Sissons et al. 1991.

GUIDELINES

This protocol explains how to mix  1000 mL of artificial saliva.

MATERIALS

Chemicals

Mucin from porcine stomach (Type III)

Trypticase peptone

Proteose peptone

Yeast Extract

KCl

NaCl

CaCl₂

Na₂HPO₄

NaHCO₃

Hemin

Menadione

Urea

L-Arginine

Equipment

1000 mL beaker

100 mL graduated cylinder

2 x 1000 mL bottle

Magnetic stirrer (with heating element)




Autoclave

pH measure


BEFORE START INSTRUCTIONS

Mix the solution under a fumehood. It can smell pretty bad.

Also have the incubator under a fumehood, if possible. I had no such luxury...



1 Add  300 mL distilled (or deionized) dH₂O to a  1000 mL beaker, with stirring and heat  60 °C .

2 Add:

-  2.5 g




Mucin from porcine stomach (Type III) Merck MilliporeSigma (Sigma-Aldrich) Catalog #M1778

-  5 g  Trypticase™ Peptone Thermo Fisher Catalog #211921

- 10 g Oxoid™ Proteose Peptone Thermo Fisher Catalog #LP0085B
- 5 g Bacto Yeast Extract Becton-Dickinson

Let the reagents completely dissolve before continuing to the next step

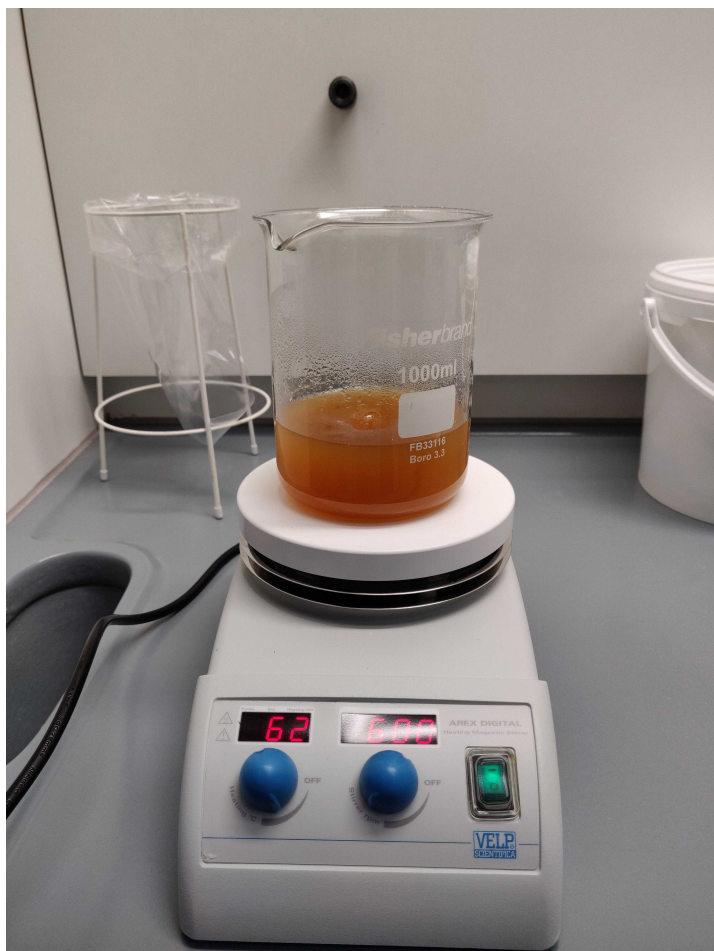

3 Add:

- 2.5 g KCl Contributed by users
- 0.35 g NaCl Contributed by users
- 0.2 g CaCl₂ Contributed by users
- 0.74 g Sodium phosphate dibasic Merck MilliporeSigma (Sigma-Aldrich) Catalog #7558-79-4
- 0.54 g NaHCO₃ Contributed by users

- 2.5 mg Hemin Contributed by users

4 Add the remaining 700 mL distilled (or deionized) H₂O



5 Adjust to pH 7 with NaOH Contributed by users and stirring

6 Transfer to two 1000 ml bottles, so around half of each bottle is filled.

15m







Autoclave at 121 °C , 1 Bar for 00:15:00 minutes

Safety information

Do NOT screw bottle caps on tightly before autoclaving.

Loosely screw the caps on the bottles or cover the tops with foil.

7 Once the solution has cooled, combine solutions to a single 1000ml bottle, and add:


-  1 mg  Menadione Contributed by users
-  0.3 g  Urea Contributed by users
-  0.17 g  L-Arginine Contributed by users Catalog #A5006

Safety information

These should be filter-sterilised as there are no more decontamination steps for the artificial saliva.

Use clean gloves and take precautions not to contaminate the artificial saliva with bacteria or starch.

8 Store in fridge at ca.  4 °C

Occasionally test the pH of the solution to ensure it stays around  7