



JUN 26, 2023

OPEN ACCESS

**DOI:**  
[dx.doi.org/10.17504/protocols.io.j8nlkobdv5r/v1](https://dx.doi.org/10.17504/protocols.io.j8nlkobdv5r/v1)

**Protocol Citation:** Justin T Savage, Juan Ramirez, Dolores Irala, Cagla Eroglu 2023. Glia-Free Cortical Neuronal Feeding Schedule - Synapse Formation.  
**protocols.io**  
<https://dx.doi.org/10.17504/protocols.io.j8nlkobdv5r/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working

**Created:** Jun 09, 2023

**Last Modified:** Jun 26, 2023

**PROTOCOL integer ID:**  
 83097

## Glia-Free Cortical Neuronal Feeding Schedule - Synapse Formation

In 1 collection


Justin T Savage<sup>1</sup>, Juan Ramirez<sup>1</sup>, Dolores Irala<sup>2</sup>, Cagla Eroglu<sup>1,2,3</sup>

<sup>1</sup>Department of Neurobiology, Duke University Medical Center, Durham, NC 27710;

<sup>2</sup>Department of Cell Biology, Duke University Medical Center, Durham, NC 27710;

<sup>3</sup>Howard Hughes Medical Institute, Duke University Medical Center, Durham, NC 27710

Eroglu\_Lab

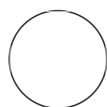

Justin T Savage

### ABSTRACT

Protocol for isolating rat cortical astrocytes and producing astrocyte-conditioned media for synaptogenesis assays.

### MATERIALS

**Penicillin-Streptomycin** (10,000 U/mL) (Gibco, 15140148)

**Sodium Pyruvate** (GIBCO 11360)

**Neurobasal plus** (Gibco, A35029-01)

**B27 plus** (Gibco, A35828-01)

**Forskolin (5 mM; 1000X)**

Add 12ml DMSO to 50 mg of Forskolin (Sigma, F6886).  
 Make 20 ul aliquots; store -20°C.

**BDNF** – (Peprotech, 450-02). Master stock should be prepared at 1mg/ml in 0.2%BSA DPBS and is stored at -80°C as 100ul aliquots.

*Preparation of 1000x Working Stocks*

Thaw on ice one master stock aliquot (100ul) of BDNF. Cool down on ice 2ml of filtered 0.2% BSA (Sigma A-8806) in DPBS solution. Add BDNF master stock aliquot to the BSA solution. Mix well but gently to avoiding foaming. Aliquot in 0.5ml tubes

20ul each, concentration 50ug/ml (1000X). Flash freeze the aliquots in liquid nitrogen and store at -80C. Working concentration: 50ng/ml, final concentration on neurons 25ng/ml.

**CNTF** – (Peprotech, 450-13). Master stock should be prepared at 1mg/ml in 0.2%BSA DPBS and is stored at -80°C as 50ul aliquots.

*Preparation of 1000x Working Stocks*

Thaw on ice one master stock aliquot (100ul) of CNTF. Cool down on ice 2.5ml of filtered 0.2% BSA (Sigma A-8806) in DPBS solution. Add CNTF master stock aliquot to the BSA solution. Mix well but gently to avoiding foaming. Aliquot in 0.5ml tubes 20ul each, concentration 20ug/ml (1000X). Flash freeze the aliquots in liquid nitrogen and store at -80C. Working concentration: 20ng/ml, final concentration on neurons 10ng/ml.

**Antibody Blocking Buffer** - 150 mM NaCl, 50 mM Tris-Base, 1% BSA, 100 mM L-lysine.

In a 1L baker add 4.383g NaCl, 3.025g Tris-Base (Fisher, Cat. No: BP152-5), 5g BSA (Sigma, Cat. No:A2153) and 9.125g L-Lysine monohydrochloride (Sigma, Cat. No: L-1137). Add 350-400ml of ddH<sub>2</sub>O and mix well. Once mixed, adjust pH with HCl to 7.5. Finally, add 5ml of Sodium Azide (NaN<sub>3</sub>). Add ddH<sub>2</sub>O to final volume 500ml. Filter through 0.22µm filter and store 4C.

**Mounting media:**

Final Composition: 20mM Tris pH8.0, 90% Glycerol, 0.5% N-propyl gallate.

In a 50ml conical tube add: 9ml Glycerol, 1ml 200mM Tris pH8.0 and 50mg N-propyl gallate. Mix on nutator in 37C incubator overnight. Store at 4C.

**PFA 16%** (Electron Microscopy Sciences, 15710)

**Normal Goat Serum** (NGS) - (ThermoFisher 01-6201)

**Triton X-100** (Roche, 11332481001)

**Anti-Basson antibody**(Enzo/Assay Designs, SAP7F07/VAM-PS003F)

**Anti-Gephyrin antibody** (Synaptic System, 147 011)

**Anti-Homer1 antibody** (Synaptic System, 160 011)

**Anti-VGAT antibody** (Synaptic Systems, 131004)

**Anti-Vglut1 antibody** (Millipore, AB5905)

**Alexa Fluor 488 goat anti-Mouse IgG2a (H+L)** (Invitrogen, A-21131)

**Alexa Fluor 568 goat anti-Mouse IgG1 (H+L)** (Invitrogen, A-21124)

**Alexa Fluor 647 goat anti-Guinea pig IgG (H+L)** (Invitrogen, A21450)

**DAPI** (Invitrogen, D1306)

## Neuronal growth media (NGM) for feeding

**1** This recipe makes 20ml of neuronal media. Make fresh per use

**2** 1. To a 50ml tube, add the following media components:

Reagent	Volume
Neurobasal plus	19ml
Pen/strep (100x)	200µl
Sodium Pyruvate (100x)	200µl
B27 plus (50x)	400µl

**3** Sterile filter these components through a syringe filter. Place media in the incubator with cap unscrewed. Allow the media to warm and equilibrate for at least 45 minutes. (Media will have an orange color and bubbles)

**4** Add growth factors right before the time to use the media.

BDNF	20µl
CNTF	20µl
Forskolin	20µl

## DIV 2 - AraC feeding

**5**

#### Note

Neurons are extremely sensitive to environmental changes, so only half of the media is replaced.

- 6 Check neurons for normal morphology under a microscope.
- 7 Prepare neuronal media with Neurobasal plus + B27 plus according to recipe fresh + add AraC (1:10000). Let equilibrate for at least 45 min in the incubator with the cap unscrewed.
- 8 Remove 230µl of old neuronal media.
- 9 Add 250µl of fresh media.
- 10 Store in the incubator

### DIV 3 - remove AraC from media

- 11 Check neurons for normal morphology under a microscope.
- 12 Prepare neuronal media with Neurobasal plus + B27 plus according to recipe fresh. Let equilibrate for at least 45 min in the incubator with the cap unscrewed.

**13** Remove 230µl of old neuronal media.

**14** Add 250µl of fresh media.

**15** Store in the incubator.

## **DIV 6 - Media feeding**

**16** Check neurons for normal morphology under microscope.

**17** Prepare neuronal media with Neurobasal plus + B27 according to recipe fresh. Let equilibrate for at least 45 min in the incubator with the cap unscrewed.

**18** Remove 230µl of old neuronal media.

**19** Add 250µl of fresh media.

**20** Store in the incubator.

## DIV 8 - ACM feeding

- 21 Check neurons for normal morphology under a microscope.
- 22 Prepare neuronal media with Neurobasal plus + B27 according to recipe fresh. Let equilibrate for at least 45 min in the incubator with the cap unscrewed.
- 23 Place the NGM on ice to cool down.
- 24 Thaw the ACM on ice.
- 25 Add AGM to the NGM for a final concentration of 50ug/ml for excitatory synapses and 100ug/ml for inhibitory synapses.
- 26 Equilibrate and warm up the NGM with ACM in the incubator (37C, 10% CO2) before adding it to the cells.
- 27 Remove 230µl of old neuronal media.
- 28 Add 250µl of fresh media.

**29** Store in the incubator.

## **DIV 11 - ACM feeding**

**30** Check neurons for normal morphology under a microscope.

**31** Prepare neuronal media with Neurobasal plus + B27 according to recipe fresh. Let equilibrate for at least 45 min in the incubator with the cap unscrewed.

**32** Place the NGM on ice to cool down.

**33** Thaw the ACM on ice.

**34** Add AGM to the NGM for a final concentration of 50ug/ml for excitatory synapses and 100ug/ml for inhibitory synapses.

**35** Equilibrate and warm up the NGM with ACM in the incubator (37C, 10% CO<sub>2</sub>) before adding it to the cells.

**36** Remove 230µl of old neuronal media.

**37** Add 250µl of fresh media.

**38** Store in the incubator.

## **DIV 13 - Neuron culture fixation and staining**

**39** The following solutions should be prepared:

a) 4% PFA in PBS – In a 50ml conical tube add one 16% PFA ampule (10ml), 26ml of ddH<sub>2</sub>O and 4ml of 10X PBS. Mix and place at 37°C to warm up. For storing keep away from the light at 4°C.

b) Blocking and permeabilization solution (0.2% Triton) - This recipe makes 5ml of blocking solution, sufficient for 1 24-well plate.

- 2.4ml of antibody-blocking buffer
- 2.5ml of Normal Goat Serum (NGS)
- 100µl of 10% Triton.

c) Primary antibodies solution: This recipe makes for 5ml, sufficient for 1 24-well plate.

- 4500µl Antibody blocking buffer
- 500µl Normal Goat Serum (NGS)
- 10µl of anti-Bassoon antibody (presynaptic marker)
- 10µl of anti-Homer1 antibody (for excitatory synapses) or anti-Gephyrin antibody (for inhibitory synapses)
- 5µl of anti-Vglu1 antibody (for excitatory synapses) or anti-VGAT antibody (for inhibitory synapses)

**40** In a chemical fume hood, aspirate media from cells and immediately add 500µl of warm 4% PFA to each well.

**41** Fix the neurons for 7 minutes at room temperature.



- 42 Wash 3 times with PBS.
- 43 Add 200ul of blocking solution to each well and block for 30min at room temperature.
- 44 Remove the blocking solution and add 200ul of primary antibody solution to each well. Incubate at 4°C overnight.

## DIV 14 - Secondary antibody staining and mounting

- 45 a) Secondary antibodies solution: This recipe makes for 5ml, sufficient for 1 24-well plate.
- 4500ul Antibody blocking buffer
  - 500ul Normal Goat Serum (NGS)
  - 10ul of secondary antibody Alexa Fluor anti-Mouse IgG2a 488
  - 10ul of secondary antibody Alexa Fluor anti-Rabbit IgG1 564
  - 10ul of secondary antibody Alexa Fluor anti-Guinea Pig IgG 647
  - 0.5ul DAPI (1/10000)
- 46 Remove the primary antibody and wash 3 times with PBS
- 47 Add 200ul of secondary antibody solution to each well and incubate at room temperature for 2h, keeping it protected from light.
- 48 Wash 3 times with PBS.
- 49 Mount coverslips onto slides using one small drop of mounting media for each coverslip.

**50** Seal with nail polish and dry at room temp for at least 30min before storing at 4°C.