



OCT 18, 2023


## Quantification of p62-positive inclusions

Núria  
Peñuelas<sup>1</sup>

<sup>1</sup>Vall d'Hebron Research Institute

Vilalab Public

Vilalab Public

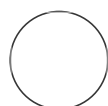

Núria Peñuelas  
Vall d'Hebron Research Institute

OPEN ACCESS



### ABSTRACT

Quantification of p62-positive inclusions in the rodent brain

**DOI:**  
[dx.doi.org/10.17504/protocols.io.8epv5x33ng1b/v1](https://dx.doi.org/10.17504/protocols.io.8epv5x33ng1b/v1)

**Protocol Citation:** Núria Peñuelas 2023.  
Quantification of p62-positive inclusions. **protocols.io**  
<https://dx.doi.org/10.17504/protocols.io.8epv5x33ng1b/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
We use this protocol and it's working

**Created:** Oct 18, 2023

**Last Modified:** Oct 18, 2023

- 1** Scan the 5um sections immunostained with guineapig anti-p62 (1:500, Progen #GP62-C), mouse anti- $\alpha$ -synuclein (1:500, BD Biosciences #610786) and TH (1:1000, Calbiochem #657012) using ZEISS LSM 980 with Airyscan 2.
  
- 2** Visualize the images with ZEN 2011 software (Zeiss, Germany) and determine one representative coronal section of the brain region of interest.
  
- 3** Determine the total number of p62-immunopositive Marinesco bodies (nuclear inclusions) and Lewy body-like inclusions (cytoplasmic inclusions) falling into each category (Syn-positive or Syn-negative) in that representative section.