



Jul 12, 2022


Primary astrocyte culture

Xiqun Chen¹, Qing Ye²

¹Department of Neurology Massachusetts General Hospital Harvard Medical School Charlestown MA 02129 USA, Aligning Science Across Parkinson's (ASAP) Collaborative Research Network Chevy Chase MD 20815 USA;

²Department of Neurology Massachusetts General Hospital Harvard Medical School Charlestown MA 02129 USA, Department of Neurology Longhua Hospital Shanghai University of Traditional Chinese Medicine Shanghai 200032 China

1 Works for me

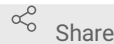

Share

dx.doi.org/10.17504/protocols.io.b4yfqxtn

Pranay Srivastava

ABSTRACT

Primary astrocytes were obtained from C57BL/6 mice embryonic day 17. The dissected cortical tissue was digested, triturated, and centrifuged. The cell pellet was resuspended in high-glucose DMEM/F12 supplemented with 10% FBS. Cells were seeded in poly-L-ornithine-coated Petri plates. Nonadherent cells were removed after 5-7 days.

DOI

dx.doi.org/10.17504/protocols.io.b4yfqxtn

PROTOCOL CITATION

Xiqun Chen, Qing Ye 2022. Primary astrocyte culture. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.b4yfqxtn>



KEYWORDS

ASAPCRN

LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Feb 11, 2022

LAST MODIFIED

Jul 12, 2022

PROTOCOL INTEGER ID

58087

- 1 For primary astrocyte culture - Use C57BL/6J mice at embryonic day 17
- 2 Anesthetized pregnant mice (1% sodium pentobarbital, 80mg/kg), dissect their embryos and collect the cortex.
(Separate and remove the soft membrane and blood vessels, rinse the cerebral cortex in PBS, and use the ophthalmic scissor to cut pieces of the cortex)
- 3 Collect the cortices in PBS in a 50 ml tube on ice
(The 50 ml tube contains 30 ml of PBS) ⚡ On ice
- 4 Transfer the cortices to 15 ml tubes containing 1.5 ml trypsin–EDTA (0.25%) and incubate it at ^{15m}
⚡ 37 °C for ⌚ 00:15:00 Dissociate the cortices by triturating with a 10 mL serological pipette 10 – 15 times
- 5 Centrifuge the dissociated cortices (⚙ 1500 rpm , ⌚ 00:05:00) and resuspend the pellet in ^{5m} 10ml




⚙ DMEM, high glucose Thermo Fisher

Scientific Catalog #11965092

medium
 supplemented with 10%

⚙ Fetal Bovine Serum, certified, heat inactivated, United States Thermo

Fisher Catalog #10082147
- 6 Triturate the cell suspension 10 times with a 1ml pipette
- 7 Count the cells and plate them in a density of 50,000 cells/cm² into PLL coated T-75 flask.
- 8 After 3 h, change the culture medium once and check the growth and survival of cells under the microscope.
- 9 Put the culture flask in a ⚡ 37 °C -cell incubator for ⌚ 24:00:00 1d
- 10 Change the media after ⌚ 24:00:00 and then every ⌚ 72:00:00 . A confluent layer of the 4d
astrocytes grows at the bottom in 5-7 days followed by a top layer of microglia.

- 11 After 5-7 days, put the flask on the shaker for 2 h at 200 rpm at 37°C (incubator). Collect and discard the supernatant (consisting of microglia and some oligodendrocytes).
- 12 Wash the adherent cells with PBS and add 0.25% trypsin at 37°C to detach the cells.
Immediately add  10 mL DMEM medium containing 10% FBS
- 13 Centrifuge the cell suspension at  1500 rpm min for  00:05:00 5m

Resuspend the cells in a T75 flask or 24 well plate with high glucose DMEM supplemented with 10% FBS.
- 14 Identification of astrocytes by immunofluorescence staining: Use GFAP antibody (astrocytic marker) to identify the population of astrocytes
- 15 **Transduction with BRAF** (Optional)
Astrocytes were transduced with BRAFV^{600E}, or BRAF^{WT}, or vector lentivirus plus 8 µg/ml polybrene for 24 h.
- 16 After transduction, the cells were cultured for 120 h in high-glucose DMEM/F12 (400ul/well) with or without FBS, and the medium as well as the cells were collected for subsequent experiments.