




Aug 03, 2020

SOLUTION- 10 - TEFF/TREG isolation buffer

Marco Cosentino¹, Elisa Storelli¹, Alessandra Luini¹, Massimiliano LM Legnaro¹, Emanuela Rasini¹, Marco Ferrari¹, Franca Marino¹

¹Center for Research in Medical Pharmacology, University of Insubria (Varese, Italy)

1 Works for me dx.doi.org/10.17504/protocols.io.bjadkia6

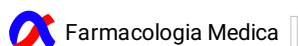

DOI

dx.doi.org/10.17504/protocols.io.bjadkia6

DOCUMENT CITATION

Marco Cosentino, Elisa Storelli, Alessandra Luini, Massimiliano LM Legnaro, Emanuela Rasini, Marco Ferrari, Franca Marino 2020. SOLUTION- 10 - TEFF/TREG isolation buffer. **protocols.io**
dx.doi.org/10.17504/protocols.io.bjadkia6

LICENSE

This is an open access document distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Aug 03, 2020

LAST MODIFIED

Aug 03, 2020

DOCUMENT INTEGER ID

39973

This buffer is composed of PBS, 0.5% BSA and 2mM EDTA.

It must be **cold (2-8°C) and sterile filtered before use.**

BSA 1.250 g,

EDTA 0.2235 g

Dissolve the following powders in **240 mL PBS 1X**, adjust the pH to 7.2 with HCl or NaOH, and then add PBS 1X to 250 mL.

Storage:  **4 °C** (Fridge - 2 Room PS03)

BSA: catalog number A2153 , Sigma, Italy

EDTA: catalog number ED2SS, Sigma, Italy