



Version 2 ▼

Jul 07, 2022

Extraction of bacterial DNA using ChargeSwitch gDNA Mini Bacteria kit on KingFisher Flex robot V.2

Leyi Wang¹, Carol Maddox¹, Melanie Prarat², Yan Zhang², Lifang Yan³, Akhilesh Ramachandran⁴, Sai Sankara Narayanan⁴, Girish Patil⁴


¹University of Illinois Veterinary Diagnostic Laboratory;

²Ohio Department of Agriculture Animal Disease Diagnostic Laboratory;

³Mississippi Veterinary Research and Diagnostic Lab;

⁴Oklahoma Animal Disease Diagnostic Laboratory

1 Works for me

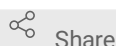

Share

dx.doi.org/10.17504/protocols.io.6qpvrdezgmk/v2

Vet LIRN



Sarah Nemser

DISCLAIMER

Reference to any commercial materials, equipment, or process does not in any way constitute approval, endorsement, or recommendation by the Food and Drug Administration.

ABSTRACT

This procedure is used to extract genome DNA of bacteria isolates using the ChargeSwitch gDNA Mini Bacteria kit on the KingFisher Flex Robot. The process consists of Sample processing, Prepare plates for Robot, and Operate on KingFisher Flex machine.

DOI

dx.doi.org/10.17504/protocols.io.6qpvrdezgmk/v2

PROTOCOL CITATION

Leyi Wang, Carol Maddox, Melanie Prarat, Yan Zhang, Lifang Yan, Akhilesh Ramachandran, Sai Sankara Narayanan, Girish Patil 2022. Extraction of bacterial DNA using ChargeSwitch gDNA Mini Bacteria kit on KingFisher Flex robot.

protocols.io

<https://dx.doi.org/10.17504/protocols.io.6qpvrdezgmk/v2>

Version created by Sarah Nemser

FUNDERS ACKNOWLEDGEMENT



Vet-LIRN

Grant ID: 1U18FD06983-01

KEYWORDS

robot, KingFisher, ChargeSwitch, DNA Extraction

LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Jul 07, 2022

LAST MODIFIED

Jul 07, 2022

PROTOCOL INTEGER ID

66209

GUIDELINES

In the whole procedure, need wear PPE

MATERIALS TEXT

Bacterial culture isolates

ChargeSwitch™ gDNA Mini Bacteria Kit (CS11301)

SAFETY WARNINGS

Need work in the biosafety level 2 hood to scoop the bacteria in the buffer tube

DISCLAIMER:

Reference to any commercial materials, equipment, or process does not in any way constitute approval, endorsement, or recommendation by the Food and Drug Administration.

BEFORE STARTING

At the first time use, add RNase A to Resuspension Buffer (R4).

1 Sample processing

1.1

Bacterial Broth Culture

If using broth, pick single colony and grow in 5 ml Trypticase Soy Broth (TSB), overnight grow at 37 degree without shaking. Use 0.4 ml of the bacterial culture and centrifuge for 4min at 8000 rpm, discard supernatant and resuspend the cell pellet in 100 μ L of Resuspension Buffer (R4) containing RNase A, then add 5 μ L of lysozyme solution (50 mg/ML). Mix by vortexing. After adding RNase A to Resuspension Buffer (R4), store the buffer at 4 degrees Celsius.

1.2

Bacterial Colonies in Plate

If using bacterial colonies, scoop 3 colonies into 100 μ L of Resuspension Buffer (R4) containing RNase A, and then add 5 μ L of lysozyme solution (50 mg/ML). Mix by vortexing.

2 Prepare plates for Robot

2.1

Load all resuspension solution prepared in 1.1. or 1.2 into 2.4 mL 96 deep well Sample plate.

2.2

Prepare Elution Plate by adding 100 μ L of Elution buffer to each well of the 0.5ml 96-well plate.

2.3

Prepare 3 Wash Plates (Wash Plate 1, Wash Plate 2 and Wash Plate 3) by adding 650 μ L Wash Buffer (W12) to each well of 2.4 mL 96 deep well plates.

3 Operate on KingFisher Flex machine

3.1

Install script of ChargeSwitch gDNA Mini Bacteria kit on KingFisher Flex by requesting the script through ThermoFisher and installing it using BindIt software

- 3.2 Select the right script/method and then click start.
- 3.3 Load 0.5ml 96-well plate containing the Tip comb.
- 3.4 Load Elution plate.
- 3.5 Load Wash Plate 1, Wash Plate 2, and Wash Plate 3 one by one.
- 3.6 Load Sample plate. The instrument will automatically start processing the sample plate.
- (While the instrument is processing the sample plate, prepare a premix of 500 ul Lysis Buffer (L14) and 10 ul of Proteinase K solution).
- 3.7 6) After 10 mins of sample processing, there is the **first pause**. Once the Instrument pauses, pick up the sample plate from the instrument and add 500 ul of the premix prepared in step 5 (500 ul Lysis Buffer/Proteinase K mixture) to each sample well of Sample plate.
- 3.8 Then load Sample plate back into the instrument and click Start.
- 3.9 After 10 mins, there is a **second pause**, take the sample plate off the instrument and add 40 ul of Charge Switch Magnetic Beads and 300 ul Binding Buffer (B8) to each sample well, and then load the Sample plate back into the instrument and click Start. The machine will do the remaining steps until it is done.