



FEB 16, 2024

🌐 Immunofluorescence staining, vibratome sections

eduard.bentea¹, María Sanchiz Calvo¹, Veerle Baekelandt¹

¹Katholieke Universiteit Leuven



eduard.bentea

ABSTRACT

Protocol for performing immunofluorescence staining on free-floating vibratome cut brain sections from rats or mice.

MATERIALS

Antigen retrieval solution (10 mM citrate buffer pH 6.0)

Prepare stock solutions:

Solution A : 0.1 M citric acid solution

Solution B : 0.1 M sodium citrate solution

Prepare working solution:

Add 9 mL of Stock solution A and 41 mL of Stock solution B to 400 mL of AD

Adjust pH to 6.0

Fill until 500 mL with AD

OPEN ACCESS



DOI:

dx.doi.org/10.17504/protocols.io.3byl4qoxrvo5/v1

Protocol Citation: eduard.bentea, María Sanchiz Calvo, Veerle Baekelandt 2024.

Immunofluorescence staining, vibratome sections. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.3byl4qoxrvo5/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

We use this protocol and it's working

Created: Feb 16, 2024

Last Modified: Feb 16, 2024

PROTOCOL integer ID: 95341




Keywords: ASAPCRN



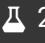

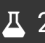
Funders Acknowledgement:

ASAP (Aligning Science Across Parkinson's)

Day 1


2h


- 1 Transfer sections to be stained to a new 24-well plate filled with 1X PBS. All incubations are with  500 μ L solution per well unless stated otherwise.
- 2 1X PBS rinse.
- 3 Antigen retrieval step (using oven).
 - 3.1 Incubate sections in antigen retrieval solution (10 mM citrate buffer pH 6.0 ; see recipe in Materials) at 80°C for  00:30:00 . 30m
 - 3.2 Place the well plate in the same buffer on ice for  00:20:00 . 20m
- 4 1X PBS rinse.

- 5 Wash sections in 1X PBS for  00:05:00 at room temperature on wobbler. 5m
- 6 Wash sections in 1X PBS for  00:05:00 at room temperature on wobbler. 5m
- 7 Incubate sections with  250 μ L blocking solution (PBS-T + 10% donkey serum) per well for  01:00:00 at room temperature on wobbler. (PBS-T = PBS + 0.1% Tergitol) 1h
- 8 Dilute primary antibodies at the required concentration in PBS-T + 10% donkey serum.
- 9 Incubate sections with  250 μ L primary antibodies overnight at room temperature.



Day 2

2h


- 10 1X PBS-T rinse.
- 11 Wash sections in 1X PBS-T for  00:05:00 at room temperature on wobbler. 5m


- 12 Wash sections in 1X PBS-T for  00:05:00 at room temperature on wobbler. 5m

- 13 Dilute secondary antibodies at the required concentration in PBS-T.

- 14 Incubate sections with  250 µL secondary antibodies for  02:00:00 at room temperature (in the dark). 2h

- 15 1X PBS rinse.

- 16 Wash sections in 1X PBS for  00:05:00 at room temperature on wobbler. 5m

- 17 Wash sections in 1X PBS for  00:05:00 at room temperature on wobbler. 5m

- 18 Briefly rinse sections in 1/2 PBS + 1/2 AD and allow to dry.

- 19 Mount with Mowiol.

