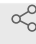




Oct 17, 2022

🌐 Drug feeding protocol

Mel Feany¹¹Brigham and Women's Hospital and Harvard Medical School1 *Works for me* Sharedx.doi.org/10.17504/protocols.io.5jyl8j9n9g2w/v1**Daniel's workspace****Daniel El Kodsi**

ABSTRACT

This protocol covers the steps for drug feeding in fly models.

DOI

dx.doi.org/10.17504/protocols.io.5jyl8j9n9g2w/v1

PROTOCOL CITATION

Mel Feany 2022. Drug feeding protocol. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.5jyl8j9n9g2w/v1>



LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Oct 17, 2022

LAST MODIFIED


Oct 17, 2022


PROTOCOL INTEGER ID

71452

- 1 Collect flies for feeding:
 - Clear bottles the day before collecting
 - Select appropriate genotype for feeding

- 2 Prepare drugs to be fed:

- Thaw drugs
- Add  3 mL water to each vial
- Add appropriate number of ul drug to water, mix by swirling

- 3 Make fly food in vials that contain water and drug:
 - Place the label that matches the drug you are using on the vial
 - Add dry fly food to make a thick slurry, water will be absorbed in minutes
 - Turn on gas and knock out flies. Put up to 40 flies into vial.
 - Put the flies in a box in the  25 °C incubator
- 4 Flies are fed 3 doses of each drug. Flies are fed for a total of 10 days.
- 5 Harvest flies at 10 days and fix in formalin for histological analysis or process appropriately for other studies as needed.