



1 ▼

Dec 15, 2021

Processing human frontal cortex brain tissue for population-scale Oxford Nanopore long-read DNA sequencing SOP V.1

Kimberley J Billingsley^{1,2}, Ramita Dewan², Laksh Malik¹,
Pilar Alvarez Jerez^{1,2}, Stith Kiley¹, Cornelis Blauwendraat^{1,2},
on behalf of the CARD Long-read Team¹

¹Center for Alzheimer's and Related Dementias, National Institute on Aging, Bethesda, Maryland, US A;

²Laboratory of Neurogenetics, National Institute on Aging, Bethesda, Maryland, USA

Kimberley J Billingsley: *These authors contributed equally to this work;

Ramita Dewan: *These authors contributed equally to this work

1



dx.doi.org/10.17504/protocols.io.b2ucqesw

NIH_CARD

Kimberley Billingsley

In development

We are still developing and optimizing this protocol.

Processing human frontal cortex brain tissue for population-scale Oxford Nanopore long-read DNA sequencing SOP

At the NIH's Center for Alzheimer's and Related Dementias (CARD) <https://card.nih.gov/research-programs/long-read-sequencing> we will generate long-read sequencing data from roughly 4000 patients with Alzheimer's disease, frontotemporal dementia, Lewy body dementia, and healthy subjects. With this research, we will build a public resource consisting of long-read genome sequencing data from a large number of confirmed people with Alzheimer's disease and related dementias and healthy individuals. To generate this large-scale nanopore sequencing data we have developed a protocol for processing and long-read sequencing human frontal cortex brain tissue, targeting an N50 of ~30kb and ~30X coverage.

†Correspondence to: Kimberley Billingsley billingsleykj@nih.gov and Cornelis Blauwendraat cornelis.blauwendraat@nih.gov

Acknowledgements:

We would like to thank the Nanopore team (Androo Markham & Hannah Lucio), Circulomics Inc team (Jeffrey Burke, Michelle Kim, Duncan Kilburn & Kelvin Liu) and

the whole CARD long-read team listed below => UCSC: Benedict Paten, Mikhail Kolmogorov, Miten Jain, Kishwar Shafin, Trevor Pesout; NHGRI: Adam Phillippy, Arang Rhie; Baylor: Fritz Sedlazeck; JHU: Winston Timp; NINDS: Sonja Scholz; NIA: Cornelis Blauwendraat, Kimberley Billingsley, Frank Grenn, Pilar Alvarez Jerez, Bryan Traynor, Shannon Ballard, Caroline Pantazis; CZI: Paolo Carnevali.

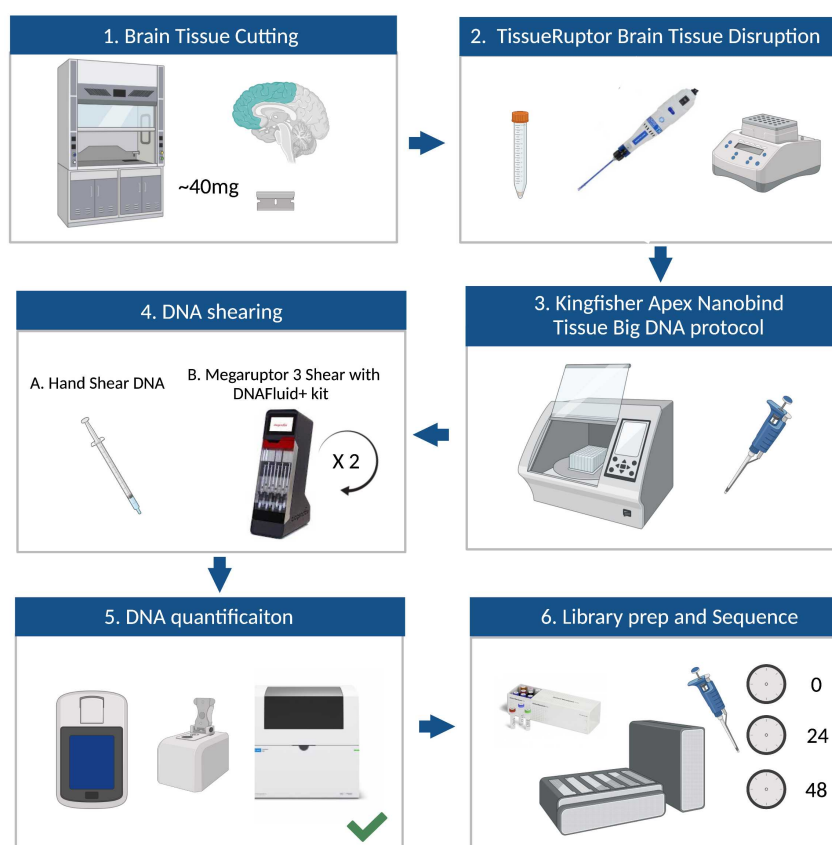


Figure 1. Overview of the HMW Brain DNA extraction and ONT sequencing protocol

Untitled-2.pdf image.png

DOI

dx.doi.org/10.17504/protocols.io.b2ucqesw

Kimberley J Billingsley, Ramita Dewan, Laksh Malik, Pilar Alvarez Jerez, Stith Kiley, Cornelis Blauwendraat, on behalf of the CARD Long-read Team 2021. Processing human frontal cortex brain tissue for population-scale Oxford Nanopore long-read DNA sequencing SOP. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.b2ucqesw>



protocol ,

Dec 14, 2021

Dec 15, 2021

55908

List of reagents/equipment needed:

Consumables:

Description	Supplier
1.5mL DNA LoBind tubes	Qiagen
1.5mL Protein LoBind tubes	Qiagen
Sterile Weigh Boat	
Razor Blade	
Cooling Block	
15mL conical-bottom tube	
TissueRuptor Disposable Tips	
2mL Protein LoBind tubes	
3mm Nanobind Disks	Circulomics -Tissue big DNA kit
1mL Luer-lock Syringes	
1.5" Needles	
DNA Fluid+ Kit	Diagenode
0.2mL thin-wall PCR tubes	
PromethION Flow Cells	Oxford Nanopore Technologies

Reagents:

Description	Supplier
Nanobind Tissue Big DNA Kit	Circulomics
TE pH 8	Millipore Sigma
Isopropyl Alcohol	
Agencourt AMPure XP beads	Agencourt
NEBNext® Companion Module for Oxford Nanopore Technologies® Ligation Sequencing	New England BioLabs
Ligation Sequencing Kit (SQK-LSK110)	Oxford Nanopore Technologies
Flow Cell Wash Kit (EXP-WSH004)	Oxford Nanopore Technologies
Ultra Pure Water	KD Medical
1x dsDNA BR Working Solution	Thermo Fisher
Absolute Ethanol	
Genomic DNA Screen Tape	Agilent

Equipment:



Description	Supplier
Chemical Fume Hood	
Weigh Scale	
Centrifuge 5425R	Eppendorf
ThermoMixer C	Eppendorf
TissueRuptor II	Qiagen
KingFisher Apex	Thermo Fisher
Megaruptor	Diagenode
C1000 Touch Thermocycler	BioRad
NanoDrop 8000	Thermo Fisher
Qubit 4	Thermo Fisher
Tapestation 4200	Agilent
Vortex Genie-T	Scientific Industries
Microfuge	
Magnetic separator (suitable for 1.5mL Eppendorf tubes)	
PromethION 48	Oxford Nanopore Technologies

:

In development

We are still developing and optimizing this protocol.














Part 1: Brain Tissue Cutting (~2.5 hours for 16 samples)

- 1.1 Add dry ice to a ice bucket
- 1.2 Place supplies (sterile weigh boat, razor blade, labeled empty 2mL protein^{5m} LoBind microcentrifuge tubes and cooling block) in dry ice and allow to chill for ~  **00:05:00**
- 1.3 Obtain tissue samples from  **-80 °C** freezer and place in dry ice
- 1.4 Wear all necessary protective equipment (lab coat, face shield, double gloved) and complete the following steps within a chemical fume hood
- 1.5 Weigh labeled empty tube to tare scale, ensuring tube is centered
- 1.6 Working in the chemical fume hood, place the chilled weighing boat on top of the cooling tissue (use right hand to grip blade and cut tissue), use left hand just above, in between tissue and in front of hood, to shield any flying pieces)
- 1.7 Using the chilled razor blade, gently lift the cut tissue piece and transfer to the chilled labeled empty tube and weigh immediately, place back in dry ice immediately to avoid tissue thawing
- 1.8 Add or remove tissue using the method outlined above as required by input specifications ~40mg for frontal cortex, DNA recovery varies based on amount of gray matter vs. white matter)
- 1.9 Dispose of used weighing boat in a burn box and razor blade in the sharps waste container between each sample to prevent inter-sample contamination

and keep all tissue samples on dry ice when not in use

2 Part 2: TissueRuptor Brain Tissue Disruption (~3 hours)

- 2.1 Place 15mL round tubes and cold Buffer CT on ice, chill centrifuge to **4 °C**, and warm ThermoMixer to **55 °C**
- 2.2 Transfer brain tissue from previous steps to 15mL round tubes (keep on ice during the entire disruption process)
- 2.3 Add **750 µL** of cold Buffer CT
- 2.4 Submerge TissueRuptor probe tip in buffer and blend at max speed for **10s** **00:00:10** (place probe tip off to side to be cleaned later)
- 2.5 Transfer homogenate to a 2 mL Protein LoBind microcentrifuge tube including all undisrupted tissue chunks and any foam that forms
- 2.6 Pellet homogenate by centrifuging at 6,000 x g and **4 °C** for **00:05:00**^{5m}. Discard supernatant (pellet may not be visible, so pipette carefully and avoid pipetting from the bottom of tube)
- 2.7 Add **1 mL** of cold Buffer CT and pipette mix 10X with a wide bore P1000 pipette to resuspend tissue
- 2.8 Pellet homogenate by centrifuging at 6,000 x g and **4 °C** for **00:02:00**^{2m}. Discard supernatant (pellet may not be visible, so pipette carefully and avoid pipetting from the bottom of tube)
- 2.9 Pulse vortex for **00:00:01** x 5 times (max setting) to dislodge pellet **1s**

- 2.10 Add  **20 µL** of Proteinase K to the previous pellet
- 2.11 Add  **50 µL** 1X TE pH 8
- 2.12 Add  **60 µL** Buffer CS
- 2.13 Add  **100 µL** Buffer CLE3 and pipette mix 15X with a wide bore P200 pipette
- 2.14 Incubate for  **01:00:00** on a ThermoMixer at  **55 °C** and 900 rpm 1h
- 2.15 Spin on a mini-centrifuge for  **00:00:02** to remove liquid from the cap 2s
- 2.16 Add  **20 µL** of RNaseA and pipette mix 3X with a wide bore P200 pipette
- 2.17 Incubate for  **00:30:00** on a ThermoMixer at  **55 °C** and 900 rpm 30m
- 2.18 Spin the tube on a mini-centrifuge for  **00:00:02** to remove liquid from the cap 2s
- 2.19 Add  **50 µL** Buffer SB and vortex for  **00:00:10** at maximum speed 10s
(transfer any foam that appears)

3 Part 3: KingFisher Apex Nanobind Tissue Big DNA protocol (~2 hours)

3.1 Prepare KingFisher Apex plates as follows:


Plate 1 - Lysis Binding: Sample +  50 µL BL3

Plate 2 - Nanobind Storage: one 3mm Nanobind disk

Plate 3 - CW1 Wash 1:  600 µL Buffer CW1

Plate 4 - CW1 Wash 2:  600 µL Buffer CW1

Plate 5 - CW2 Wash 1:  600 µL Buffer CW2

Plate 6 - CW2 Wash 2:  600 µL Buffer CW2

Plate 7 - Elution:  100 µL Buffer EB

Plate 8 - Tip: KingFisher Flex 96-Tip Combo

3.2 Run KingFisher Apex program “210804_nanobind_tissue_kf_apex_v2.kfx” (KF script available by request from Circulomics Inc)

3.3 After 00:12:00 when the program pauses, add 300 µL IPA 12m

3.4 Collect DNA by transferring eluate to a new 1.5 mL microcentrifuge tube OR if sample is to be sheared, the sample can be transferred to a DNA Fluid+ tube.

3.5 Let the sample rest at room temperature overnight to allow DNA to solubilize

4 Part 4: DNA Shearing (8 hours per 48 samples)

4.1 Hand shear DNA




Hand-shear 10X with 1mL Luer-Lock syringes and 1.5” needles (bringing sample up into needle and depressing plunger counts as 1 cycle)

4.2 Quantify:

1. Quantify the samples with the nanodrop and Qubit and size using the Agilent Tapestation 4200
2. Upload the tapestation reports
3. Add the quantifications
4. Only take forward samples that are > 4.5ug. If a sample does not reach this requirement then repeat DNA extraction from Part 1.

4.3 Megaruptor 3 Shear with DNAFluid+ kit:

2h

1. Transfer all of sample to DNA Fluid+ tubes if not in already
2. Make sample up to  150 µL with nuclease free water
3. Attach the DNA Fluid+ needle onto the tube and push the entire item into the Megaruptor 3 slots until it fits snugly. If running fewer than 8 samples, put the tubes in the 1st and/or 8th slots, working your way in. Samples should always be balanced, if running an odd number of samples, samples can be balanced with an empty corresponding tube. Shear at speed 45 (takes around ~  01:00:00)
4. Once the MR3 shearing is finished, repeat the run by navigating back to the main menu and select speed 45 again (takes around ~  01:00:00)
5. Avoid any vortexing of DNA from this point on to avoid any unnecessary further shearing, instead mix by gently flicking the tube and spin down

5 Part 5: DNA Quantification

5.1 Quantify the samples with the nanodrop and Qubit and size using the Agilent Tapestation 4200

5.2 Upload the tapestation reports and quantifications




5.3 For Step 6 the starting material must be > 4.5ug. If a sample does not reach this requirement then repeat DNA extraction from Part 1.

5.4 DNA can be stored at 4° C for up to four weeks, or  -80 °C indefinitely

Part 6: Library Prep (~6 hours, not including flushing and returning cells)

6.1 A. DNA Repair and End-Prep

20m 30s

1. Put all the necessary reagents on ice to thaw and the Agencourt AMPure XP beads out at room temperature
2. Prepare the following in a 0.2 mL thin-walled PCR tube:
 -  **48 µL** DNA (load 4.5 ug, this will likely be over 48 µL but that is fine, just adjust the amount of beads to match the total volume of this mixture)
 -  **3.5 µL** NEBNext FFPE DNA Repair Buffer (vortex)
 -  **3.5 µL** Ultra II End-prep reaction buffer (vortex)
 -  **3 µL** Ultra II End-prep enzyme mix (do not vortex)
 -  **2 µL** NEBNext FFPE DNA Repair Mix (do not vortex)
3. Mix thoroughly by gently flicking tube or pipetting up and down 10X, and then spin down
4. Using a thermal cycler, incubate samples at  **20 °C** for  **00:05:00** and  **65 °C** for  **00:05:00**
 - Start and pause thermal cycler to allow lid to come to  **85 °C** before putting samples in
5. Allow Thermocycler to cool to  **4 °C** and then remove your samples.
6. Resuspend the AMPure XP beads by vortexing
7. Transfer DNA samples to clean 1.5 mL Eppendorf DNA LoBind tube
8. Add  **60 µL** (or equivalent volume, see step 2) of resuspended beads to the reaction and mix by pipetting up and down 10X
9. Incubate for  **00:05:00** at room temperature
10. Prepare  **500 µL** per sample of fresh 75% ethanol in Nuclease-free water
11. Pellet sample on magnet until eluate is clear and colorless, about  **00:02:00**
12. Pipette off the supernatant and retain, just in case the following quant is uncharacteristically or surprisingly low
13. With the samples remaining on the magnet, wash the beads with  **200 µL** of the ethanol, pipetting on the opposite wall (the goal here is to make sure the beads are fully covered) and making sure not to disturb the

- pellet, count to 3 and remove and discard ethanol
14. Repeat previous step
 15. Spin down and place the tube back on magnet, pipetting off any residual ethanol
 16. Allow to dry for ~ 🕒00:00:30 but do not over-dry
 17. Remove the tube from the magnetic rack and resuspend the pellet in
 🧪61 µL Nuclease-free water, incubate for 🕒00:03:00 at room temperature
 (hold in hands and gently flick every so often)
 18. Pellet the samples on a magnet until eluate is clear and colorless
 19. Remove and retain 🧪61 µL of eluate into a clean 1.5 mL Eppendorf DNA LoBind tube
 20. The sample concentration must be > 40ng/ul. If the sample does not reach this requirement restart from Part 6.
 21. It is possible to store samples at 🌡4 °C overnight at this step if needed

6.2 B. Adapter Ligation and Clean-Up

35m 30s



1. Spin down the AMX-F, Quick T4 ligase, and LNB, then return to ice
 - Do not allow AMX-F to remain at room temperature for too long
2. Thaw LNB at RT and mix by pipetting up and down (vortexing is ineffective due to viscosity)
3. Thaw EB at RT, mix by vortexing, spin down, and place on ice
4. Thaw SFB at RT, mix by vortexing, spin down, and keep at RT
5. In a 1.5 mL Eppendorf DNA LoBind tube, mix the following in order:
 - 🧪60 µL DNA sample (if not 60, add water until it is)
 - 🧪25 µL LNB
 - 🧪10 µL Quick T4
 - 🧪5 µL AMX-F
6. Mix gently by flicking the tube and spin down
7. Incubate the reaction for 🕒00:10:00 at RT
8. During this time, put flow cells out at RT
9. Resuspend beads by vortexing
10. Add 🧪45 µL of resuspended beads to the reaction and mix by flicking
11. Incubate for 🕒00:05:00 at RT
12. Spin down sample and pellet on magnet
13. Pipette off the supernatant and retain, just in case the final elution quant is uncharacteristically or surprisingly low
14. Wash the beads with 🧪250 µL SFB, flick to resuspend and repellent, remove and discard supernatant
15. Repeat previous step
16. Spin down and place the tube back on magnet, pipetting off any residual SFB

17. Allow to dry for ~ 🕒 00:00:30 , but do not over-dry
18. Remove the tube from magnet and resuspend pellet in 25 μL EB, spin down, and incubate for 🕒 00:20:00 at 🌡 37 °C
19. During this time, QC the flow cells (only use flow cells with >6000 pores)
20. Pellet the beads on magnet until eluate is clear and colorless
21. Remove and retain 📦 25 μL of eluate into a clean 1.5 mL Eppendorf DNA LoBind tube (this is the DNA library)
22. Quantify samples on Qubit
23. Reprep library from Part 6 if < 1200 ng.
24. Keep libraries on ice until ready to load on flow cell



6.3 C. Priming and Loading Flow Cell

1h 15m

1. Thaw SBII, FLT, and FB, vortex, and spin down
2. Thaw LBII
3. Add 📦 30 μL of thawed and mixed FLT directly to tube of FB and vortex
4. Expose inlet port on flow cell and draw back a small volume to remove any air bubbles (usually about 20-30 μL , just until a small volume of buffer enters the pipette tip)
5. Flush 📦 500 μL of Priming Mix into the inlet port of the flow cell, being extremely careful to avoid the introduction of air bubbles at the end
6. Wait 🕒 00:05:00
7. During this time, separate the DNA library into three equal aliquots (ideally with 400 ng of DNA each). Bring each aliquot up to 📦 24 μL with EB. (i.e, if your final elution is exactly 1200 ng in 24 μL , move 📦 8 μL to three separate tubes and add 📦 16 μL of EB to each.)
8. Prepare the first library mix for loading:
 - 📦 75 μL SBII (vortex)
 - 📦 51 μL LB (pipette up and down immediately before use)
 - 📦 24 μL DNA library in EB (400 ng)
9. Immediately load all 📦 150 μL of the library mix
10. Close valve to seal inlet port and close PromethION lid
11. Wait 🕒 00:10:00 and then initiate sequencing
12. Ideally, the library quant yielded at least 1200 ng to allow for 3X 400 ng loads, the latter 2 loaded approximately after 24 and 48 hours. However this will vary slightly depending on pore usage, data generated, as well as other factors i.e. if after 24 hours there are still +3000 pores then the sample does not need to be reloaded until 48 hours.
 - To wash and reload a flow cell, begin by thawing Wash Mix (WMX) on ice and Wash Diluent (DIL) at RT (DIL should be vortexed, WMX should NOT be vortexed, only spun)
 - Add 📦 2 μL WMX into 📦 398 μL DIL and pipette mix

- Pause the PromethION runs and export .pdf reports
- Rotate the inlet port cover to reveal inlet port 1
- Using a P1000, insert tip into inlet port and draw back a small volume using the wheel to remove any air (usually around 20-30 µL)
- Load  400 µL Flow Cell Wash Mix into the inlet port, avoiding any introduction of air
- Wait  01:00:00
- Repeat priming steps and reload samples (steps 1 - 12)

6.4 D. Flushing and Recycling Flow Cells (~15 minutes per set of 4 flow cells)

1. Following the completion of the sequencing, flow cells may be removed from the sequencer
2. Place enough absorbent material to take up approximately 4 mL of flush waste
3. Rotate valve to reveal inlet port 1
4. Place flow cell at a 45° angle on the absorbent material and, using a P1000, flush  1 mL of DI water into the inlet port
5. Repeat 3 more times for a total of  4 mL
6. Once complete, close the input port cover and remove all liquid from the waste port
7. Dispose of absorbent material as local biological waste guidelines dictate
8. Return flow cells to clear plastic tray in which it was shipped, making sure to record the flow cell IDs
9. Seal the tray with the sticker provided in the packaging
10. Put the clear plastic lid back on the tray
11. Place the tray back in the packaging
12. Place packaged cells in the returns box (large box can hold 80)
13. Once returns box is filled, follow the instructions [here](#) and follow the prompts to request the box to be sent back to Nanopore

Sequencing results:

7

1. Following 72 hours of sequencing the sample should yield an N50 \geq 30kb with a data output ~ 100-160GB.

Read Length Histogram Basecalled Bases

Estimated N50: 31.85 kb

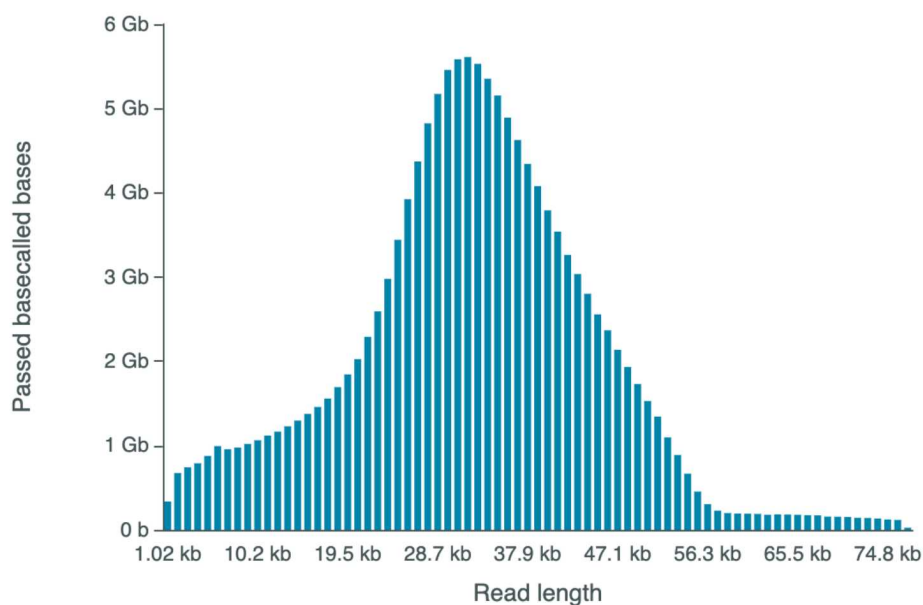


Figure 2. Expected Read Length Histogram:

Mux Scan Grouped

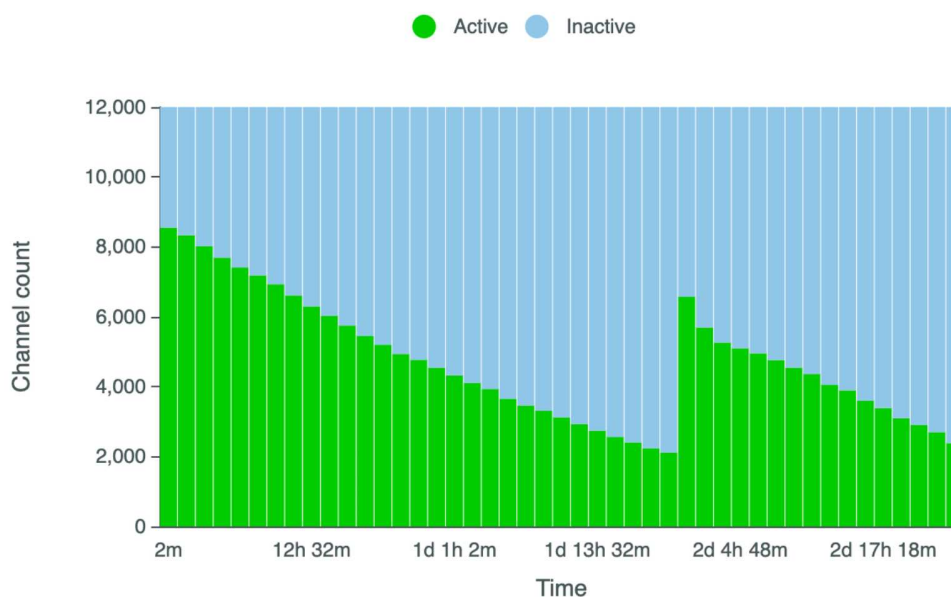


Figure 3. Expected Muxscan (from 2 loads):