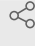




Sep 06, 2022

🌐 Preparing samples for NGS

 In 1 collectionHanqin Li¹, Yogendra Verma¹, Dirk Hockemeyer¹, Frank Soldner²¹University of California, Berkeley; ²Albert Einstein College of Medicine1 *Works for me* Sharedx.doi.org/10.17504/protocols.io.b4nwqvfe

Devin E Snyder

ABSTRACT

This protocol describes a standard procedure used to prepare PCR samples for Next Generation Sequencing (NGS)

DOI

dx.doi.org/10.17504/protocols.io.b4nwqvfe

PROTOCOL CITATION

Hanqin Li, Yogendra Verma, Dirk Hockemeyer, Frank Soldner 2022. Preparing samples for NGS. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.b4nwqvfe>



FUNDERS ACKNOWLEDGEMENT

Aligning Science Across Parkinson's
Grant ID: ASAP-000486

COLLECTIONS

**Genotyping by next generation sequencing**

KEYWORDS

ASAPCRN

LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Feb 03, 2022

LAST MODIFIED

Sep 06, 2022

PROTOCOL INTEGER ID

57782

PARENT PROTOCOLS

Part of collection

[Genotyping by next generation sequencing](#)

MATERIALS TEXT

Item	Vendor	Catalog #
Wizard SV Gel and PCR clean-up system	Promega	A9301
AMPure XP beads alternative	UC Berkeley	
Microseal® 'B' Adhesive Seals	Biorad	MSB1001
96-well PCR plate	GeneMate	T-3152-1

- 1 Run 5 µl PCR product in a 1% agarose gel to confirm successful amplification
- 2 Purify the rest of the PCR product using column based or SPRI-beads based method, like Wizard SV Gel and PCR clean-up system or AMPure XP beads
- 3 Following the instruction of the NGS facility, dilute the purified amplicon to a proper concentration, usually 10 ng/µl
- 4 Send the purified samples to a NGS facility for multiplexing, library prep and NGS.

