



VERSION 8  
SEP 10, 2023

## DNA extraction (BOMB) V.8

Yin-Tse

Huang<sup>1</sup>,

Tsu-Chun Hung<sup>1</sup>

<sup>1</sup>KMU

Mycology Lab



Yin-Tse Huang

### ABSTRACT

DNA extraction (BOMB)

OPEN  ACCESS



#### DOI:

[dx.doi.org/10.17504/protocols.io.n2bvj6mdnlk5/v8](https://dx.doi.org/10.17504/protocols.io.n2bvj6mdnlk5/v8)

**Protocol Citation:** Yin-Tse Huang, Tsu-Chun Hung 2023. DNA extraction (BOMB).

**protocols.io**

<https://dx.doi.org/10.17504/protocols.io.n2bvj6mdnlk5/v8> Version created by Yin-Tse Huang

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
We use this protocol and it's working

**Created:** Sep 09, 2023

**Last Modified:** Sep 10, 2023

**PROTOCOL integer ID:**  
87598

## MATERIALS

### 1. Lysis master mix (870 uL/sample)

A	B
TE buffer	225 uL
Lysis buffer	375 uL
Ammonium acetate	270 uL

### 2. TE buffer


A	B
Tris HCl pH8.0	10mM
EDTA	1mM

### 3. Lysis buffer

A	B
GITC	4M
Tris HCl pH8.0	50mM
SDS	0.5g
EDTA	20mM

## Sample Collection

3m

1 Add  200 µL of **0.5 mm beads** to 2mL screw tube


30s



- 2 Add  200  $\mu\text{L}$  of **1 mm beads** to 2mL screw tube

30s



- 3 Add  870  $\mu\text{L}$  Lysis master mix to 2mL screw tube. The final look:

30s



#### Note

In 11F, 4°C fridge

Lysis master mix: **225  $\mu$ L** of TE buffer + **375  $\mu$ L** of lysis buffer + **270  $\mu$ L** of 10M ammonium acetate

- 4 Collect  20-50 mg of **sample** to 2mL screw tube


1m

#### Note

You can collect up to 100 mg of sample if you can until you bump into the low DNA quality or PCR success rate; by then it means too many inhibitors in the sample and you have to lower the input.

## Sample crush

4m

- 5 Put the 2mL screw tube in mixmill for sample crush, at 3200 rpm  00:04:00


4m

## Note

Remember to balance if you have odd number of samples

## Centrifugation

3m

- 6 Put 2mL screw tube in centrifuge for centrifugation, at this condition:  10 x g, 25°C, 00:03:00

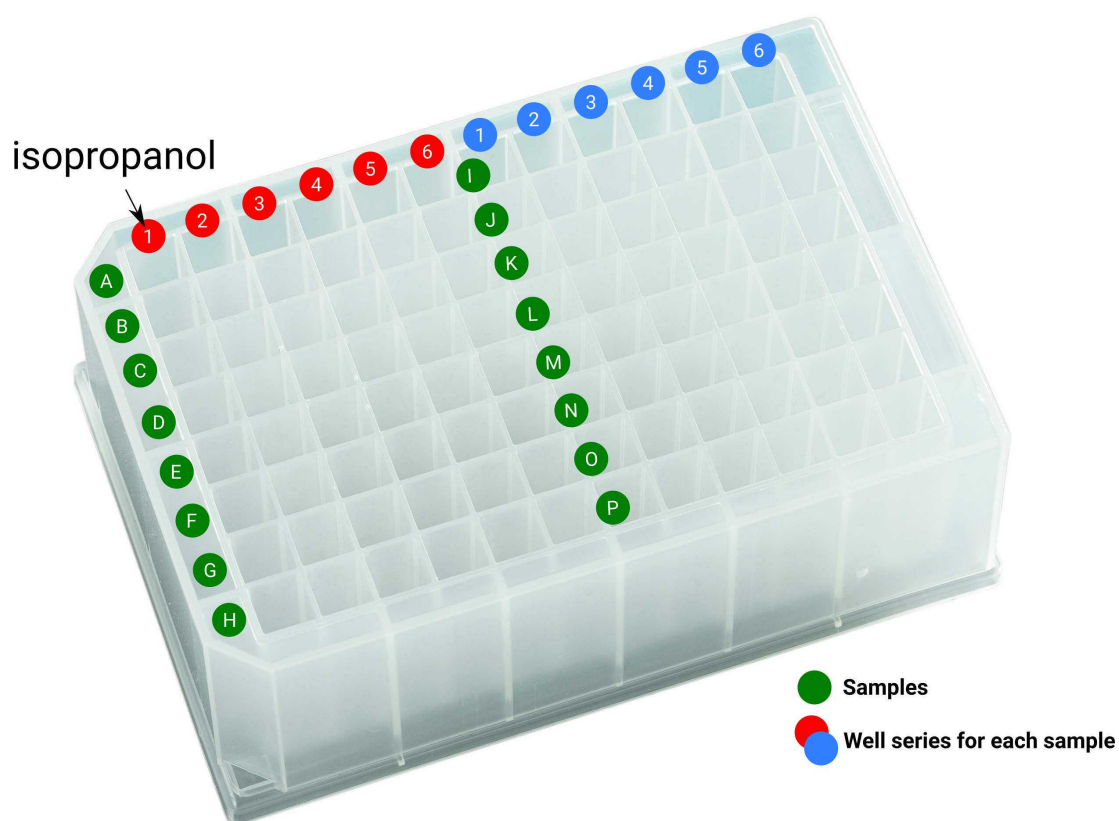
3m


## DNA purification


37m 30s

- 7 Add  350  $\mu$ L of **isopropanol** to the 1st well of 96 well plate

30s




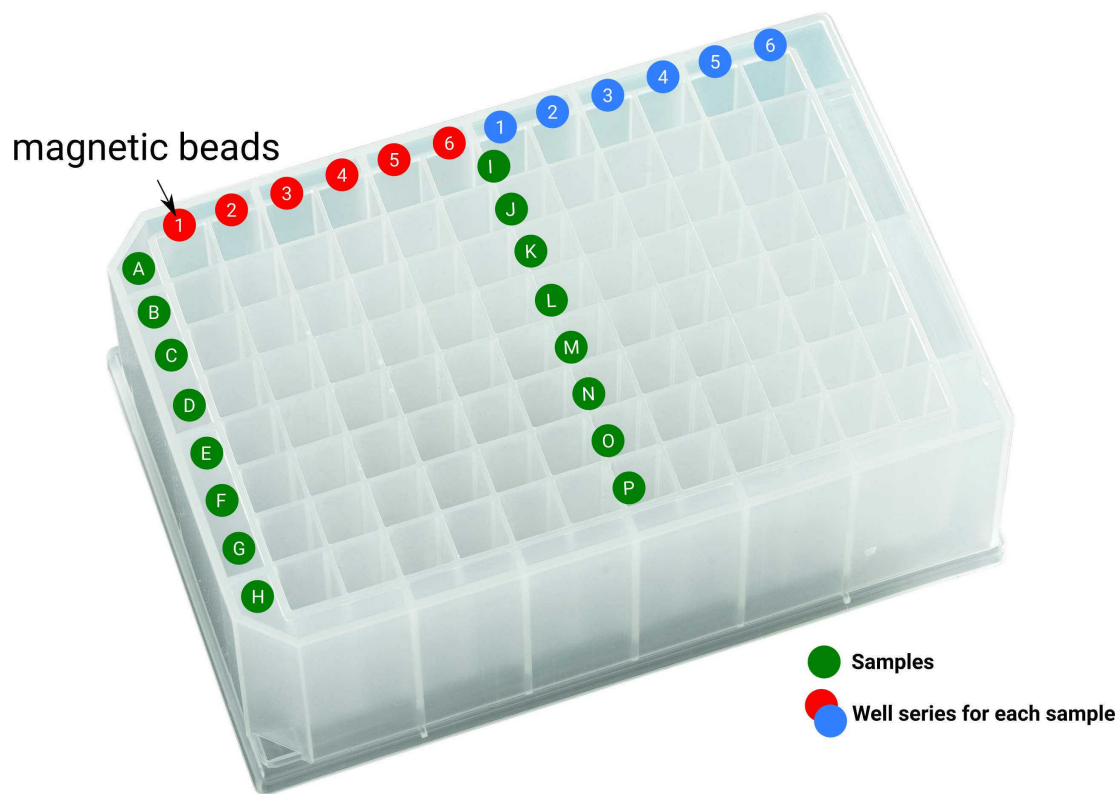- 8 Add  50  $\mu$ L of **magnetic beads (10mg/ml)** to the 1st well of 96 deep well plate




9 Add  400  $\mu\text{L}$  of **isopropanol** to the 2nd well of 96 deep well plate

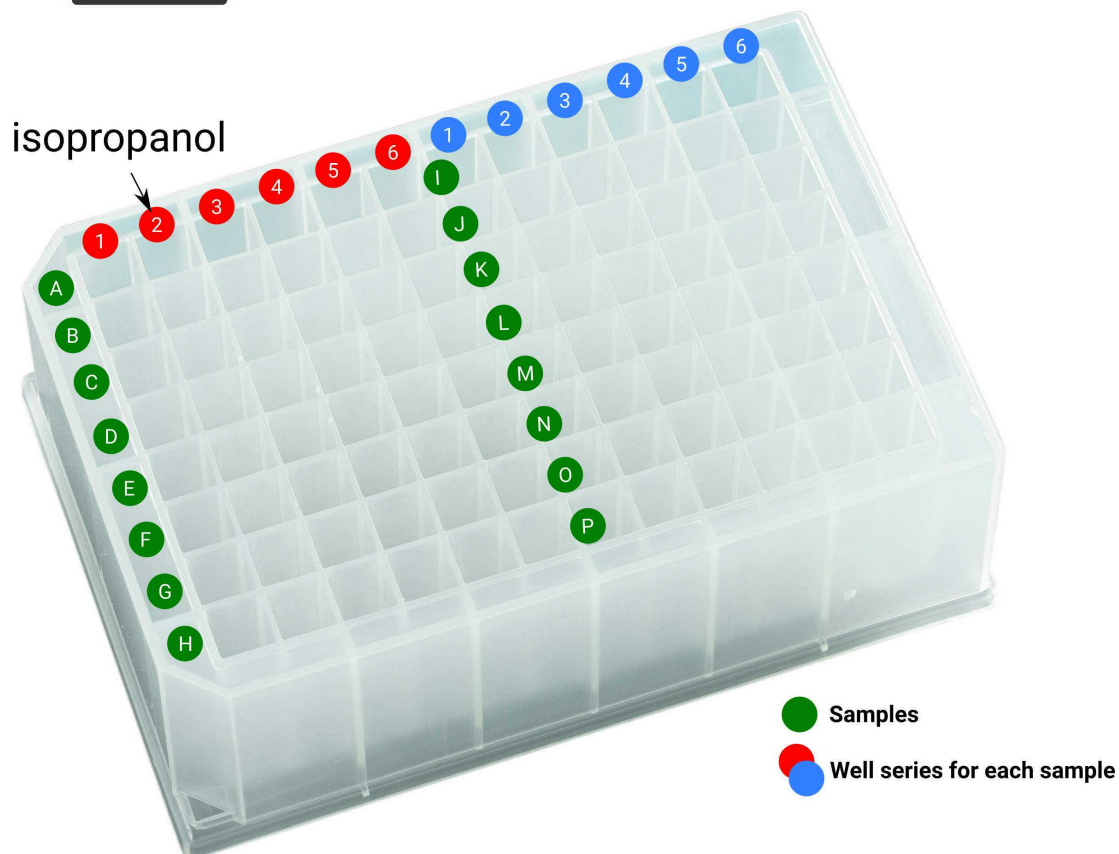30s




10 Add  300  $\mu$ L of **80% ethanol** to the 3rd well of 96 deep well plate

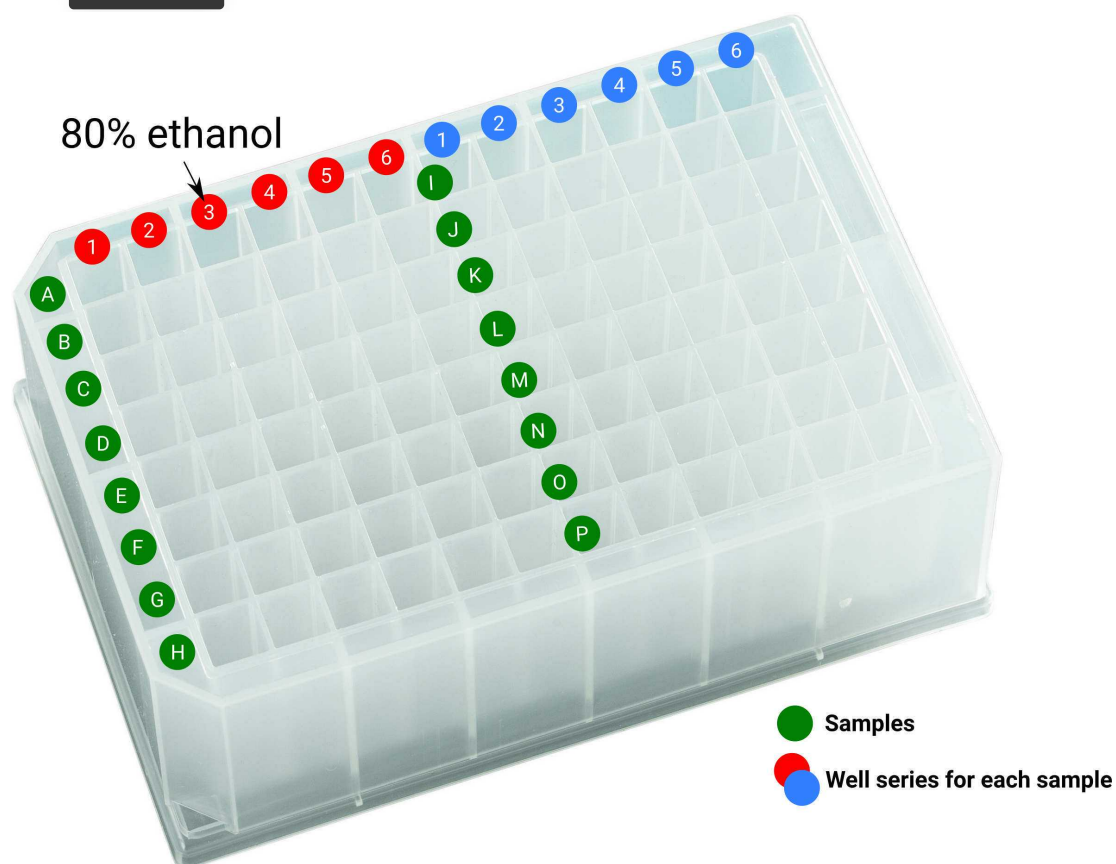
30s




11 Add  300  $\mu$ L of **95% ethanol** to the 4th well of 96 deep well plate

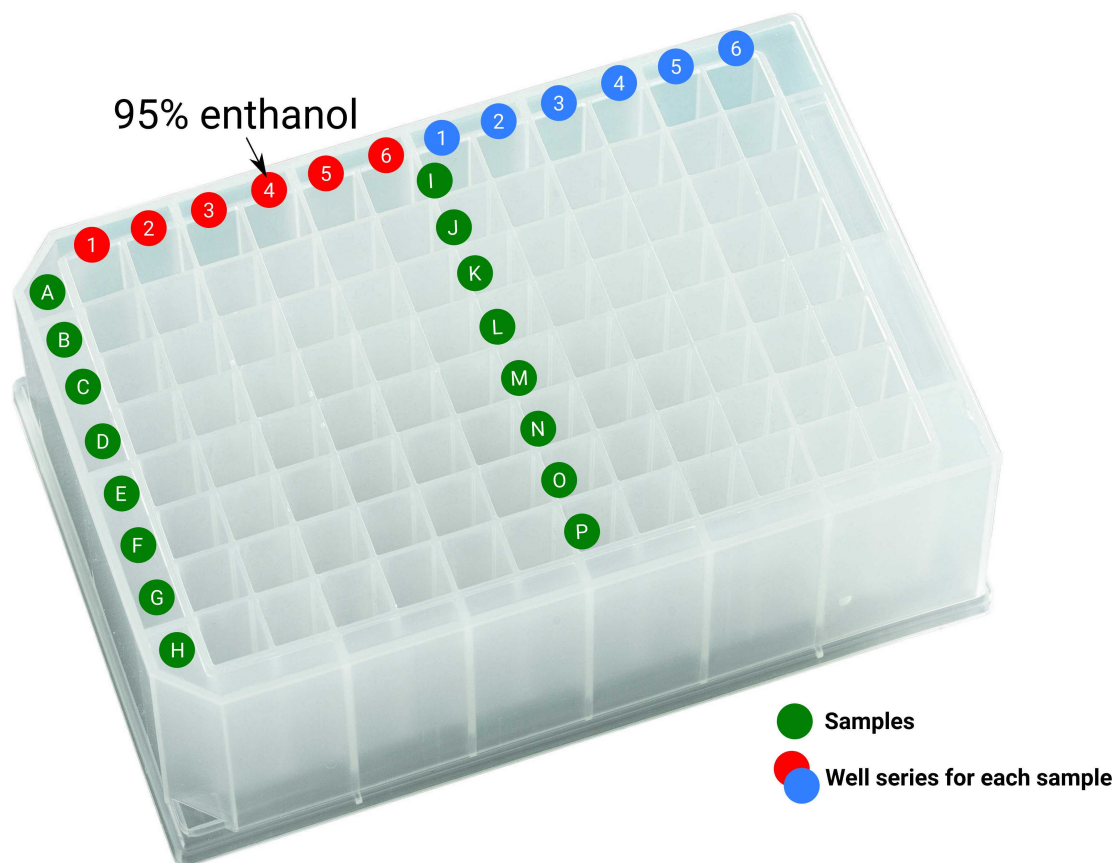30s



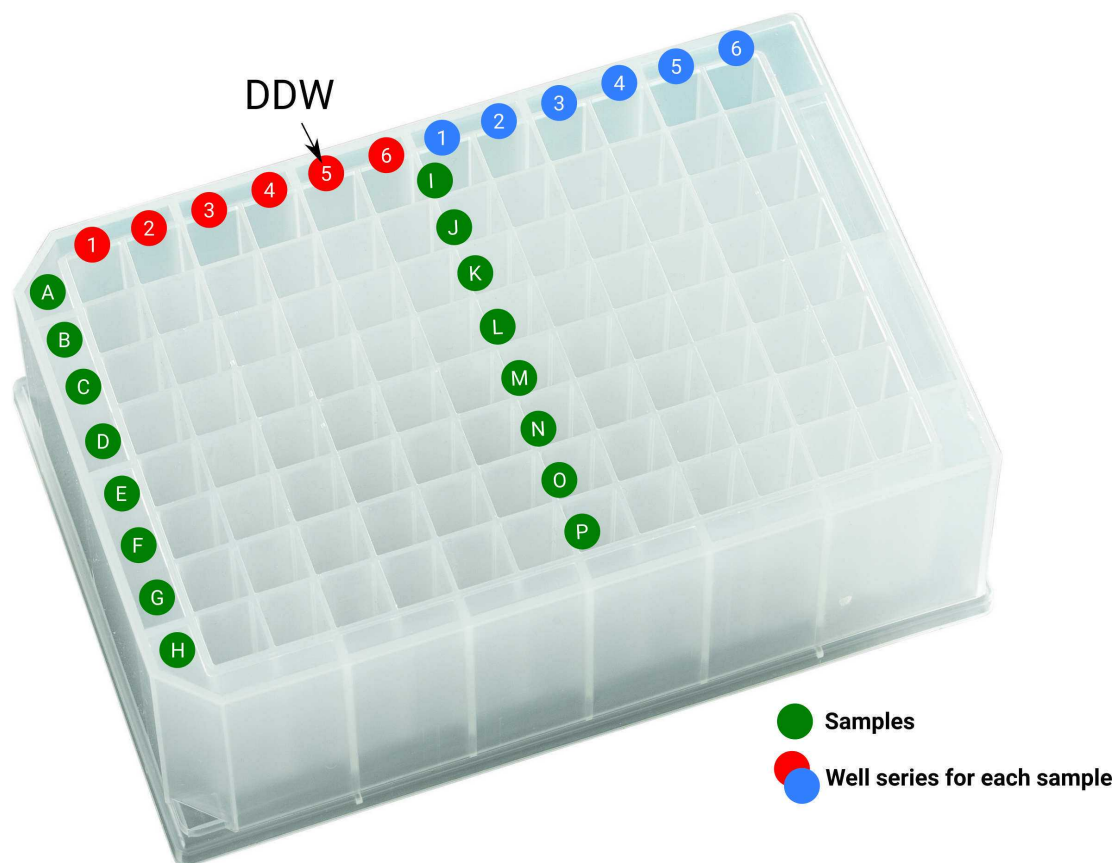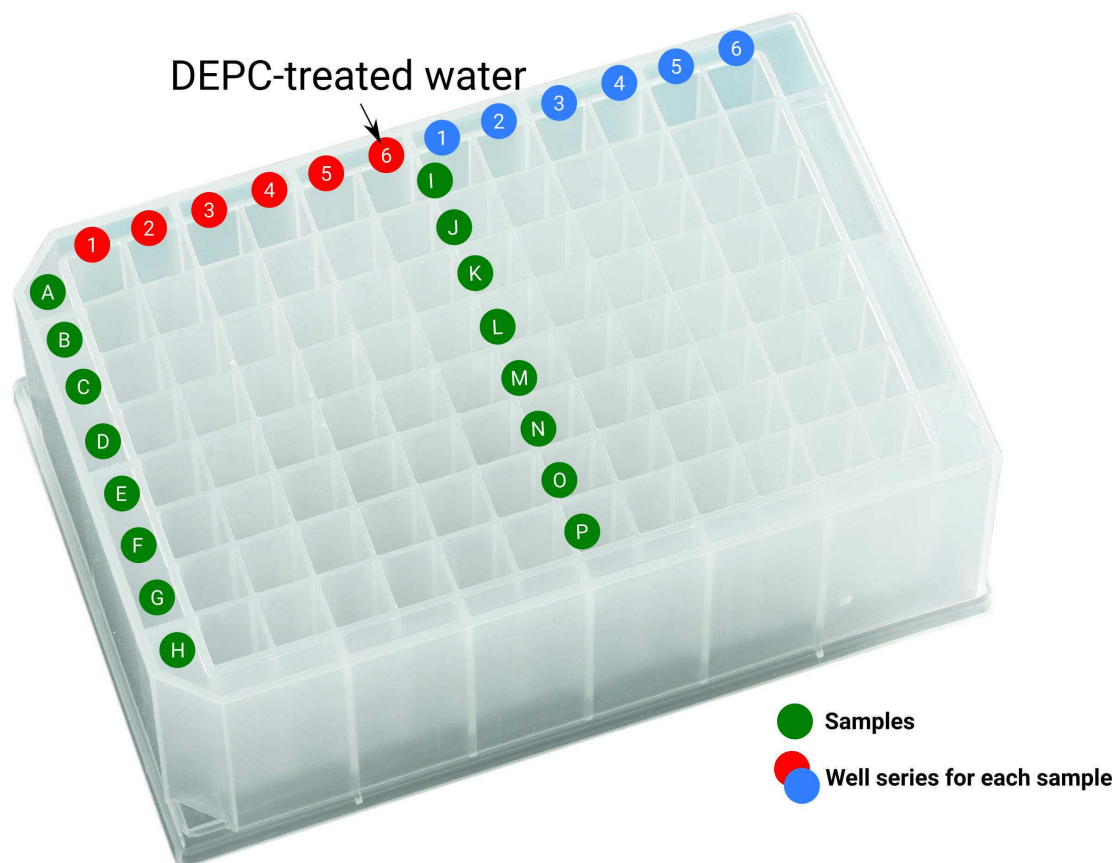
12 Add  300  $\mu\text{L}$  of **DDW** to the 5th well of 96 deep well plate




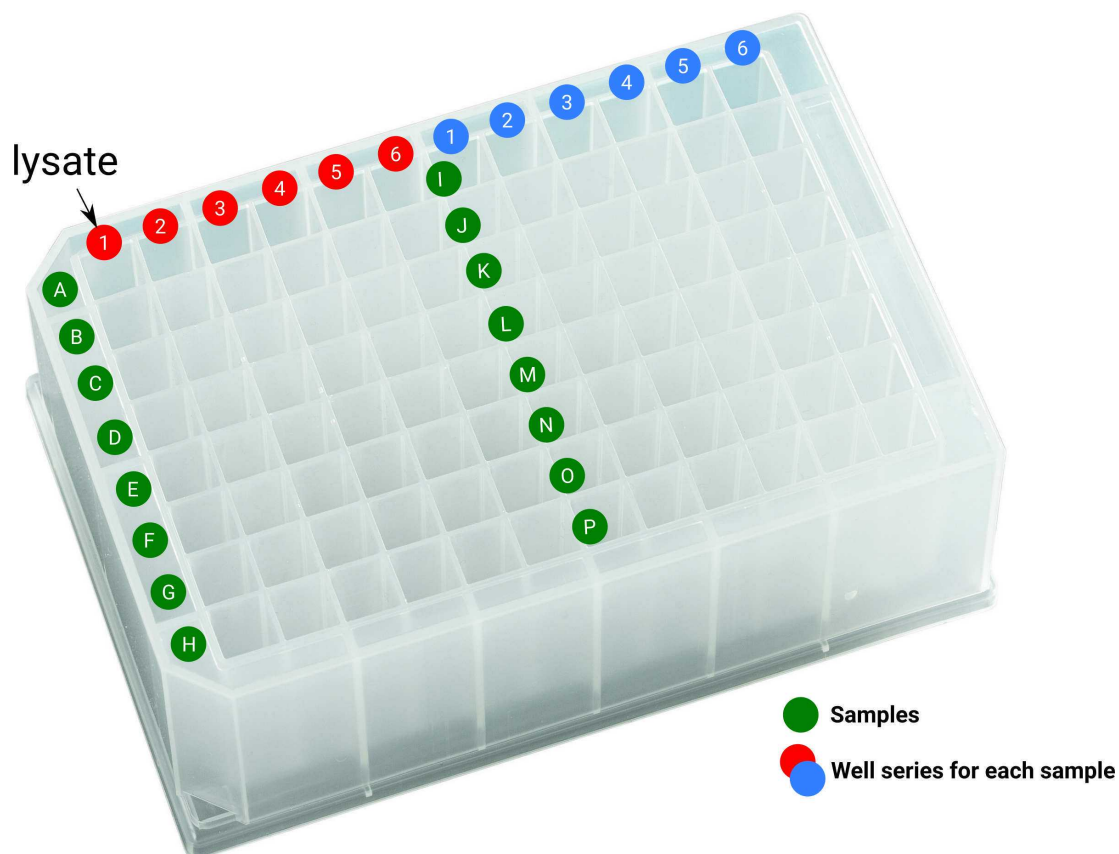


13 Add  100  $\mu\text{L}$  of **DEPC-treated water** to the 6th well of 96 deep well plate

30s



- 14 Add  300-500  $\mu\text{L}$  of the **sample (lysate)** from the 1.5mL centrifuged tube to the 1st well of 96 deep well plate 30s





#### Note

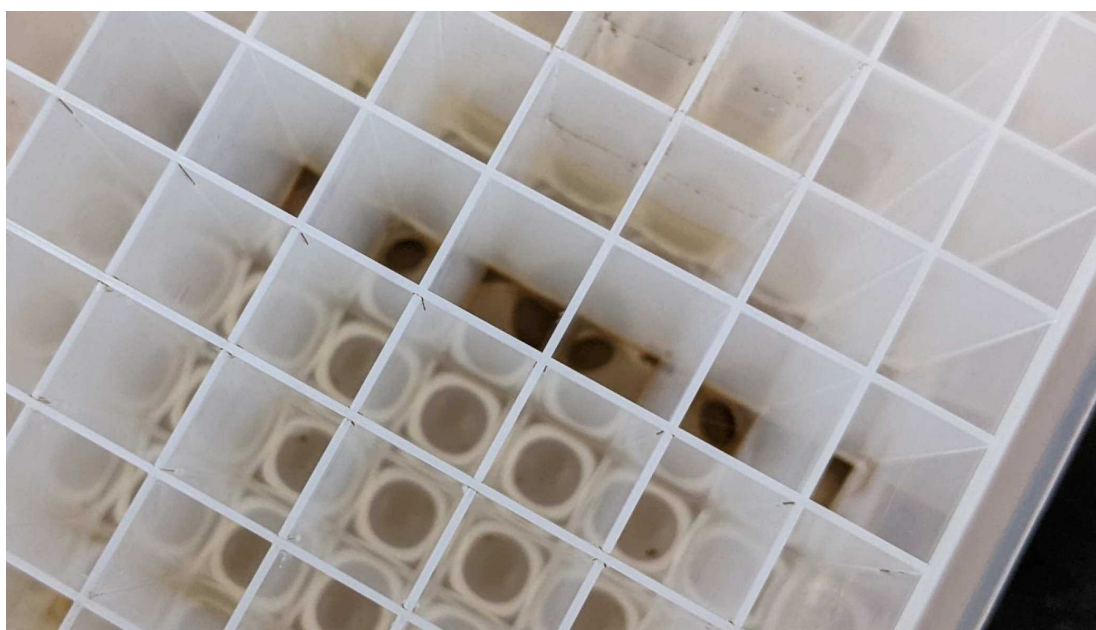Pipetting **as many lysate as you can**, as long as it's free of any cell debris (no solids in your tip)

- 15 Put the prepared 96 deep well plate in the automated DNA extraction machine and select the BOMB protocol
- 16 After the extraction is done, put on the 96 magnetic plate to pellet the magnetic bead residues.

34m



- 17 Collect  100  $\mu\text{L}$  of the **eluted sample** (avoid getting magnetic bead) as the DNA template for downstream experiments

