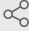





Sep 09, 2022

Preparation of LRRK2 RCKW cryo-EM grids

 Forked from [Preparation of LRRK2 RCKW trimer cryo-EM grids](#)**Mariusz Matyszewski**¹¹Department of Cellular and Molecular Medicine, University of California, San Diego, La Jolla, CA 920931 *Works for me* Sharedx.doi.org/10.17504/protocols.io.kxygxp7ywl8j/v1 **Mariusz Matyszewski**

ABSTRACT

This is Leschziner's Lab updated protocol for making cryo-EM grids for LRRK2 RCKW. This protocol, when using lower protein concentration, results in better monomer and dimer formation than the old protocol.

DOI

dx.doi.org/10.17504/protocols.io.kxygxp7ywl8j/v1

PROTOCOL CITATION

Mariusz Matyszewski 2022. Preparation of LRRK2 RCKW cryo-EM grids.
protocols.io
<https://protocols.io/view/preparation-of-lrrk2-rckw-cryo-em-grids-bryqm7vw>



FUNDERS ACKNOWLEDGEMENT

ASAP
Grant ID: ASAP-000519

FORK NOTE

FORK FROM

Forked from [Preparation of LRRK2 RCKW trimer cryo-EM grids, Mariusz Matyszewski](#)

KEYWORDS

cryo-EM, LRRK2, structural biology, ASAPCRN

LICENSE

_____ This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Jan 31, 2021

LAST MODIFIED

Sep 09, 2022

PROTOCOL INTEGER ID

46832

MATERIALS TEXT

LRRK2 Buffer:

- **[M]20 millimolar (mM) HEPES pH 7.4**
- **[M]80 millimolar (mM) NaCl**
- **[M]0.5 millimolar (mM) TCEP**
- **[M]5 % volume Glycerol**
- **[M]2.5 millimolar (mM) MgCl₂**
- **[M]20 micromolar (μM) GDP**

Note: please change salt as needed to maintain final salt of 80 mM NaCl

SAFETY WARNINGS

For hazard information and safety warnings, please refer to the SDS (Safety Data Sheet).

Take proper precautions while freezing grids.

BEFORE STARTING

Decide which protein concentration to use, and create the proper LRRK2 buffers in order to obtain the right salt concentration (80 mM NaCl).

Preparing Sample

10m

- 1 Spin down purified LRRK2 RCKW. 🌀 **10000 rcf, 4°C, 00:10:00** , (can be faster)
Leave protein on ice afterwards.

For best results, reduce the amount of time between spinning and freezing samples.

Freezing Grids

20s

- 2 Plasma clean grids.
We used UltrAuFoil Holey Gold 1.2/1.3 300 mesh grids and plasma cleaned them in the Solarus II (Gatan) using the QuantiFoil Au preset.
- 3 Dilute samples to desired concentration in the **LRRK2 buffer**. Make sure final salt is at 80 mM NaCl.
For best results, make 📏 **10 µL** samples, good for freezing 2 grids. This is to minimize time spent outside of storage buffer, reducing aggregation.

High concentrations favor dimers and trimers, while lower concentrations favor monomers.

3.1 (Optional) If using inhibitors, let them incubate on ice in the final LRRK2 buffer before plunge freezing.

- 4 Apply protein to grids and plunge freeze.
We used a Vitrobot (FEI) to blot away excess sample and plunge freeze
- 5 Store grids in liquid nitrogen until ready for imaging.