


Apr 13, 2026

# 🌐 Extraction of DNA from canine feces for detection of *Campylobacter jejuni* using the MagMax™ CORE Nucleic Acid Purification Kit on KingFisher™ Flex Instrument

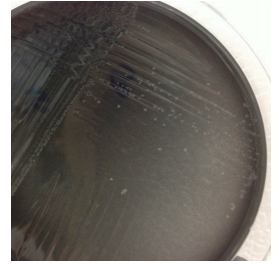

DOI

<https://dx.doi.org/10.17504/protocols.io.8epv5j1b5l1b/v1>

Sara D. Lawhon<sup>1</sup>, Ching-Yuan Yang<sup>1</sup>, Melanie Prarat<sup>2</sup>, Jing Wu<sup>1</sup>

<sup>1</sup>Texas A&M University - College Station; <sup>2</sup>Ohio Animal Disease Diagnostic Laboratory

Vet LIRN



Sara D. Lawhon

Texas A&M University - College Station

OPEN  ACCESS



**DOI:** <https://dx.doi.org/10.17504/protocols.io.8epv5j1b5l1b/v1>

**Protocol Citation:** Sara D. Lawhon, Ching-Yuan Yang, Melanie Prarat, Jing Wu 2026. Extraction of DNA from canine feces for detection of *Campylobacter jejuni* using the MagMax™ CORE Nucleic Acid Purification Kit on KingFisher™ Flex Instrument. protocols.io <https://dx.doi.org/10.17504/protocols.io.8epv5j1b5l1b/v1>

**Manuscript citation:**

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working

**We use this protocol and it's working**

**Created:** September 27, 2022

**Last Modified:** April 13, 2026

**Protocol Integer ID:** 70577

**Keywords:** campylobacter jejuni, dna extraction, manual dna extraction method, extraction of dna, kingfishertm flex instrument campylobacter jejuni, feces for the detection, feces for detection, campylobacter, rapid detection of animal, foodborne illness, extraction, significant cause of foodborne illness, cause of foodborne illness, magmaxtm core nucleic acid purification kit, rapid detection, veterinary patient, fece

**Funders Acknowledgements:**

U.S. Food & Drug Administration Veterinary Laboratory Investigation and Response Network  
Grant ID: U18FD005013

U.S. Food & Drug Administration Veterinary Laboratory Investigation and Response Network  
Grant ID: U18FD006171

U.S. Food & Drug Administration Veterinary Laboratory Investigation and Response Network  
Grant ID: U18FD006446

U.S. Food & Drug Administration Veterinary Laboratory Investigation and Response Network  
Grant ID: U18FD006664

U.S. Food & Drug Administration Veterinary Laboratory Investigation and Response Network  
Grant ID: U18FD007242

## Disclaimer

Reference to any commercial materials, equipment, or process does not in any way constitute approval, endorsement, or recommendation by the Food and Drug Administration.

The authors thank the members of the U.S. Food & Drug Administration Veterinary Laboratory Investigation and Response Network who reviewed and tested the method during multiple interlaboratory comparison exercises and blinded method tests.

## Abstract

*Campylobacter jejuni* is a significant cause of foodborne illness in humans and animals. Rapid detection of animals shedding *C. jejuni* is a key step in limiting the scope of outbreaks and reducing the potential exposure of humans and other animals. Manual DNA extraction methods, although effective, can be labor-intensive and time-consuming. The utilization of automated or semi-automated methods can significantly reduce the time needed for DNA extraction and significantly increase the number of samples that can be processed in a given amount of time. This method was developed to provide a standard method for semi-automated DNA extraction from feces for the detection of *C. jejuni* in veterinary patients.

Validation data (in-house and by an independent laboratory via collaborative study such as Randomized Blinded Method Test) are available upon request.

## Image Attribution

Sara Lawhon



## Guidelines

All manufacturers' intended use, safety information, storage guidelines, and recommendations should be followed.

*Campylobacter jejuni* is infectious to humans and is a biological safety level 2 organism. All work with specimens that potentially contain *C. jejuni* should follow the best practices outlined in the most current version of Biosafety in Microbiological and Biomedical Laboratories (BMBL) and institutional guidance.

Reporting of positive test results should follow all local, state, and federal regulations.

Follow the MagMax<sup>TM</sup> CORE Nucleic Acid Purification Kit manufacturer's recommendations, warnings and procedural guidelines.

## Materials

### For DNA Extraction

1X PBS (Life Technologies cat#10010-023)

Applied Biosystems™ MagMAX™ CORE Nucleic Acid Purification Kit (ThermoFisher Scientific, cat#A32700/A32702)

VetMAX™ Xeno™ Internal Positive Control DNA (ThermoFisher cat#s A29762/A29764)

Adhesive PCR Plate Foils, or equivalent (ThermoFisher Scientific, cat# AB0626)

Plasticware - please confirm that plasticware works with your equipment.

KingFisher™ 96 tip comb for DW magnets, 100 combs (ThermoFisher Scientific, cat#97002534)

KingFisher™ Flex Microtiter Deepwell 96 plates (ThermoFisher Scientific, cat#95040460)

KingFisher™ 96 KF microplates (200 µL), (ThermoFisher Scientific, cat#97002540)\*

\*Note: Only works with KingFisher™ Flex 96 Heating Block (Cat. No. 24075420). Not to be used with KingFisher™ Flex 96 Deep-Well Heating Block (Cat. No. 24075430).

### From any manufacturer:

Microcentrifuge tubes - 1.5-mL and 2-mL tubes sterile

50 mL conical tubes - sterile

Pipette tips (1 mL, 200 µL, 20 µL and 10 µL)

Serological pipette, 25 mL

### For PCR

TaqMan Fast Virus 1-step Master Mix (ThermoFisher Cat. Nos. 4444432/4444434/4444436)

VetMAX™ Xeno™ Internal Positive Control - VIC™ Assay (ThermoFisher Scientific A29765 or A29767)

Positive control DNA from *C. jejuni* ATCC 33291, ATCC 33292, or ATCC 33560

Negative control DNA from *E. coli* ATCC 25922

|  | A          | B                                                       |
|--|------------|---------------------------------------------------------|
|  | Oligo Name | 5'-3' Sequence                                          |
|  | gyrA-F     | AAGATACGGTCGATTTTGTTC                                   |
|  | gyrA-R     | CTACAGCTATACCACTTGAACCATTTAATA                          |
|  | gyrA-P     | 5'-<br>[FAM]TGATGGTTCAGAAAGCGAACCTGATGTTTT[BHQ<br>1]-3' |

Primer and probe sequences.

Primers and probe are diluted to a working concentration of 10 µM with a final concentration in the total reaction volume of 0.5 µM.

### Equipment needed:



KingFisher™ Flex Magnetic Particle Processor (ThermoFisher Scientific, cat#5400110)

Biotang Inc Microplate Shaker, or equivalent titer plate shaker (for mixing beads with samples; Fisher Scientific™ cat#50-751-4965)

Applied Biosystems™ 7500 Fast Real-Time PCR System or similar real-time PCR machine

Scale for weighing 8 g of feces

Pipettes for appropriate volumes

Centrifuge for microcentrifuge tubes capable of speeds of up to 16,000 X g

Benchtop centrifuge with plate adaptors (for lysate preparation in plates)

Laboratory mixer, Vortex or Equivalent

Optional: Qiagen (MoBio) Vortex Adapter (Cat#13000-V1-24)

### Reference for primers and probe

Iijima Y, Asako NT, Aihara M, Hayashi K. (2004) Improvement in the detection rate of diarrhoeagenic bacteria in human stool specimens by a rapid real-time PCR assay. J Med Microbiol. 2004 Jul;53(Pt 7):617-22. PMID: 15184531

### Protocol materials

⊗ VetMAX Xeno Internal Positive Control ND **Thermo Fisher Scientific Catalog #A29762/A29764**

⊗ 1X PBS **Life Technologies Catalog #10010-023**

⊗ MagMAX™ CORE Nucleic Acid Purification Kit **Thermo Fisher Scientific Catalog #A32700/A32702**

### Safety warnings

! All manufacturers' intended use, safety information, storage guidelines, and recommendations should be followed.

Appropriate safety and personal protective equipment should be worn while working with chemicals including a suitable lab coat, disposable gloves, and protective eyewear.

*Campylobacter jejuni* is infectious to humans and is a biological safety level 2 organism. Work with specimens that may contain this organism should follow all applicable institutional safety guidelines.

## Background

- 1 Fecal samples should be stored at  4-8 °C refrigerated until processed.
- 2 Samples should be **processed within 72 hours of receipt at the laboratory.**
- 3 If performing bacterial culture, bacterial culture plates and enrichment broth should be incubated in a microaerophilic environment at  42 °C + or - 2 .
- 4 Remaining samples from the associated culture method can be used for molecular detection.
- 5 Before using the MagMax™ CORE Nucleic Acid Purification Kit, invert bottles of solutions and buffers to ensure thorough mixing.
- 6 In this protocol, you will mix samples with reagents from the MagMax™ CORE Nucleic Acid Purification Kit by either using a plate shaker or by pipetting up and down.  
  
If you are using a plate shaker, you will want to determine the maximum setting of your plate shaker prior to running samples.
- 6.1 First verify that the plate fits securely on the plate shaker.
- 6.2 Test the maximum setting that you can use on the plate shaker by adding  1 mL of water to each well of a test plate and then covering the plate with sealing foil.  
  
Next, determine the maximum setting that you can use on your shaker without any of the water splashing onto the sealing foil.
- 7 To prevent cross-contamination, 1) cover the plate or tube strip during the incubation and shaking steps, to prevent spill-over and 2) carefully pipet reagents and samples to avoid splashing.
- 8 To prevent nuclease contamination, wear laboratory gloves throughout all procedures, use nucleic acid-free pipette tips to handle reagents, and decontaminate the lab benches and pipettes before use.

9 This protocol uses

 VetMAX Xeno Internal Positive Control ND Thermo Fisher Scientific Catalog #A29762/A29764

as a control for PCR inhibition.

#### Note

Note: This reagent will be added to the MagMAX™ CORE Lysis Solution as a DNA extraction control and/or after DNA extraction is complete to serve as the PCR inhibition control. If the laboratory has a system in place for DNA extraction and PCR inhibition, that method can be substituted in this protocol.


10 Prior to using this protocol, the script for the MagMAX™ CORE kit should be loaded onto the instrument. The script is: **KingFisher Flex heated script: MagMAX CORE\_Flex.bdz**

## Sample preparation



2m 30s

11 Use  8 g of feces in 50 mL conical tubes

12 Add  32 mL  1X PBS Life Technologies Catalog #10010-023 to  8 g of feces. The feces in PBS will hereafter be called the feces-PBS slurry.

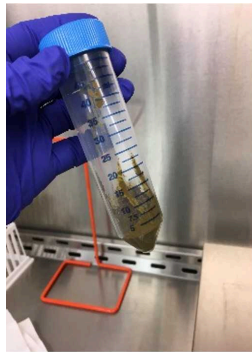

8 g feces

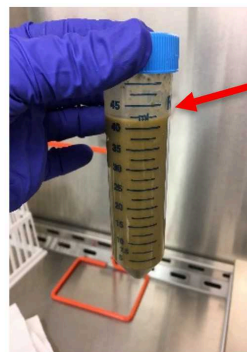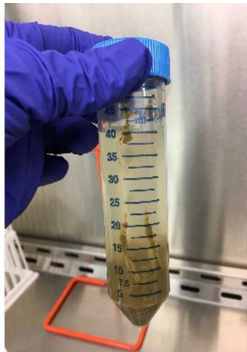

Add **32 mL PBS** to 8 g feces to prepare the feces-PBS slurry and allow the tube to sit for 2 minutes.

Vortex slurry for 30 seconds. Use slurry for direct PCR.

When transferring sample from the feces-PBS slurry pipette from the meniscus at the red arrow.

Making the Feces-PBS slurry.

12.1 Allow samples to sit for  00:02:00 after adding the  32 mL PBS

2m

12.2 Vortex samples for  00:00:30 to prepare fecal slurry

30s

## PCR Method - DNA Extraction

1m

13 Extract DNA from the feces using the

 MagMAX™ CORE Nucleic Acid Purification Kit **Thermo Fisher**  
Scientific Catalog #A32700/A32702

following the manufacturer's protocol (outlined below).

14 Before processing the fecal slurry, prepare the Processing Plates.

|  | A            | B              | C          | D                                                                                                                               | E               |
|--|--------------|----------------|------------|---------------------------------------------------------------------------------------------------------------------------------|-----------------|
|  | Plate ID     | Plate Position | Plate Type | Reagent                                                                                                                         | Volume per well |
|  | Sample Plate | 1              | Deep Well  | *The sample plate will be set up in subsequent protocol steps and will be loaded on the KingFisher™ Flex Instrument at Step 24. |                 |
|  | Wash Plate 1 | 2              | Deep Well  | MagMAX™ CORE Wash Solution 1                                                                                                    | 500 µL          |
|  | Wash Plate 2 | 3              | Deep Well  | MagMAX™ CORE Wash Solution 2                                                                                                    | 500 µL          |
|  | Elution      | 4              | Standard   | MagMAX™ CORE Elution Buffer                                                                                                     | 90 µL           |
|  | Tip Comb     | 5              | Standard   | Place a tip comb in the plate                                                                                                   |                 |

Table 1. Plate setup: KingFisher™ Flex instrument

15 Before processing the fecal slurry, prepare the Bead-Proteinase K mix.



- 15.1 Vortex the MagMAX™ CORE Magnetic Beads thoroughly to ensure that the beads are fully resuspended.
- 15.2 Combine the following components for the required number of samples plus 10%.

| A                           | B                 |
|-----------------------------|-------------------|
| Component                   | Volume per sample |
| MagMAX™ CORE Magnetic Beads | 20 µL             |
| MagMAX™ CORE Proteinase K   | 10 µL             |
| Total Bead/PK Mix           | 30 µL             |

Table 2. Preparation of the Magnetic Beads and Proteinase K

#### Note

The manufacturer recommends that you prepare new Bead/PK mix for each processing run. If necessary, the Bed PK mix can be stored at  4 °C for up to 1 week.

- 15.3 Add the Bead/PK Mix to the wells in Sample Plate

- 16 Prepare the CORE Lysis Solution.

- 16.1 Combine the following components for the required number of samples plus 10%.

| A                                                      | B                 |
|--------------------------------------------------------|-------------------|
| Component                                              | Volume per sample |
| MagMAX™ CORE Lysis Solution                            | 450 µL            |
| VetMAX™ Xeno™ Internal Positive Control DNA            | 4 µL              |
| Total Lysis Solution and Internal Positive Control DNA | 454 µL            |

Table 2. Preparation of the Lysis Solution

### Note


Validation of the protocol used  $4\ \mu\text{L}$  VetMAX™ Xeno™ Internal Positive Control DNA per reaction instead of the  $2\ \mu\text{L}$  as noted in the manufacturer's protocol for consistency with prior manual extraction methods. Laboratories may choose to validate the use of  $2\ \mu\text{L}$  VetMAX™ Xeno™ Internal Positive Control DNA per reaction internally.

16.2 Mix by inverting the tube or bottle at least 10 times.

### Note

(Optional) Store Lysis Solution at room temperature for up to 24 hours.

17 To a microcentrifuge tube, add  $1\ \text{mL}$  of feces-PBS slurry from Sample Preparation Step 12.

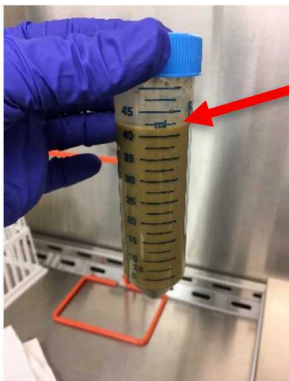

When transferring sample from the feces-PBS slurry pipette from the meniscus at the red arrow.

Pipetting the Feces-PBS slurry to a microcentrifuge tube.

18 Centrifuge the sample at 100 x g, Room temperature, 00:01:00

1m



19 Samples can be processed either in tubes (Step 20) or plates (Step 21)

20 Processing in tubes

20.1 Add 454  $\mu\text{L}$  of the Lysis Solution to a new 2-mL tube.

20.2 Transfer 200  $\mu\text{L}$  of the supernatant of the centrifuged fecal slurry to the tube.

20.3 Vortex vigorously for 00:03:00 or until the sample is suspended.

3m

#### Note

Optional - The sample in the microcentrifuge tube can be placed in a vortexer with Qiagen (MoBio) Vortex Adapter (Cat#13000-V1-24)

20.4 Centrifuge the sample at  15000 x g, Room temperature, 00:02:00

2m



20.5 Without disturbing the pellet, transfer  200  $\mu\text{L}$  of the the supernatant (clarified lysate) to the Sample Plate containing the Bead/PK Mix from Step 15 and proceed to Step 22.

21

Processing samples in Plates

21.1 Add  454  $\mu\text{L}$  of the Lysis Solution to the appropriate wells of a deep-well plate.

21.2 Transfer  200  $\mu\text{L}$  of the supernatant of the fecal slurry to the tube.

21.3 Seal the plate with sealing foil.

21.4 Shake the plate at moderate speed for 5 minutes.

21.5 Centrifuge at  3000 x g, Room temperature, 00:05:00 .

5m

21.6 Remove the supernatant (clarified lysate) without disturbing the pellet and transfer  200  $\mu\text{L}$  to the Sample Plate containing the Bead/PK Mix from Step 15 and proceed



to Step 22.

22

Mix the clarified lysate (sample) and bead mix for  00:02:00 at

 Room temperature .

2m

23

Add  350  $\mu$ L of MagMAX™ CORE Binding Solution to the clarified lysate-bead mixture in wells of the Sample Plate.

24

Immediately proceed to process the samples on the KingFisher™ Flex instrument

25

Select the appropriate script on the instrument. **KingFisher Flex heated script: MagMAX CORE\_Flex.bdz**

 00:27:00

27m

#### Note

Processing on the machine takes approximately 27 minutes.

## PCR Method - Setting up and running the PCR reactions

26

The probe is labeled with FAM with Black Hole Quencher 1 (BHQ1) so detection parameters appropriate for FAM or SYBR should be used. The developing laboratory uses Applied Biosystems™ VetMAX™ Xeno™ Internal Positive Control DNA (ThermoFisher Scientific A29762 or A29764) and the VetMAX™ Xeno™ Internal Positive Control - VIC™ Assay (ThermoFisher Scientific A29765 or A29767) to monitor for PCR inhibition.

## PCR Method - DNA Extraction

1m

27

After samples have finished processing, proceed with PCR or store samples at  4 °C until ready to perform PCR.

## PCR Method - Setting up and running the PCR reactions

28 Create run on ABI 7500 Fast machine. Setting up the ABI 7500 Fast Thermocycler: Refer to the manual provided for running the machine. Reaction times are as follows in Table 1.

Equipment

7500 Fast Real-Time PCR System, desktop

NAME

Applied Biosystems

BRAND

4351107

SKU

|  | A       | B                      | C           | D              | E                | F      |
|--|---------|------------------------|-------------|----------------|------------------|--------|
|  | Stage   | Function               | Temperature | Duration (sec) | Repeats          | Optics |
|  |         | Reverse Transcription  | 50°C        | 300            |                  |        |
|  | Stage 1 | Initial PCR Activation | 95.0°C      | 20             |                  | Off    |
|  | Stage 2 | Annealing              | 95.0°C      | 3              | Repeats 45 times | Off    |
|  | Stage 2 | Extension              | 60°C        | 60             |                  | On     |

**Table 4a. PCR amplification conditions in ABI 7500 under standard conditions are as follows:**

|  | A       | B                      | C           | D              | E       | F      |
|--|---------|------------------------|-------------|----------------|---------|--------|
|  | Stage   | Function               | Temperature | Duration (sec) | Repeats | Optics |
|  |         | Reverse Transcription  | 50°C        | 300            |         |        |
|  | Stage 1 | Initial PCR Activation | 95.0°C      | 20             |         | Off    |



|  | A       | B         | C      | D  | E                | F   |
|--|---------|-----------|--------|----|------------------|-----|
|  | Stage 2 | Annealing | 95.0°C | 3  | Repeats 45 times | Off |
|  | Stage 2 | Extension | 60°C   | 60 |                  | On  |

**Table 4b. PCR amplification conditions in ABI 7500 under FAST conditions are as follows:**

Note

**Passive reference dye should be set as ROX**

- 29 Calculate the number of PCR reactions to be performed. Using Table 5 as a guide, determine the quantity of each reagent needed. Include a positive (*C. jejuni* ATCC 33291, ATCC 33292, or ATCC 33560) and negative (*E. coli* ATCC 25922) control and a “No Template” negative control.

For calculations, reference **Table 5**. PCR Worksheet for Reagent Volumes for Master Mix Preparation for *Campylobacter jejuni* PCR.

 Table 5. PCR Worksheet for Reage...

- 30 Combine reagents except for the DNA template in a 1.5ml or 2ml microcentrifuge tube and mix well.
- 31 Aliquot  18 µL into each PCR tube or reaction well if using a plate.
- 32 Add samples and positive (*C. jejuni* ATCC 33291, ATCC 33292, or ATCC 33560 DNA) and negative (*E. coli* ATCC 25922 DNA) controls along with PCR inhibition controls
-  VetMAX Xeno Internal Positive Control ND **Thermo Fisher Scientific Catalog #A29762/A29764**
- 33 Place Samples in thermocycler and perform run.
- 34 Following run, analyze the results. Results should be reported with a Ct value and interpretation.



- 35 Samples with a Ct value of 40 or less should be interpreted as positive. Samples with a Ct value >40 or an undetermined Ct value should be considered negative.

|  | A      | B        | C              |
|--|--------|----------|----------------|
|  | Target | CT Value | Interpretation |
|  | gyrA   | <=40     | Detected       |
|  |        | >40      | Not detected   |

**Table 6.** Interpretation of Ct Value

| A          | B                                               |
|------------|-------------------------------------------------|
| Oligo Name | 5'-3' Sequence                                  |
| gyrA-F     | AAGATACGGTCGATTTTGTTCCTA                        |
| gyrA-R     | CTACAGCTATACCACTTGAACCATTTAATA                  |
| gyrA-P     | 5'-[FAM]TGATGGTTCAGAAAGCGAACCTGATGTTTT[BHQ1] 3' |

**Table 7.** Primers and Probe

## Protocol references

Life Technologies Corporation. 2017. MagMAX™ CORE Mechanical Lysis Module USER GUIDE. Publication Number MAN0015945 Revision A.0. Austin, TX. <https://www.thermofisher.com/document-connect/document-connect.html?url=https://assets.thermofisher.com/TFS-Assets%2FSLG%2Fmanuals%2FMAN0015945-MagMAXCORE-Mech-Lysis-Module-UG.pdf>