

## Histological staining of fish gonadal tissue

Jhennifer Gomes Cordeiro<sup>1</sup>, Caio Maximino<sup>1</sup>, Diógenes Henrique de Siqueira-Silva<sup>1</sup>

<sup>1</sup>Universidade Federal do Sul e Sudeste do Pará

Jun 10,  
2019

Working dx.doi.org/10.17504/protocols.io.ygvftw6

Fish behavior and physiology Medicinal Plants Southeastern Pará Research Group



Caio Maximino

Universidade Federal do Sul e Sudeste do Pará



### ABSTRACT

This protocol describes the dissection, fixation, and histological processing of gonadal tissue of small and medium-sized fish. It has been used in our lab to describe the microanatomy of the gonads of Amazonian Characidae, as well as of ornamental fish such as zebrafish.

### MATERIALS

#### NAME ▾

#### CATALOG # ▾

#### VENDOR ▾

Ethanol

scalpel blades

20 mg Eugenol

orb104769

biorbyt

Tweezers

Sterile glass slides

Toluidine blue O

TB0962.SIZE.5g

Bio Basic Inc.

Xylene

XC9800.SIZE.1L

Bio Basic Inc.

Embedding base molds

22-363-553

Fisher Scientific

Standard square cover slips 18×18 mm

Surgical scissors (Iris - straight)

SCT-I 528

Braintree Scientific

Hematoxylin

[View](#)

Sakura Finetek

Eosin

[View](#)

Sakura Finetek

Bouins Fixative

[View](#)

Polysciences

25 ml glass bottles

[View](#)

Glycol methacrylato resin

[View](#)

Resin hardener powder

[View](#)

Epoxy glue

[View](#)

Tekbond

Tabletop drying oven

[View](#)

Thermo Fisher Scientific

Rotary microtome

[View](#)

Leica Biosystems

Permount mounting medium

[View](#)


Fisher Scientific

## SAFETY WARNINGS

Some of the reagents used, including Bouin's Fixative, may cause serious skin and eye damage, and may cause cancer if inhaled. Wear eye protection, face protection, protective clothing, and protective gloves when manipulating these reagents.

### Collection and fixation of gonads

- 1 Anesthetize and euthanize the animals
- 2 Perform a ventro-longitudinal incision to expose the gonads, Bouin's fixative must be immediately applied on the gonads aiming the prefixation. Then, gonads are removed and immediately placed in the fixative

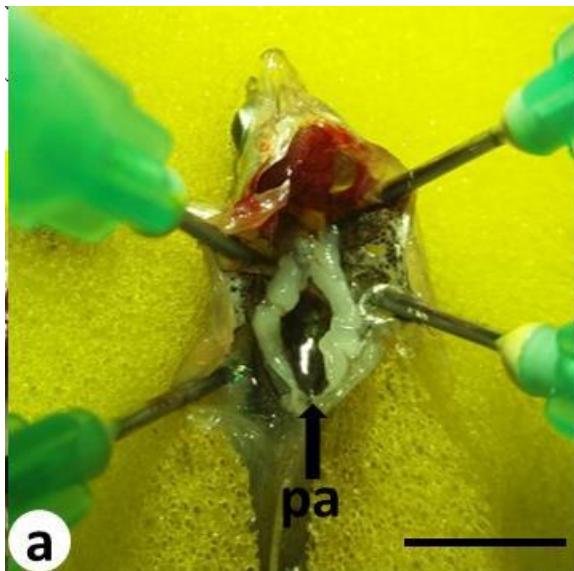

- 3 Leave gonads overnight in Bouin's fixative for 24h ⏱ **24:00:00**
- 4 After 24h, the fixative solution must be replaced with 70% ethanol, that should be renewed every day until there is no trace of fixative

### Dehydration

- 5 Before dehydration the tissue must be sectioned in the desired shape and fit into the molds
- 6 Samples need to be placed in glass vials, to which the solutions will be added



7 Transfer the samples to three fast changes of 70%, 80%, and 90% ethanol **00:03:00** for each change

8 Leave the samples at 95% ethanol for 40 min. **00:40:00**

#### Inclusion

9 Include samples in 1:1 solution of 95% alcohol and infiltration resin (Methacrylate glycol) for 40 minutes **00:40:00**

10 Leave samples in infiltration resin (methacrylate glycol) "overnight" (or at least 12 hours) **12:00:00**

11 Prepare the inclusion resin solution (Methacrylate glycol) with the hardener (for each 1.0 mL of resin use 0.07 mL of hardener)

12 Embedding the material: place the material on the mold and place the resin + hardener solution on it by completing the mold

13 After the resin polymerizes, leave embedded samples in a drying stove for 24 h **24:00:00**

14 Unblock and glue on microtome support stick with epoxy resin glue



- 15 After the glue has completely dried, the sample should be taken to the microtome
- 16 The sample should be sectioned to the thickness of 3 to 5  $\mu\text{m}$
- 17 The cuts should be deposited on properly identified slides
- 18 The slides should be kept in the drying stove for 1 day and then followed for staining procedures. If staining with H&E, follow steps 19-27. If staining with toluidine blue, follow steps 28-33 **⌚24:00:00**




#### Staining: Hematoxylin and eosin

- 19 Moisturize cuts in running tap water for 10 minutes ⌚ 00:10:00
- 20 Leave in Harri's hematoxylin solution for 15 minutes ⌚ 00:15:00
- 21 Wash in running tap water for 10 minutes ⌚ 00:10:00
- 22 Dehydrate in 90% alcohol with fast change (if eosin is alcoholic) ⌚ 00:03:00
- 23 Leave in eosin solution for 15 minutes ⌚ 00:15:00
- 24 Rinse well under running water
- 25 Transfer the samples to three fast changes of 70%, 80%, and 95% ethanol, followed by fast passes in absolute ethanol, and two changes of xylene (3 min for each station) for clearing before coverslipping ⌚ 00:03:00 for each change
- 26 Transfer the sample to a third change of xylene ⌚ 00:03:00
- 27 Mount and coverslip in Permount (xylene-based mounting medium)


28

Check slides on the microscope and take photographs as needed.

Example slides:

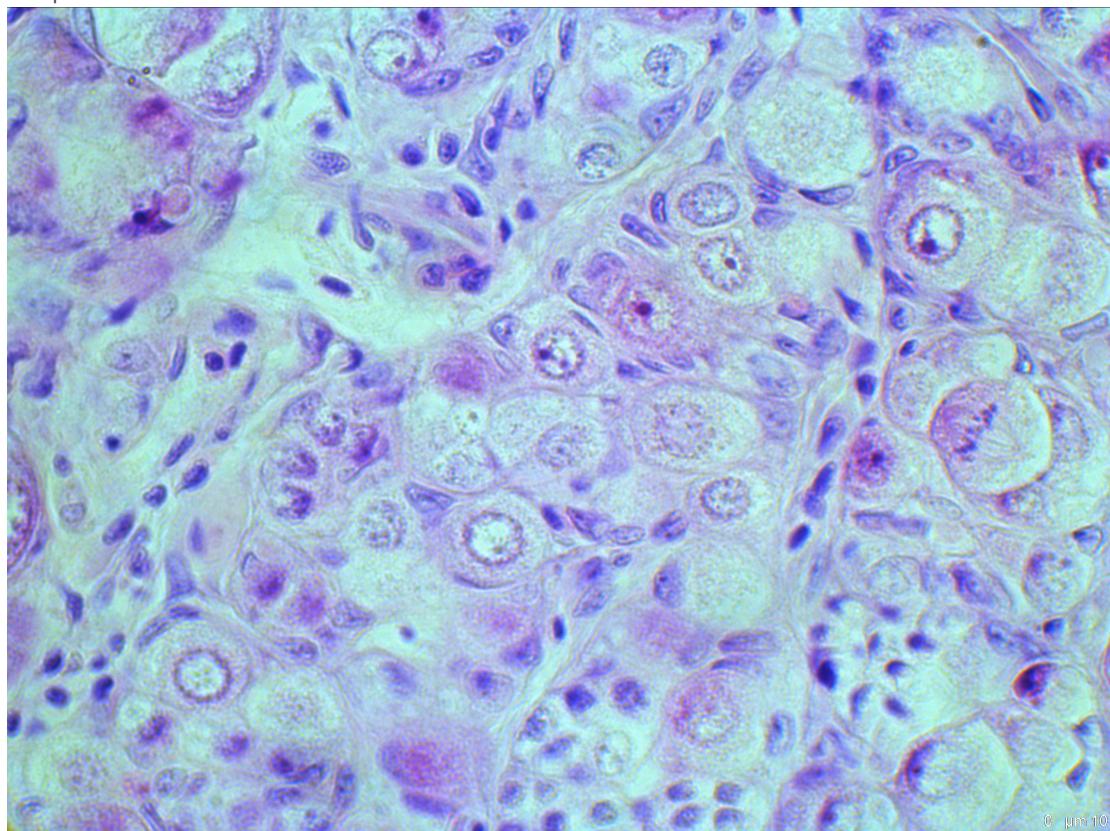

Immature female *Astyanax bimaculatus* gonads

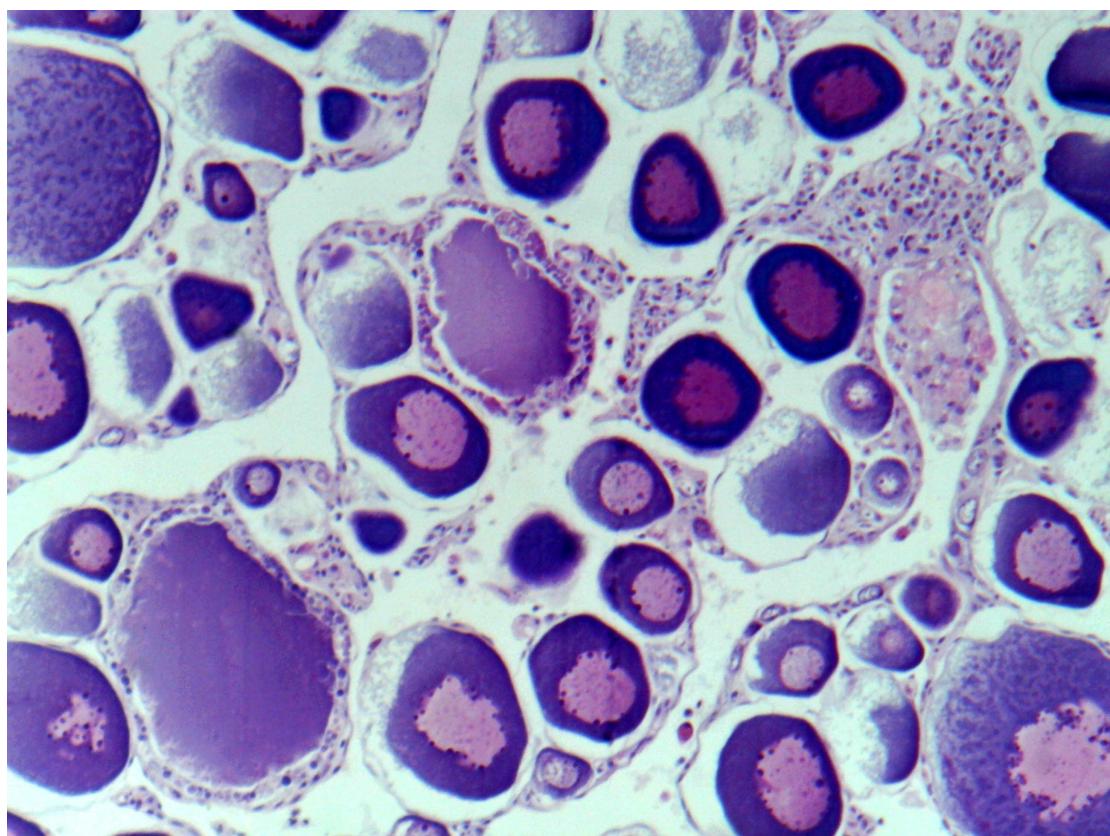

Initial gonadal development in male *Astyanax bimaculatus*

Staining: Toluidine blue

- 29 Moisturize cuts in running water for 10 minutes **⌚00:10:00**
- 30 Leave in toluidine blue solution for 15 minutes **⌚00:15:00**
- 31 Wash in running water for 10 minutes **⌚00:10:00**
- 32 Transfer the samples to three fast passes of 70%, 80%, and 95% ethanol, followed by fast passes in absolute ethanol, and two changes of xylene (3 min for each station) for clearing before coverslipping **⌚00:03:00 for each change**
- 33 Transfer the sample to a third change of xylene **⌚00:03:00**
- 34 Mount and coverslip in Permount (xylene-based mounting medium)

 This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited