



Oct 13, 2022

Movie_Timepoint_CopytoClipboard_Tool

Condon ND¹¹University of Queensland

1 Works for me

 Sharedx.doi.org/10.17504/protocols.io.n92ldp9e9l5b/v1 Condon ND

ABSTRACT

This protocol describes how to use the Movie-Timepoint_CopytoClipboard_Tool.ijm program made by Nicholas Condon at the Institute for Molecular Bioscience's Microscopy core facility at The University of Queensland, Brisbane, Australia.

This tool allows the user to quickly generate figures in programs external to FIJI for time-series data. It allows the user to choose which channels of the open image they would like incorporated, as well as an optional Merged panel. Users then navigate to each frame and a series of popups will guide them through pasting the image frames into their program of choice. This tool records the frame numbers chosen and the corresponding time value (if known from the image metadata file). There is also the option to save out the results into a user-chosen directory for better record keeping. A log with the frame number and timepoint information is also saved.

DOI

dx.doi.org/10.17504/protocols.io.n92ldp9e9l5b/v1

PROTOCOL CITATION

Condon ND 2022. Movie_Timepoint_CopytoClipboard_Tool. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.n92ldp9e9l5b/v1>



KEYWORDS

FIJI Macro, Figure making, automation, script

LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Oct 13, 2022

LAST MODIFIED

Oct 13, 2022

PROTOCOL INTEGER ID



71283

Running the Script


1 Launch FIJI




2 Open your time-series dataset into FIJI

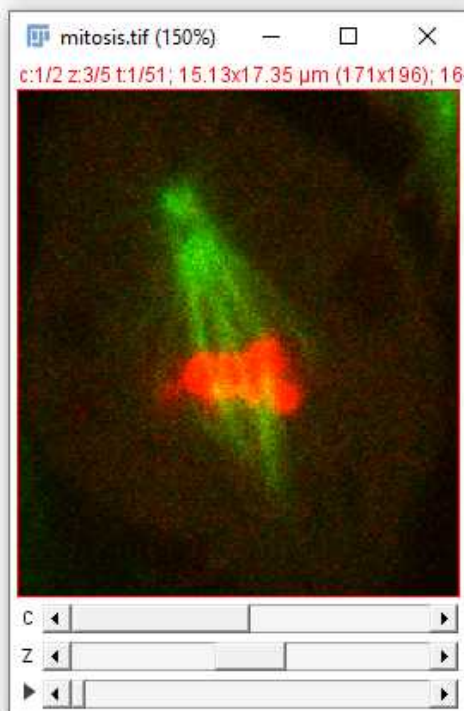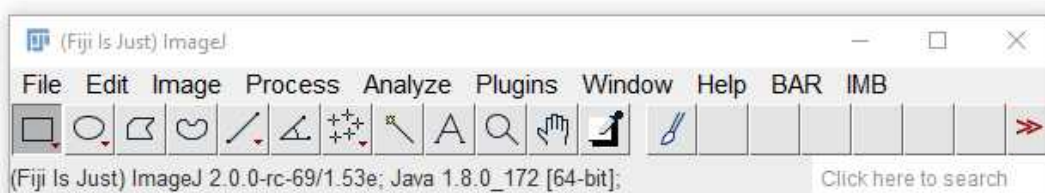


- 3 Open the **Movie Timepoint CopytoClipboard Tool** script
[Plugins > Macros > Run...]

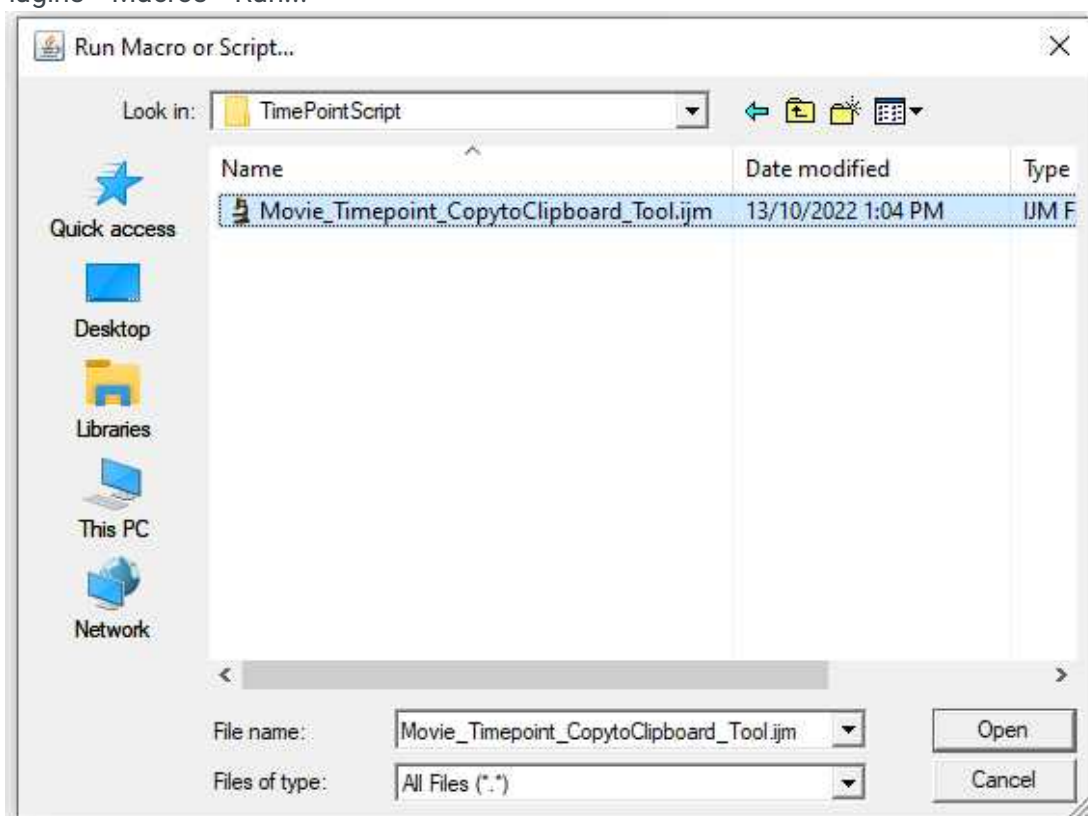


- 4 Read the instructions & acknowledgement window

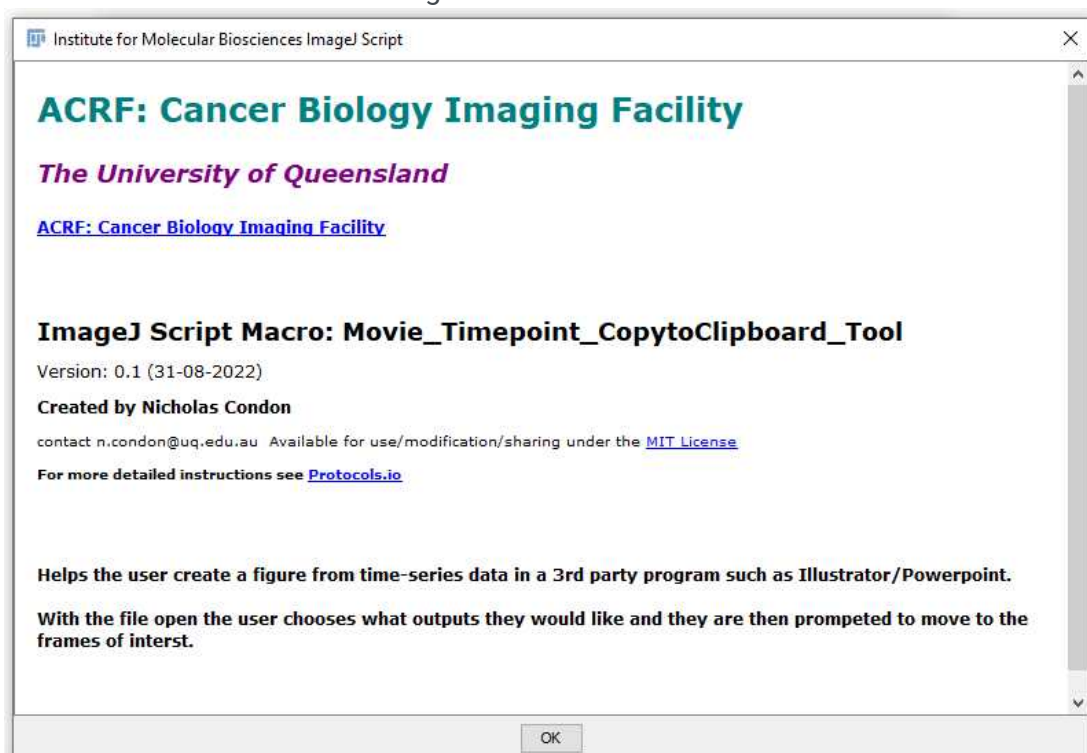

- 5 Ensure a time-series file is open, and is at "Top Level" (i.e. select the window) before clicking ok.

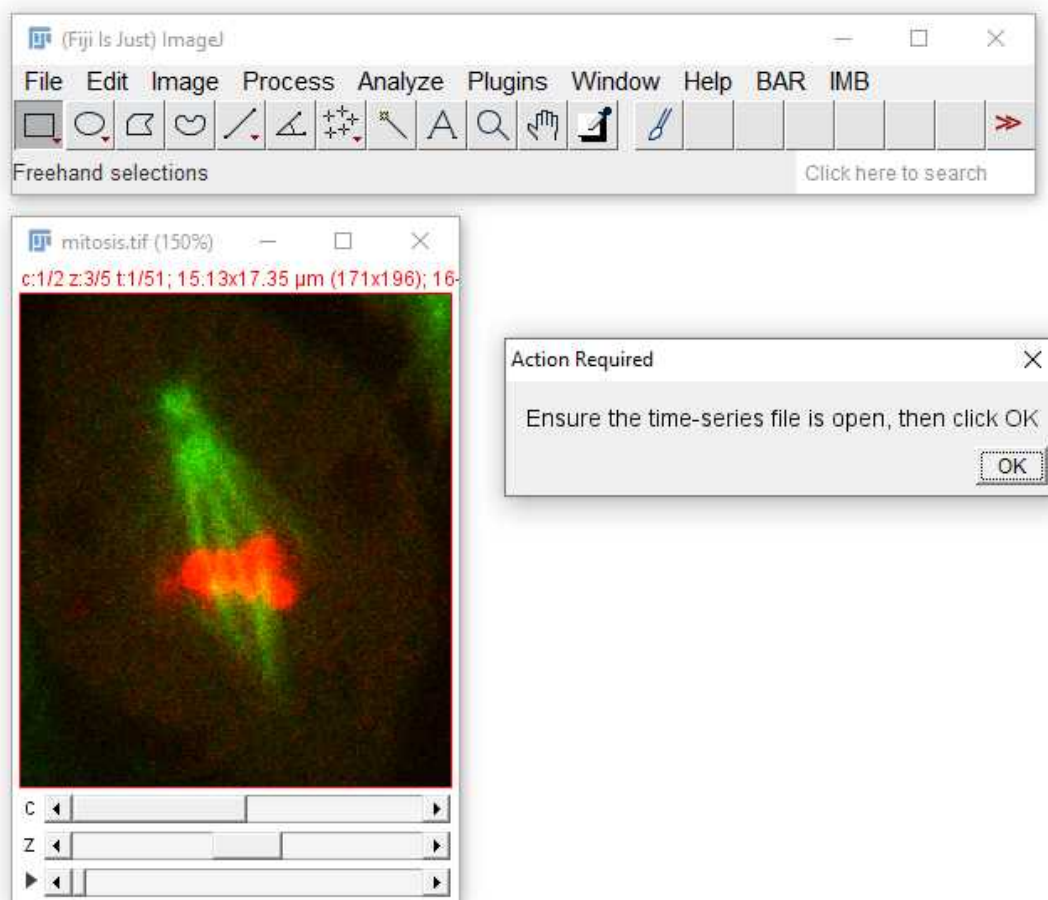

- 6 Select which Channel(s) you would like to include in the run. There is also the option to include a merged channel (ticked by default).
Images added to the clipboard can also be saved in a folder for later recalling if needed.

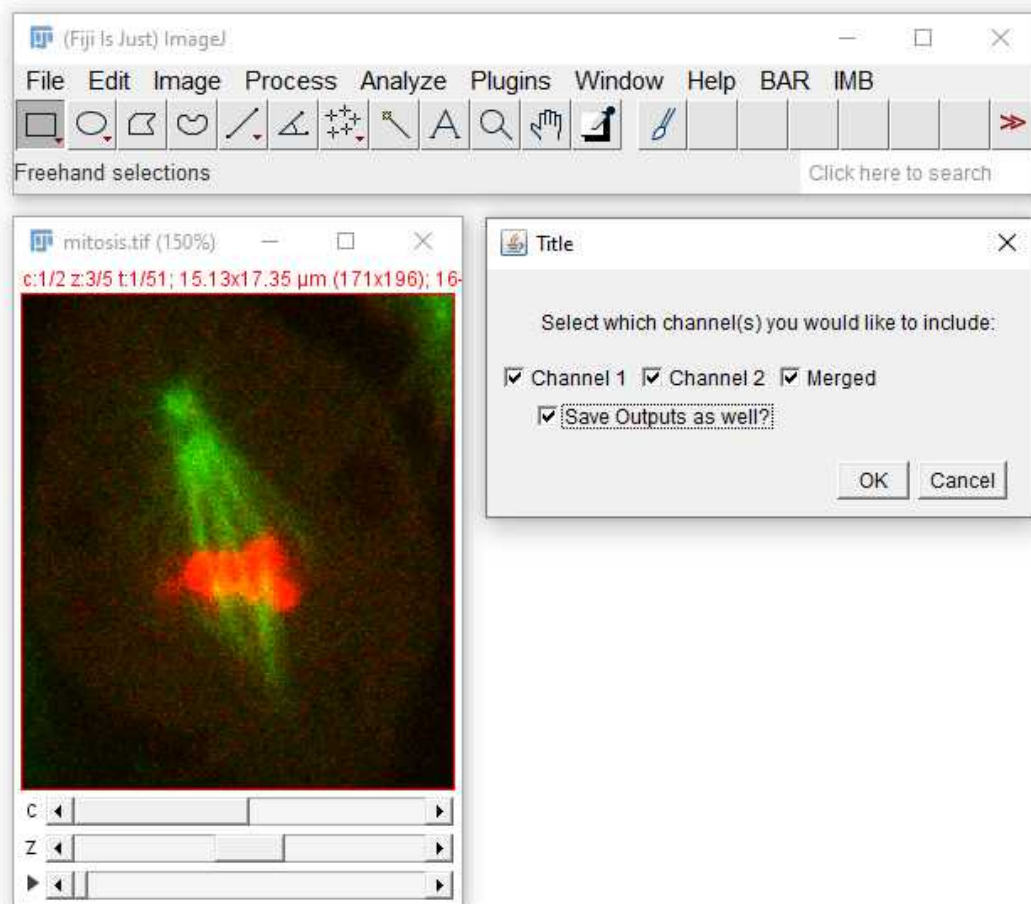

Step 6 includes a Step case.


Save Outputs as well? (Enabled)

Save Outputs as well? (disabled)


step case

Save Outputs as well? (Enabled)


7 Choose an Output Directory location.

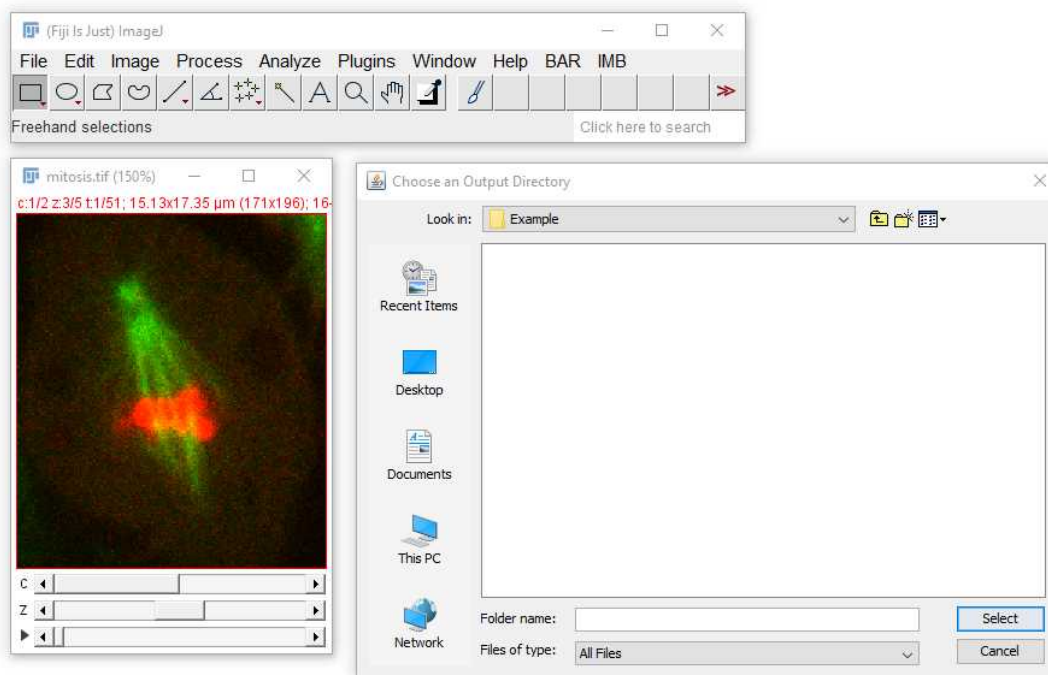

8 Move to the desired timepoint (in the image window), and then click OK.

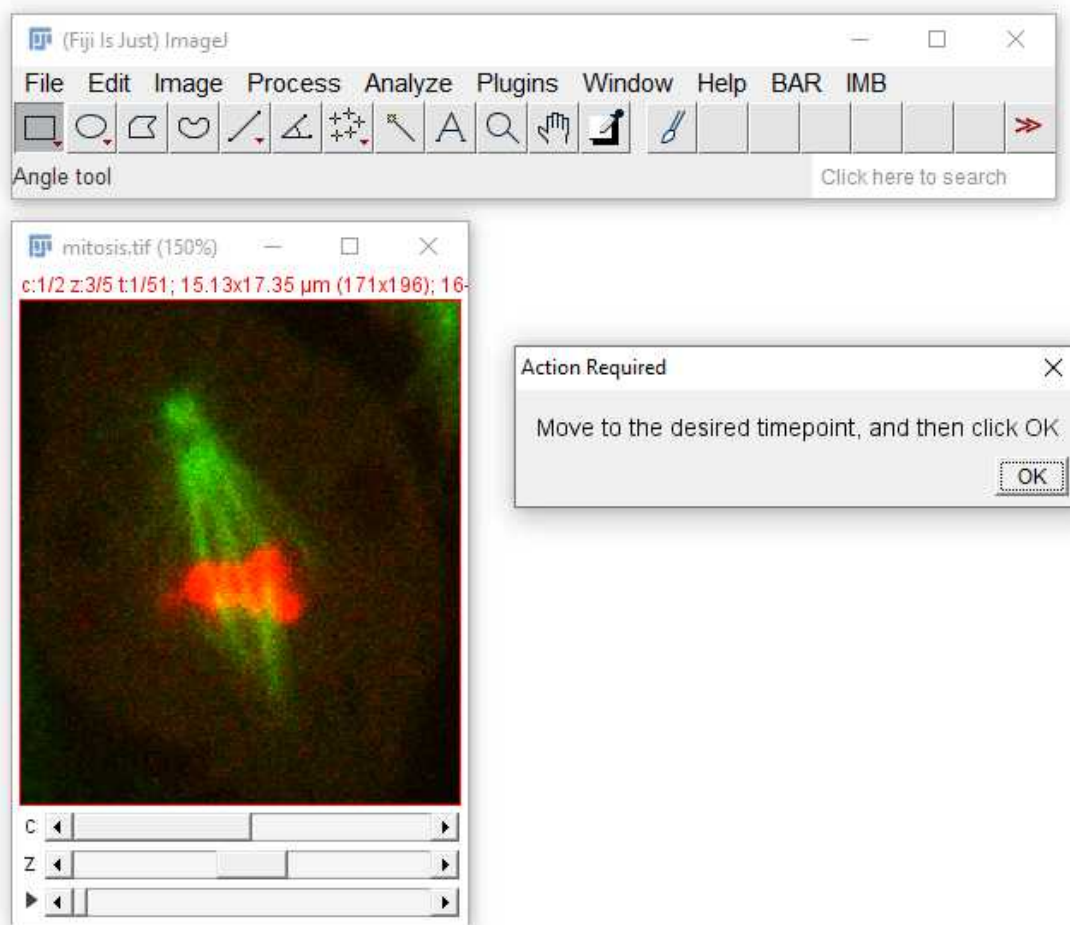

- 9 Merged Channel now copied, return when pasted and click OK.

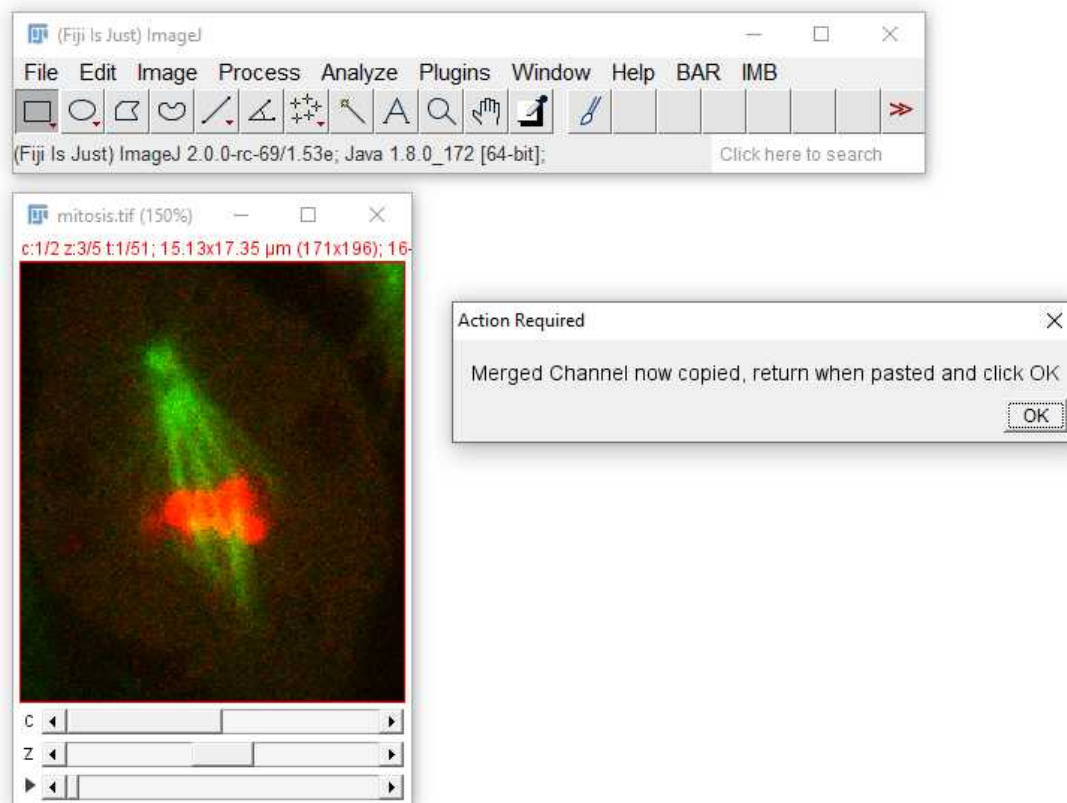


Navigate to your program of choice (e.g. Microsoft PowerPoint or Adobe Illustrator/In Design) and paste the image from the clipboard. Then click OK.


- 10 Channel 1 now copied, return when pasted and click OK.

Still in your program of choice, you can paste in the image file without going back to FIJI.


Then click OK.

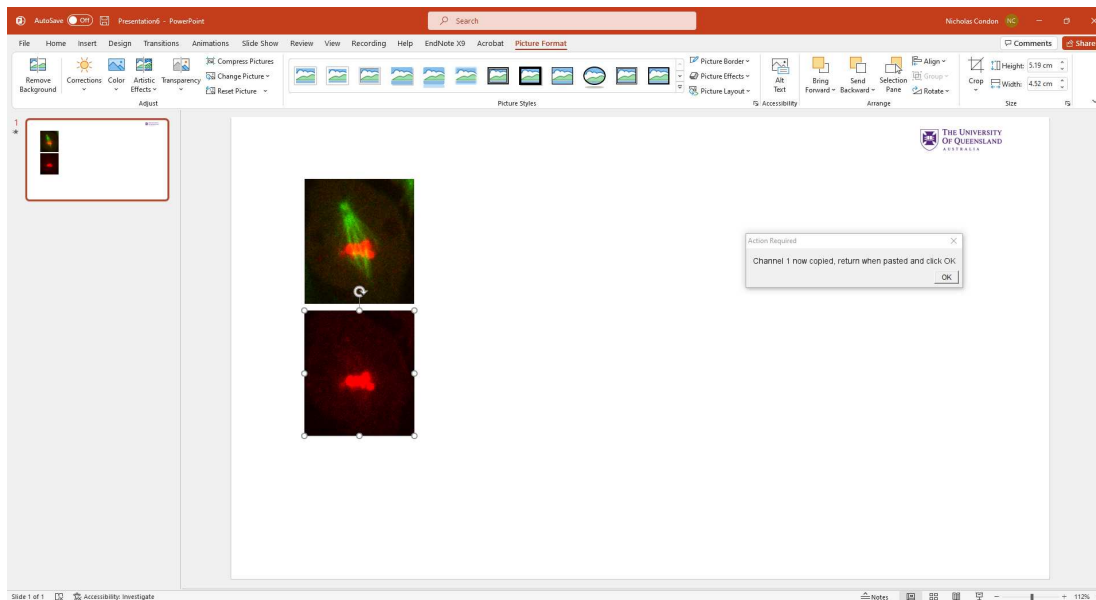


11 Repeat for however many channels were selected in the window in Step 6.

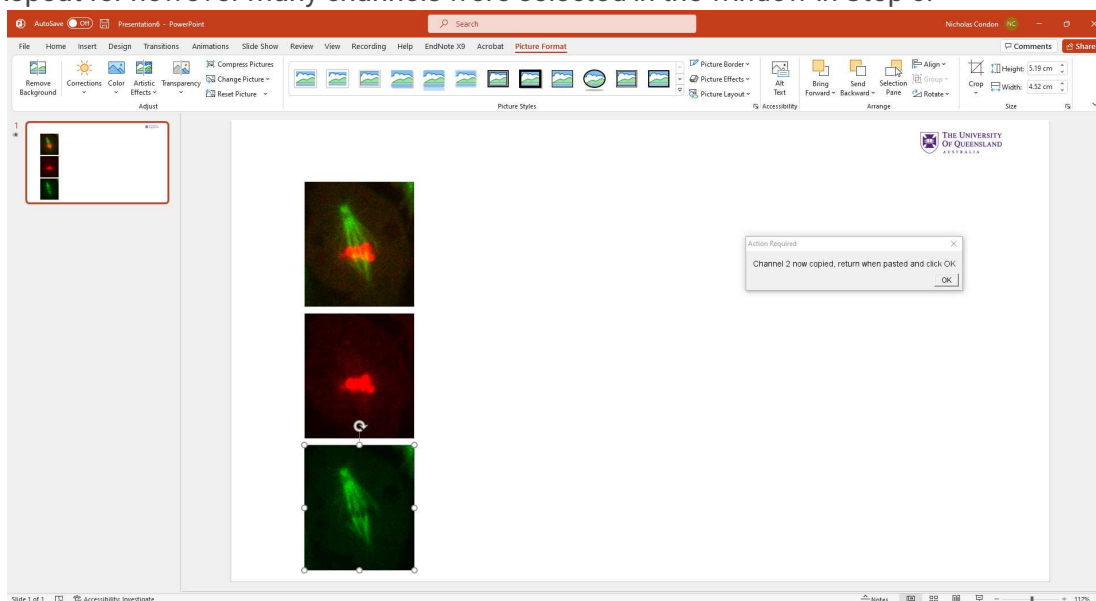


12 Do you have more timepoints?
If you require more timepoints to be selected, leave the Checkbox selected and click OK.

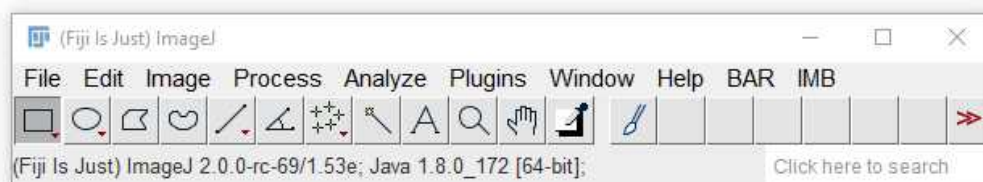


13 Move to the desired timepoint (using the sliders in the image window), and then click OK.

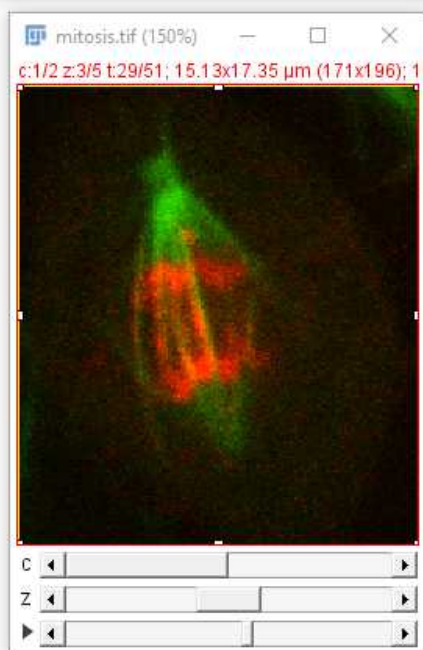


14 Repeat Steps 9-11 by pasting the images into your program of choice.

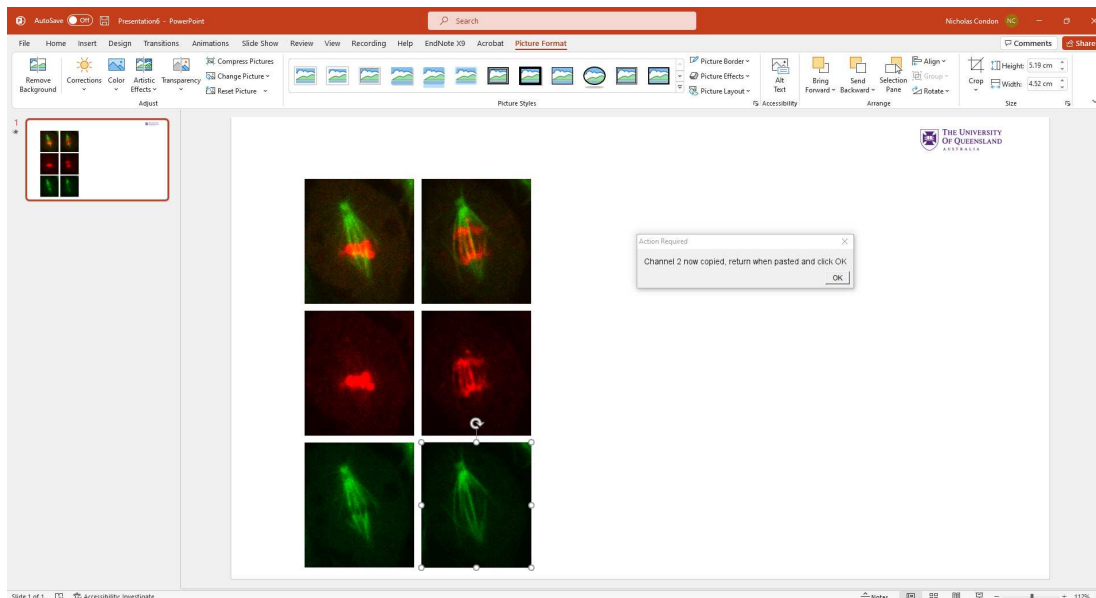

- 15 Repeat Steps 8-14 for as many timepoints as you desire.
When finished uncheck the box on the Do you have more timepoints? window.

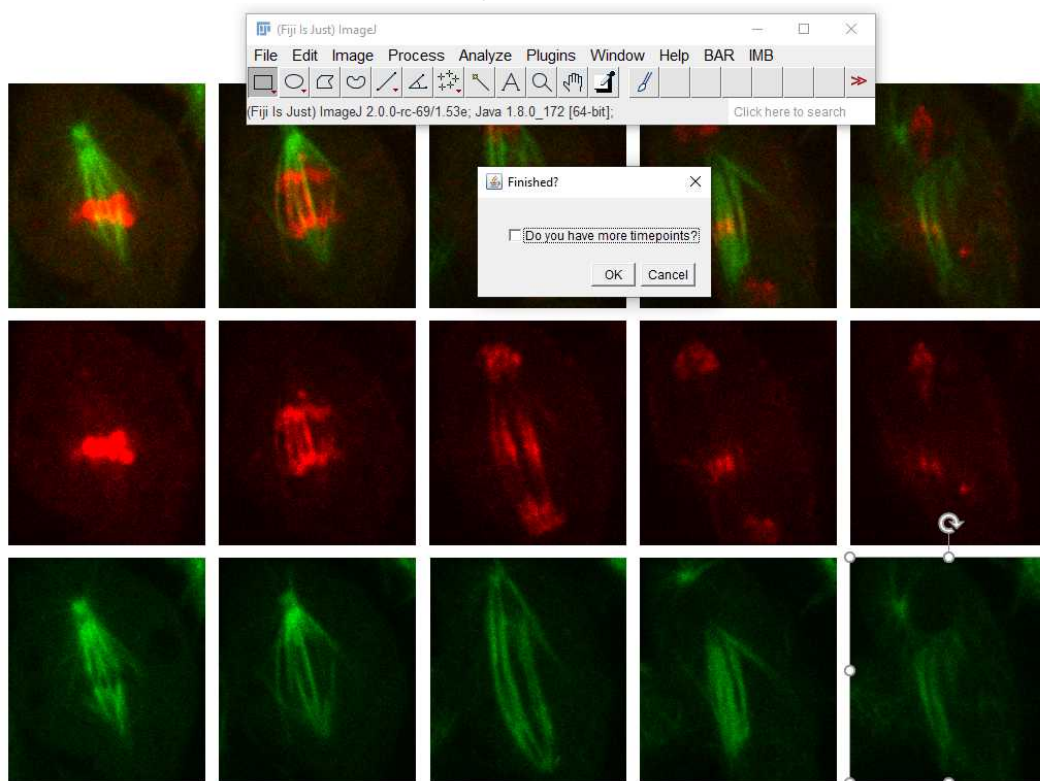

- 16 Batch Completed. Put down that coffee! Your analysis is finished, dialog will appear to inform you of the safe completion of the program. Note the log file also records which frames and the corresponding time value for that frame to help annotating the figure.

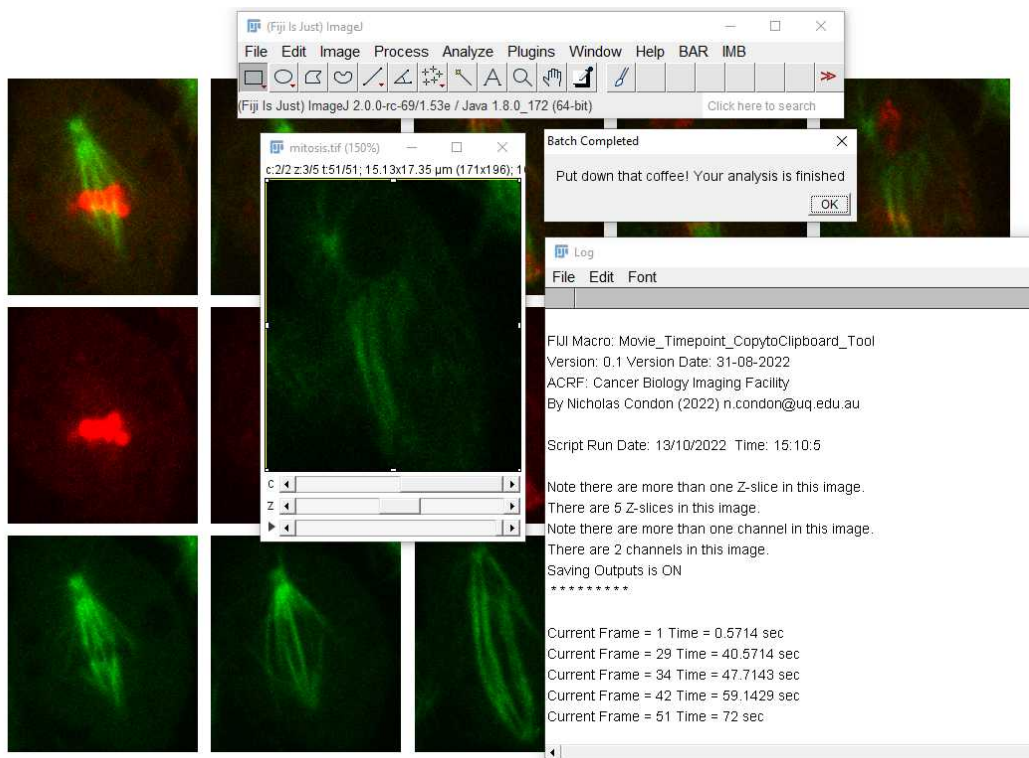

- 17 Navigate to the output directory location to access your saved files.
Note the log is also saved with the frame number and timepoint information.

