



NOV 24, 2022

WORKS FOR ME

1

# Cell treatments for time-lapse experiments


COMMENTS 0

DOI

[dx.doi.org/10.17504/protocols.io.j8nlkw8bxi5r/v1](https://dx.doi.org/10.17504/protocols.io.j8nlkw8bxi5r/v1)

Electra Brunialti<sup>1</sup>,  
Alessandro Maria Villa<sup>1</sup>, Paolo Ciana<sup>1</sup>

<sup>1</sup>Department of Health Sciences, University of Milan, Milan, Italy

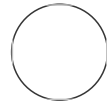

Electra Brunialti

## ABSTRACT

Cell treatments for time-lapse experiments.

DOI

[dx.doi.org/10.17504/protocols.io.j8nlkw8bxi5r/v1](https://dx.doi.org/10.17504/protocols.io.j8nlkw8bxi5r/v1)

## PROTOCOL CITATION

Electra Brunialti, Alessandro Maria Villa, Paolo Ciana 2022. Cell treatments for time-lapse experiments. **protocols.io**  
<https://dx.doi.org/10.17504/protocols.io.j8nlkw8bxi5r/v1>



## FUNDERS ACKNOWLEDGEMENT

The work was supported by the EU Joint Programme - Neurodegenerative Disease Research (JPND) project GBA-PaCTS, 01ED2005B

## LICENSE

\_\_\_\_\_ This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

## CREATED

Nov 23, 2022

## LAST MODIFIED

Nov 24, 2022

## PROTOCOL INTEGER ID

73173

- 1 For the timelapse experiments, cell cultures were treated for 6 h with vehicle (water) or lipopolysaccharide (LPS) at the final concentration of 10 µg/ml LPS O111:B4, (Cat. L2630, Merk) starting from a 0,5 mg/ml stock solution, and then subjected to timelapse microscopy for 2h. For the GCase inhibition, cells were treated for 48 h with 200 µM CBE (Cat. 234599, Merk) starting from a 10 mM stock solution, or vehicle (water), and then subjected to timelapse microscopy for 2h.