

JAN 09, 2024

Experiment Data Depot (EDD) Data Import

Forked from [Experiment Data Depot \(EDD\) Study Creation](#)

Jennifer Gin¹, Christopher J Petzold¹

¹Lawrence Berkeley National Laboratory

LBNL omics

Agile BioFoundry

[1 more workspace](#) ↓



Christopher J Petzold

Lawrence Berkeley National Laboratory

OPEN ACCESS



DOI:
dx.doi.org/10.17504/protocols.io.bp2l6bd15gqe/v1

External link: <https://public-edd.agilebiofoundry.org/>

Protocol Citation: Jennifer Gin, Christopher J Petzold 2024. Experiment Data Depot (EDD) Data Import . **protocols.io** <https://dx.doi.org/10.17504/protocols.io.bp2l6bd15gqe/v1>

MANUSCRIPT CITATION: Morrell, William C., et al. "The experiment data depot: a web-based software tool for biological experimental data storage, sharing, and visualization." ACS synthetic biology 6.12 (2017): 2248-2259.

DISCLAIMER

DISCLAIMER – FOR INFORMATIONAL PURPOSES ONLY; USE AT YOUR OWN RISK

The protocol content here is for informational purposes only and does not constitute legal, medical, clinical, or safety advice, or otherwise; content added to [protocols.io](#) is not peer reviewed and may not have undergone a formal approval of any kind. Information presented in this protocol should not substitute for independent professional judgment, advice, diagnosis, or treatment. Any action you take or refrain from taking using or relying upon the information presented here is strictly at your own risk. You agree that neither the Company nor any of the authors, contributors, administrators, or anyone else associated with [protocols.io](#), can be held responsible for your use of the information contained in or linked to this protocol or any of our Sites/Apps and Services.

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Jan 11, 2021

Last Modified: Jan 09, 2024

PROTOCOL integer ID:
46196

Keywords: EDD, Metadata, Synthetic Biology, Experiment Data Depot, Metabolic Engineering, REST API

ABSTRACT


This protocol details how to import data into an EDD Study. Data can be shared with collaborators via exported CSV files or through the REST API.

For more information about the EDD visit: [Experiment Data Depot Google Site](#).

Citation:
Morrell, William C., et al. "The experiment data depot: a web-based software tool for biological experimental data storage, sharing, and visualization." ACS synthetic biology 6.12 (2017): 2248-2259. ([link to publication](#))

Example Links:
<https://public-edd.agilebiofoundry.org/>
<https://public-edd.agilebiofoundry.org/s/example-data-quality-study-2/overview/>

MATERIALS

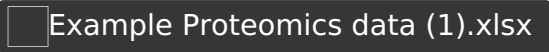
Example proteomics data import file:


BEFORE START INSTRUCTIONS

The following instructions show you how to import data for proteomics and metabolomics.

Create an EDD study ([EDD Study Creation protocol](#)) before you start the data import process.

1 Click the **Import Data** button on the top right.

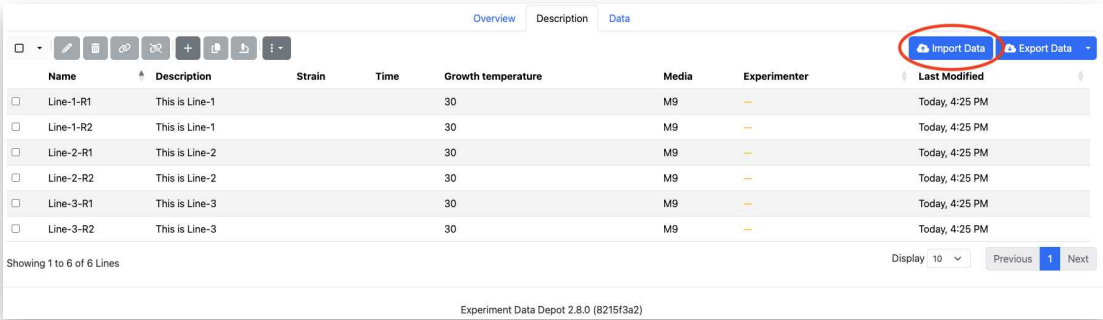


Import Data

2 Go through the set of data import steps:

Identify what protocol was used. Click **Shotgun (Discovery) Proteomics**.
Select the layout of your data. Click **Generic**.

Click **Start Import** on the bottom right.

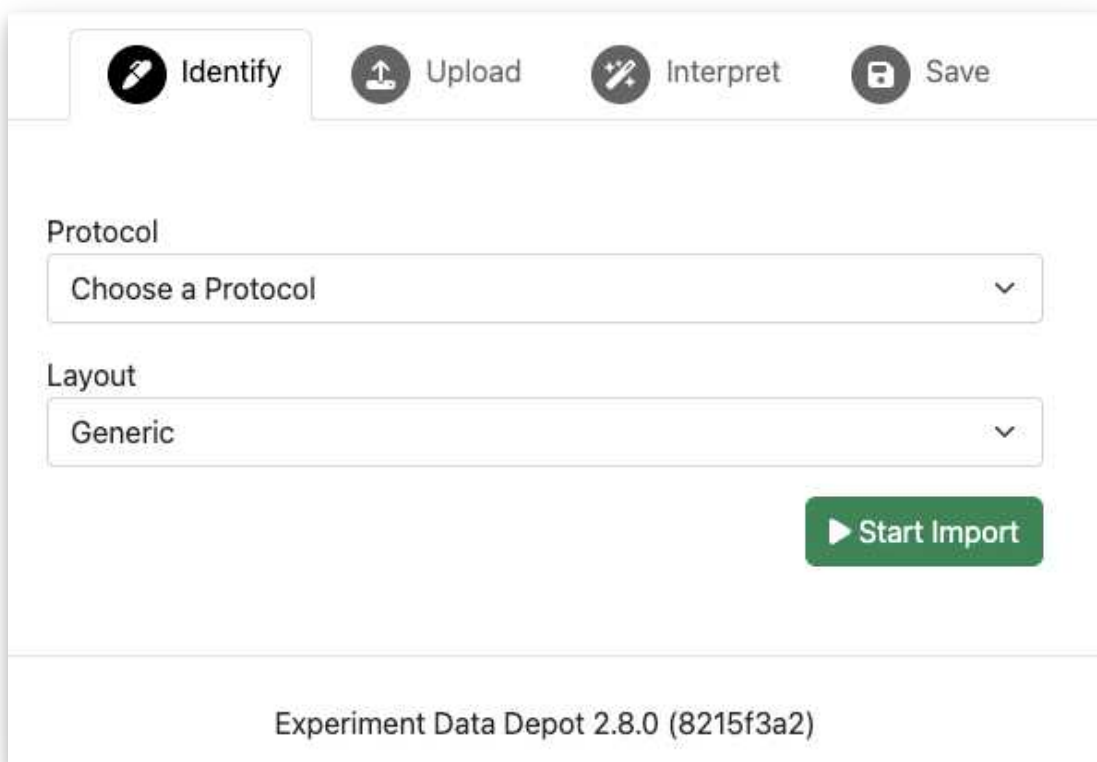

The screenshot shows the EDD interface with four tabs: Identify, Upload, Interpret, and Save. The 'Identify' tab is active. Below the tabs, there are two dropdown menus. The first is labeled 'Protocol' and has 'Choose a Protocol' selected. The second is labeled 'Layout' and has 'Generic' selected. A green 'Start Import' button is located to the right of the 'Layout' dropdown. At the bottom of the interface, the text 'Experiment Data Depot 2.8.0 (8215f3a2)' is displayed.

Note

Repeat for each data type you want to import into the EDD.

3 Example data layouts:

Here, we have measured the amount of protein "XYLA_ECODH", metabolite "D-Limonene", and extracellular "D-Glucose" for sample "181-aceF" in triplicates at time "12" hours.

A	B	C	D	E
Line Name	Measurement Type	Time	Value	Units
181-aceF-R1	sp B1X8I1 XYLA_ECODH	12	2795	counts
181-aceF-R2	sp B1X8I1 XYLA_ECODH	12	2611	counts
181-aceF-R3	sp B1X8I1 XYLA_ECODH	12	2155	counts

Proteomics Example Data File

 Example Proteomics data.xlsx

A	B	C	D	E
Line Name	Measurement Type	Time	Value	Units
181-aceF-R1	CID:440917 D-Limonene	12	0.079585069	mM
181-aceF-R2	CID:440917 D-Limonene	12	3.712638406	mM
181-aceF-R3	CID:440917 D-Limonene	12	0.416450273	mM

Metabolomics Example Data File


 Example Metabolomics data.xlsx

A	B	C	D	E
Line Name	Measurement Type	Time	Value	Units
181-aceF-R1	CID:5793 D-Glucose	12	22.07066865	g/L
181-aceF-R2	CID:5793 D-Glucose	12	21.84441193	g/L
181-aceF-R3	CID:5793 D-Glucose	12	21.45756405	g/L

HPLC Example Data File


 Example HPLC aceF.xlsx

Drag-and-drop to upload data file into the dropzone box.


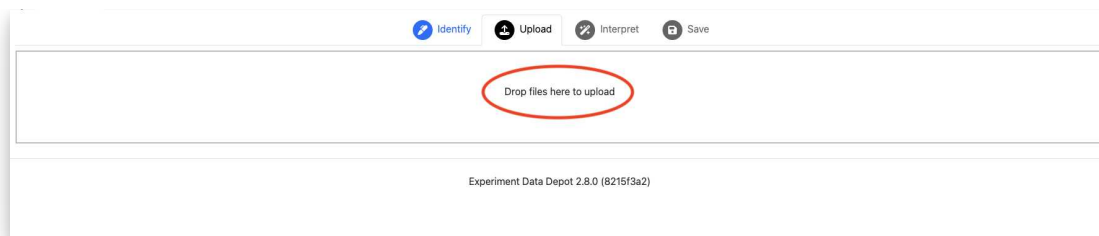


Drag data file (.xlsx) to the dropzone

- 4 After the EDD interprets the file it will calculate the number of records to save. Check that this is correct, then click the **Save** button.

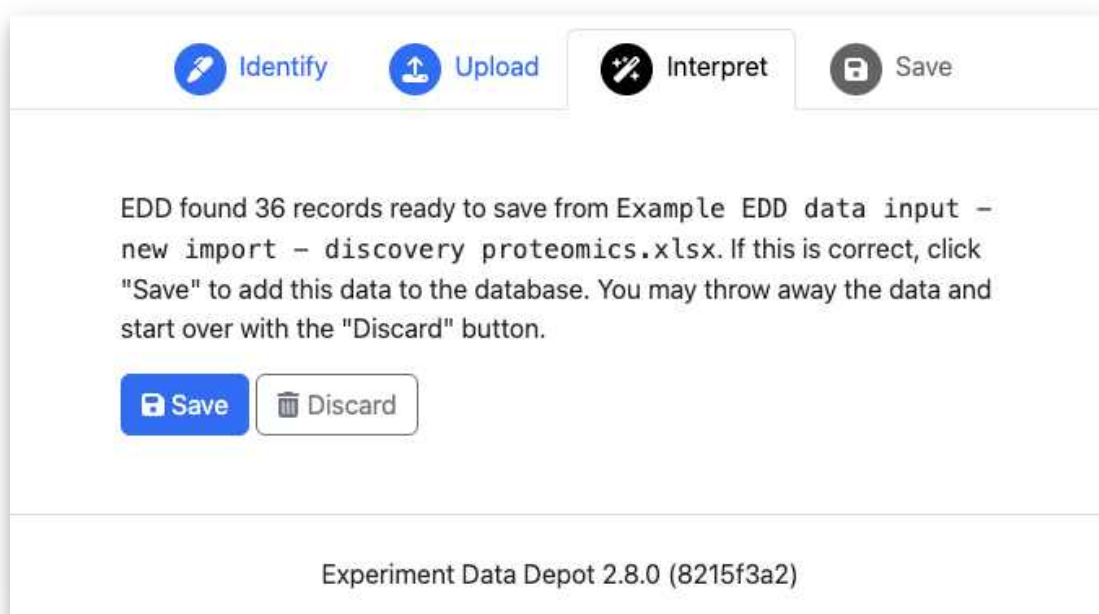

- 5 Once the data is uploaded, the EDD will bring you to the **Data** page in the **Line graph** view. You can toggle to the **Bar Graphs**.

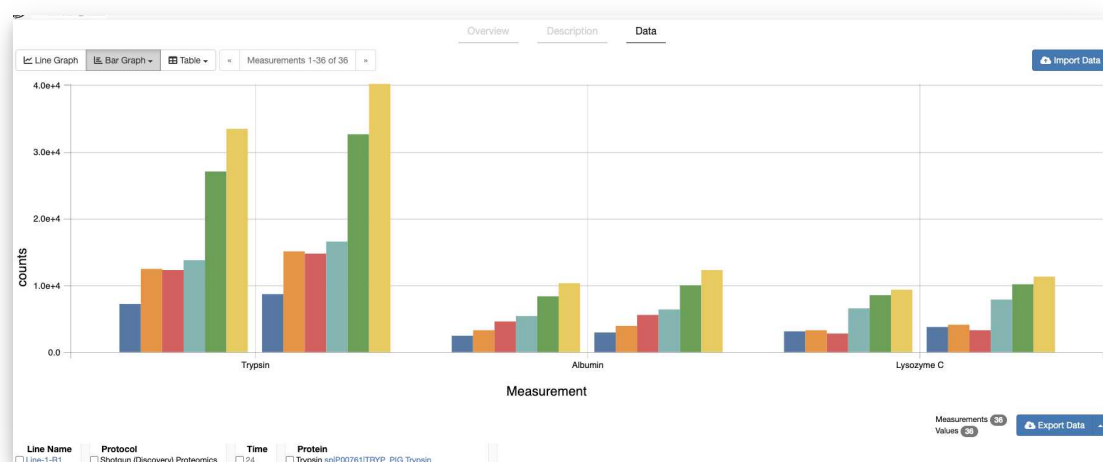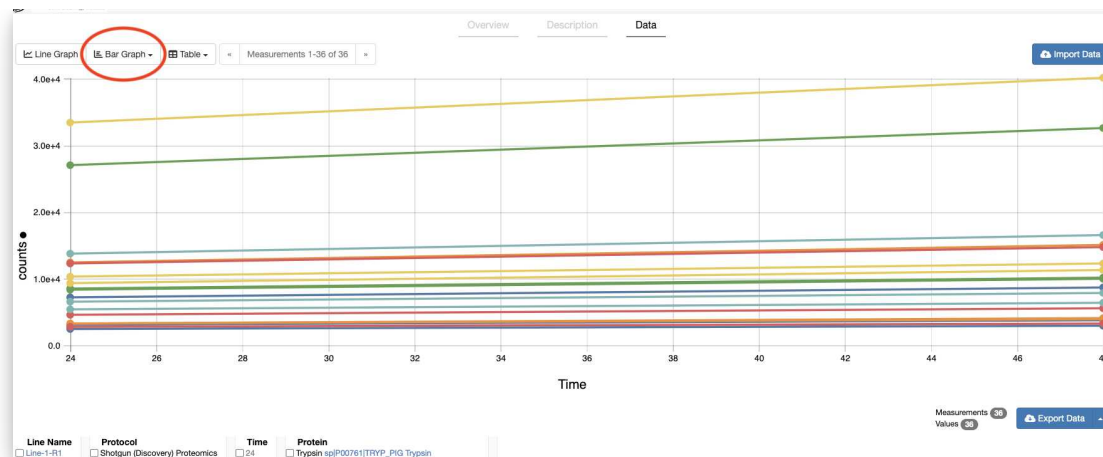

6 Repeat the Data Import process for other data types.