

JUN 09, 2023

# Pole Test to assess motor coordination in parkinsonian mice

## Zachary

natalia.lopezgonzalezdelrey<sup>1</sup>, Gaertner<sup>1</sup>

<sup>1</sup>Northwestern University, Feinberg School of Medicine, Neurology Department



natalia.lopezgonzalezdelrey

# OPEN ACCESS

#### DOI:

dx.doi.org/10.17504/protocol s.io.dm6gp3ed8vzp/v1

**Protocol Citation:** natalia.lo pezgonzalezdelrey, Zachary Gaertner 2023. Pole Test to assess motor coordination in parkinsonian mice.

#### protocols.io

https://dx.doi.org/10.17504/protocols.io.dm6gp3ed8vzp/v1

License: This is an open access protocol distributed under the terms of the Creative Commons
Attribution License, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working We use this protocol and it's working

Created: Jun 09, 2023

Last Modified: Jun 09, 2023

**PROTOCOL** integer ID:

83161

**Keywords:** ASAPCRN, behavior protocol, parkinsonian behavior, pole test

#### **ABSTRACT**

The pole test evaluates the motor coordination ability of a mouse by having them grasp a pole and then descend to its home cage.

#### **MATERIALS**

16in tall, 1/2in diameter wooden pole coated in tape to allow gripping

## **Pole Test**

1d


1 Total duration: 1 day

Three consecutive trials per animals are run.

### **Protocol**

1h

2 Place the animal cage in the testing room to allow the animals to acclimate. 01:00:00

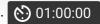

3 Mice are placed with their head oriented upward at the top of a 16in tall, 1/2in diameter wooden pole coated in tape to allow gripping. The animals normally orient themselves downward and then descend the length of the pole in order to return to their home cage.

5m

4 Return the mouse to its cage.

# **Analysis**

5 The time required for the animals to orient themselves facing in a downward direction (time to T turn) and to descend to the base of the pole (total time) is recorded.

Mean times per mouse are compared across groups.