



OCT 12, 2023


🌐 Glucose Tolerance Test

Sabina

Marciano¹, Roberta Marongiu¹

¹Weill Cornell Medicine

ASAP Collaborative Research Network

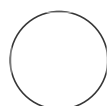

Dylan J Barbera

Weill Cornell Medicine

OPEN  ACCESS



ABSTRACT

The glucose tolerance test is performed to determine how quickly the glucose is cleared from the blood. It is used to test for [diabetes](#) or [insulin resistance](#).

DOI:

dx.doi.org/10.17504/protocols.io.261geddodv47/v1

Protocol Citation: Sabina Marciano, Roberta Marongiu 2023. Glucose Tolerance Test. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.261geddodv47/v1>


License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Aug 14, 2023

Last Modified: Oct 12, 2023

Keywords: ASAPCRN

- 1 Single cage the mice for 1 week in advance
- 2 Weight the mice
- 3 Fast mice for  06:00:00 6h
- 4 After fasting, measure glucose with a glucose meter
- 5 Perform an intraperitoneal injection of 2g / kg body weight of glucose (20% D-glucose stock solution dissolving 2g of glucose in 10ml saline and give 10ul per gram of body weight).
- 6 Collect blood samples for glucose measurement at 15, 30, 45, 60, 90 and 120 min post glucose injection