

Version 2

Jul 15, 2020

Fitting enzyme kinetics data with Solver V.2

Nithesh Chandrasekharan¹, Chris Berndsen¹¹James Madison University

1

Works for me

This protocol is published without a DOI.

Chris Berndsen
James Madison University

ABSTRACT

Reaction kinetics are a fundamental component of the biochemical characterization of a biomolecule. The V_{\max} and the K_M values derived from experiments describe how fast a biomolecule catalyzes the reaction and the approximate affinity between biomolecule and substrate. For almost a century, these properties have been used to describe the mechanism and functions of enzyme and other biomolecules. Fitting data to calculate kinetic values like the V_{\max} and the K_M values is a fundamental skill in biochemistry.

This protocol describes fitting data determined at several substrate concentrations to calculate V_{\max} and K_M values using Excel/Google Sheets and the macro with these programs Solver. This protocol is based on the work of Kemmer and Keller.

This protocol is written for students in Biochemistry at James Madison University, however the approach is general enough for any enzyme kinetics data.

THIS PROTOCOL ACCOMPANIES THE FOLLOWING PUBLICATION

Kemmer, G., and Keller, S. (2010) Nonlinear least-squares data fitting in Excel spreadsheets. Nat. Protoc. 5, 267–281.

PROTOCOL CITATION

Nithesh Chandrasekharan, Chris Berndsen 2020. Fitting enzyme kinetics data with Solver. **protocols.io**
<https://protocols.io/view/fitting-enzyme-kinetics-data-with-solver-bik7kczn>

MANUSCRIPT CITATION please remember to cite the following publication along with this protocol

Kemmer, G., and Keller, S. (2010) Nonlinear least-squares data fitting in Excel spreadsheets. Nat. Protoc. 5, 267–281.

KEYWORDS

kinetics, excel, solver, data fitting

LICENSE

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Jul 15, 2020

LAST MODIFIED

Jul 15, 2020

PROTOCOL INTEGER ID

39295

MATERIALS TEXT

Microsoft Excel or Google Sheets

Kinetics data from at least 5 concentrations of substrate

Where to start

- 1 If you have calculated rate data already skip to step 2.3. If you have time course data and need to calculate the rate, start with step 2.


Calculating initial rates

- 2 Calculate the initial rate from the slope of the product formed vs. time plot.



It is important to calculate the rate when the plot of product formed vs. time is linear with time. This ensures that the enzyme is in the **steady-state**. Typically this is in the early part of the reaction before more than 10% of the substrate has been converted to product (see figure). *The rate is the slope of the line.*

In the steady-state, the concentration of enzyme bound complexes is not changing with time and the enzyme is catalyzing the reaction at the fastest possible rate at that concentration of substrate.

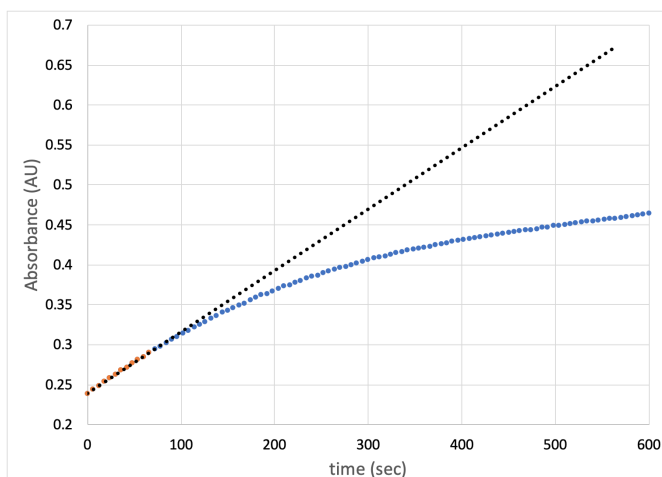

Example of fitting to calculate the initial velocity. Blue points are data. Blank dashed line shows the fitting of the trendline to the portion of the data where it is linear and presumably the reaction is still in the steady-state. The slope of the black dashed line would be the initial velocity or initial rate.

- 2.1 Repeat for each concentration of substrate.
- 2.2 Subtract the rate of the background reaction from each rate.



Background rate is typically from a reaction performed in the absence of enzyme which shows the rate of the reaction in the absence of catalysis.

- 2.3 Record the rate at each concentration of substrate in the table below. Be sure to modify the units to match the units from the experiment. Add rows to the table as needed.

Concentration of substrate (uM)	Observed Rate (uM/sec)

Plotting data to estimate kinetic values

- Plot rate (y-axis) vs Substrate concentration (x-axis) and estimate the V_{\max} and K_M values using the plot. Plot the data as a scatter plot (no lines)



The V_{\max} estimate is the y-value where the graph levels out. The K_M is the x-value at $1/2 V_{\max}$.

- Record your estimates with units of V_{\max} and K_M in the table below.

V_{\max}		
K_M		

Fitting data to calculate V_{\max} and K_M

- To draw a trendline, you'll need to calculate theoretical values for rates using the estimated V_{\max} and K_M values from 2.1 and the Michaelis-Menten Equation. This is your V Calculated (V_{calc}). This will also help you calculate the correct V_{\max} and K_M values later on.
- In Microsoft Excel or Google Sheets, set up a table with headers as shown below:

Concentration of substrate	Observed rate	Theoretical rate (V_{Calc})	Squared residual		Variables	Estimates	Units
					V_{\max}		
					K_M		
					SSR		

5.1 Input your data table from step 1.3 into columns A and B.

5.2 Input your estimates of V_{max} and K_M from step 2.1 into the appropriate box in column G

- 6 In column C, insert the following code to use the Michaelis-Menten equation to calculate the theoretical rate given the estimates for V_{max} and K_M .

```
=($G$2 * A2)/($G$3 + A2)
```



This code uses the Michaelis-Menten equation $v = \frac{V_{max} \times [S]}{K_M + [S]}$ to calculate the rate based on the estimates. The \$ locks those values in so that the equation can be "dragged" down for all concentrations of substrate without needing to re-write the equation.

- 6.1 "Drag" the equation down until the theoretical rate has been calculated for all the concentrations of substrate. After this step the table will appear similar to below:

Conc	Obs	V Calc
10	0.05628	0.05
20	0.09715	0.08571429
50	0.15269	0.15
100	0.17799	0.2
200	0.26219	0.24

- 6.2 In Column D, insert the following code to calculate the squared residual value.

```
=(B2 - C2)^2
```



The squared residual will be used to calculate error and refine the estimates of V_{max} and K_M . It is the squared difference between the real rate and the calculated rate. For a perfect match the squared residual should be zero.

- 6.3 "Drag" the equation down until the squared residual has been calculated for all the concentrations of substrate.

- 6.4 In cell G4, insert the following code to calculate the sum of squared residuals. Replace DXX with D followed by the number of the last substrate concentration (e.g. D10 if there are 9 substrate concentrations).

The resulting number is the sum of squared residuals or SSR.

```
=sum(D2:DXX)
```

- 7 Use SOLVER to fit the data by changing the V_{\max} and K_M estimates to minimize the sum of squared residuals.

7.1

If Solver has not been used before, it may not have been installed as described below. If it has been used before skip to step 6.2.


Note: For some users of Excel, add-ins are disabled, therefore Google Sheets has to be used.

FOR EXCEL USERS:

If Solver has not been installed, it can be added to Excel by going to File -> Options -> Add-ins -> Manage Excel Add-ins -> Solver Add-in

FOR GOOGLE SHEETS USERS:

If Solver has not been installed, it can be added to Sheets by going to Add-ons -> Get Add-ons, and searching for Solver.

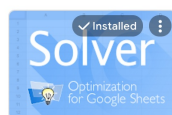

- 7.2 Load the Solver window. In Excel, Solver is found in the Data tab. In Sheets, it is found in the Add-ons menu.

- 7.3 The Solver window (example from Excel shown below) has 4 parameters to set:

- **Set Objective:** \$G\$4
- **To:** Min
- **By Changing Variable Cells:** \$G\$2, \$G\$3
- **Select a Solving Method:** GRG Nonlinear

The Objective and Variable cells can also be selected using the button at the right of the window and using the cursor to click on the appropriate cells.

Solver

Set Objective:

To:

☒ Max
☐ Min
☐ Value Of:

By Changing:

Subject To:

Add

Change

Delete

Solving Method:

Standard LSGRG Nonlinear

Reset All

Insert Example

Solve

Options

Copyright © 2019 Frontline Systems, Inc.

Solver in Google Sheets

Solver Parameters

Set Objective:

\$D\$19

To:

☐ Max
☒ Min
☐ Value Of:

0

By Changing Variable Cells:

'Example_JRKinetics_data (2)!'\$B\$10:\$B\$11

Add

Change

Delete

Reset All

Load/Save

☒ Make Unconstrained Variables Non-Negative

Select a Solving Method:

GRG Nonlinear

Options

Solving Method

Select the GRG Nonlinear engine for Solver Problems that are smooth nonlinear. Select the LP Simplex engine for linear Solver Problems, and select the Evolutionary engine for Solver problems that are non-smooth.

Close

Solve

Solver in Excel

7.4 Click Solve and wait for the fitting to finish (<10 seconds).



Solver will substitute in various values for V_{\max} and K_M to reduce the SSR value as much as possible which means that the theoretical values match the experimental values closely. At the end, the values will be changed automatically in the spreadsheet

7.5 In Sheets: Solver will run and change the values.

In Excel: A window will pop-up asking to accept the values and display if there were errors. Click OK and the estimated V_{\max} and K_M will be changed to the optimized values.

7.6 Record the values V_{\max} , K_M , and SSR below. These V_{\max} and K_M values are accurate compared to your original estimates.

Vmax	
Km	
SSR	

Visualizing the fit to the data

8 Add a second set of data to the original plot of observed rate vs. concentration of substrate.


8.1 Plot the theoretical (V Calc) rates vs. concentration of substrate.

8.2

Click on the newly plotted V Calc data and edit the data set to show only the line.

To show only the line, highlight the Vcalc data by left-clicking on the data and then right clicking to bring up the pop-up menu. Format the data series to add the line.

The plot should look similar to that shown below:

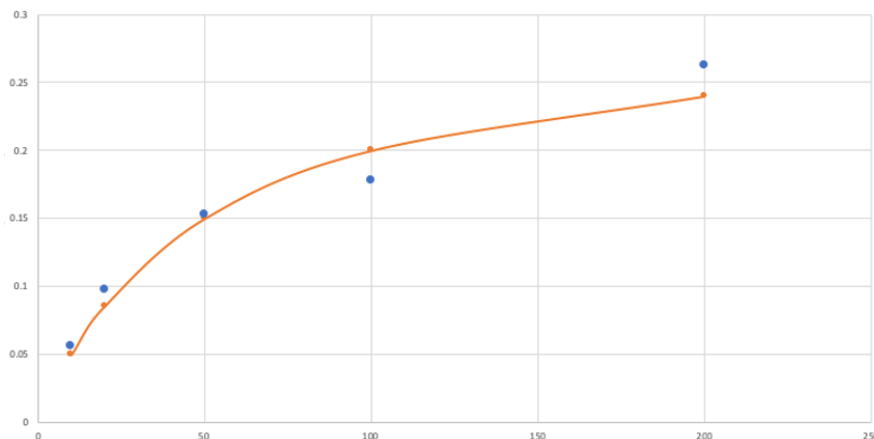

Label the plot axes with appropriate labels and units.



8.3

8.4 Save your plot as a picture and upload the image.