

JUN 05, 2023

OPEN ACCESS

DOI:
dx.doi.org/10.17504/protocols.io.81wgb6j51pk/v1

Protocol Citation: Lisa Schunk, Anja Cramer, Konstantin Bob, Ivan Calandra, Guido Heinz, Olaf Jöris, Joao Marreiros 2023. Enhancing lithic analysis: Introducing 3D-EdgeAngle as a semi-automated 3D digital method to systematically quantify stone tool edge angle and design.
protocols.io
<https://dx.doi.org/10.17504/protocols.io.81wgb6j51pk/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
 We use this protocol and it's working

Created: Aug 05, 2022

Last Modified: Jun 05, 2023

🌐 Enhancing lithic analysis: Introducing 3D-EdgeAngle as a semi-automated 3D digital method to systematically quantify stone tool edge angle and design

Lisa Schunk^{1,2,3}, Anja Cramer⁴, Konstantin Bob⁵, Ivan Calandra^{6,2},
 Olaf Guido Heinz⁴, Jöris^{7,3}, Joao Marreiros^{2,3,8}

¹Institute of Archaeology, Faculty of Historical and Pedagogical Sciences, University of Wroclaw, Poland;

²TraCEr. Laboratory for Traceology and Controlled Experiments. MONREPOS Archaeological Research Centre and Museum for Human Behavioural Evolution, LEIZA, Neuwied, Germany;

³Institute of Ancient Studies, Department of Prehistoric and Protohistoric Archaeology, Johannes Gutenberg University, Mainz, Germany;

⁴Leibniz Research Institute for Archaeology. Mainz, Germany;


⁵Scientific Computing and Bioinformatics, Institute of Computer Science, Johannes Gutenberg University, Mainz, Germany;

⁶Imaging Lab, LEIZA;

⁷MONREPOS Archaeological Research Centre and Museum for Human Behavioural Evolution, LEIZA, Schloss Monrepos, 56567 Neuwied, Germany;

⁸ICArEHB, Interdisciplinary Center for Archaeology and Evolution of Human Behaviour, University of Algarve, Faro, Portugal

TraCErMonreposLEIZA

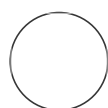

Lisa Schunk

Institute of Archaeology, Faculty of Historical and Pedagogical...

ABSTRACT

This protocol describes a method ("*3D-EdgeAngle*") to measure edge angles on 3D models. The method is semi-automated and script based. *3D-EdgeAngle* illustrates an objective way of measuring the edge angle at cross sections along the entire tool edge in defined steps and also allows measurements at different distances perpendicular to the edge. Depending on the scale of analysis, three parameters can be controlled accordingly. Within this protocol, all steps, starting from the acquisition of the 3D data to the statistical analysis, are described. Although the method is script based, *3D-EdgeAngle* can also be executed manually. The protocol illustrates both versions. The scripts applied to perform the edge angle measurements are incorporated in the protocol as Zenodo links.

Keywords: 3D model, edge
angle, lithic analysis,
standard method, tool
design, quantitative analysis

Samples

1 To demonstrate the applicability of the method, three samples have been selected:

1x experimental flake (EAP-flake)

1x Keilmessser (BU-072)

1x calibrated standard angle (WEM-60; gauge block)

1.1 Experimental flake

- made of Baltic flint (Southern Sweden, secondary deposit:



55.945852 N, 12.767851

E

knapped by Frank Moseler)


- one edge of the tool is in the distal part marginally, unifacially retouched



Ventral and dorsal view on the experimental flake (EAP-flake).

1.2 Keilmesser

- artefact from the Late Middle Palaeolithic site of Buhlen, Germany

 51.191022 N, 9.086585
E

- made of silicified schist

- the tool has one distinct, retouched active edge characterised by the scare of a lateral
tranchet blow

CITATION

Lisa Schunk (2022). Understanding Middle Palaeolithic asymmetric stone tool design and use. Functional analysis and controlled experiments to assess Neanderthal technology. Verlag des Römisch-Germanischen Zentralmuseums, Mainz.

LINK


<https://doi.org/10.11588/propylaeum.1076>



Upper, lateral and lower view on the Keilmesser (BU-072). Photos taken by I. Görner.

1.3 Calibrated standard angle

- made of certified steel
- calibrated to 60°40' (60.67°)
- coated with oil to avoid corrosion

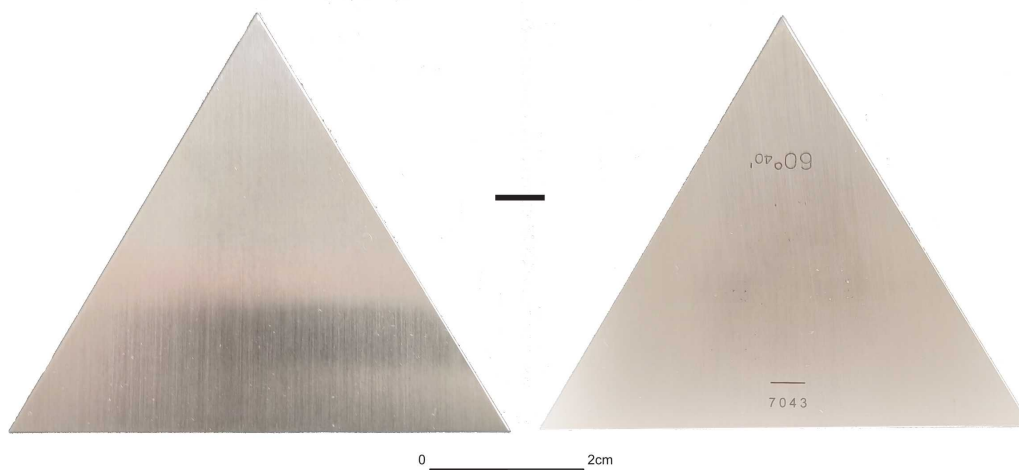


Back and front view of the calibrated standard angle (WEM-60).

Sample documentation

- 2 Sample documentation is separated into two part: generating the 3D model via scanning and processing the 3D model.

2.1 3D scanning

The three samples were scanned with an AICON smartScan-HE R8 from the manufacturer Hexagon (software version OptoCat 2018R1), featuring a blue light LED and two black and white cameras with 8 megapixels each. The S-150 FOV (field of view) used has a point-to-point distance of 33 μm .

Equipment	
smartScan-HE R8	NAME
3D structured light scanner	TYPE
AICON	BRAND
-	SKU
https://www.hexagonmi.com/products/structured-light-scanners/aicon-smartscan	LINK
S-150 FOV, resolution of 33 μm	SPECIFICATIONS
	

Software	
OptoCat	NAME
Hexagon Manufacturing Intelligence Software	DEVELOPER

The following settings have been used to scan the samples:

I. Turning on the system


Before using the AICON smartScan-HE R8 the system had to be warmed up for one hour in advance.

II. Calibration


The scanner was calibrated with the appropriate calibration plate. Accuracy was reached when the deviation of the calibration was $\leq 7.5 \mu\text{m}$ for the FOV 150.

III. Placing the sample

Each sample was placed in the middle of the turntable. Polystyrene was used to keep the samples in place. All three samples were scanned in an upright position. The sample holder (polystyrene) was covered with a black cloth.

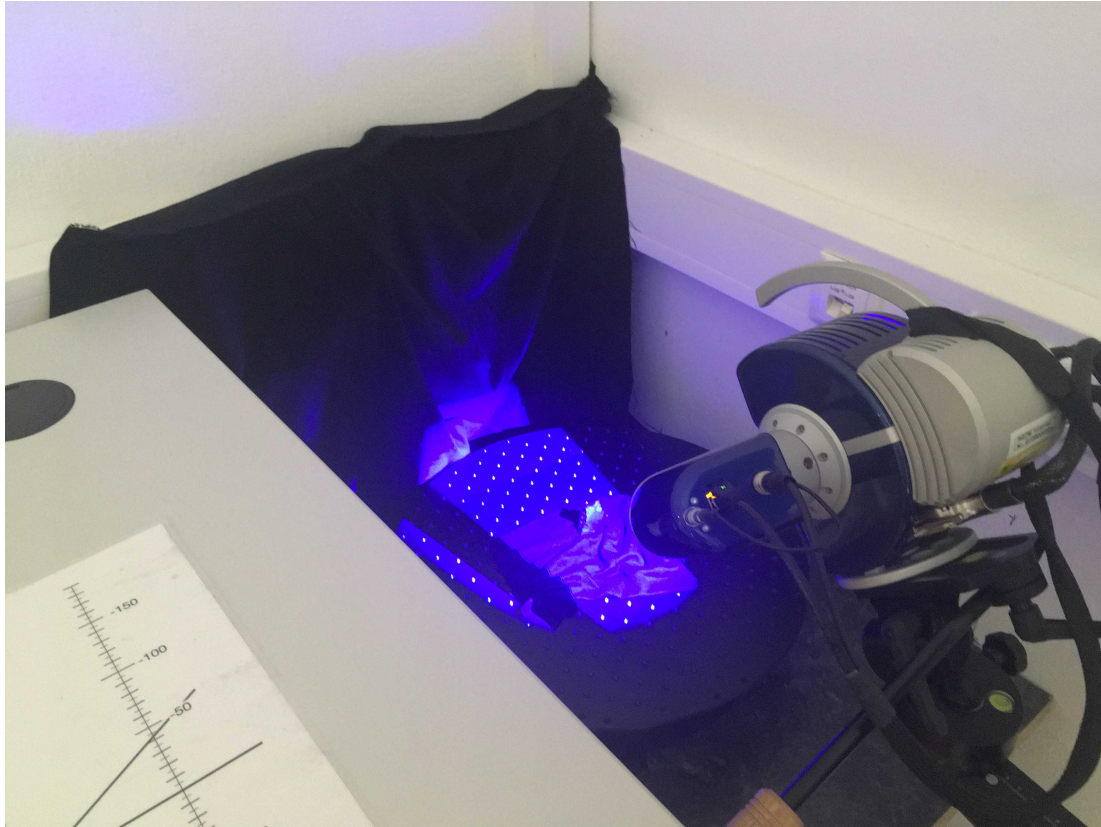


AICON smartScan-HE R8 with the FOV S-150 and the turntable. The sample (here Keilmesser) is placed in the middle of the turntable.

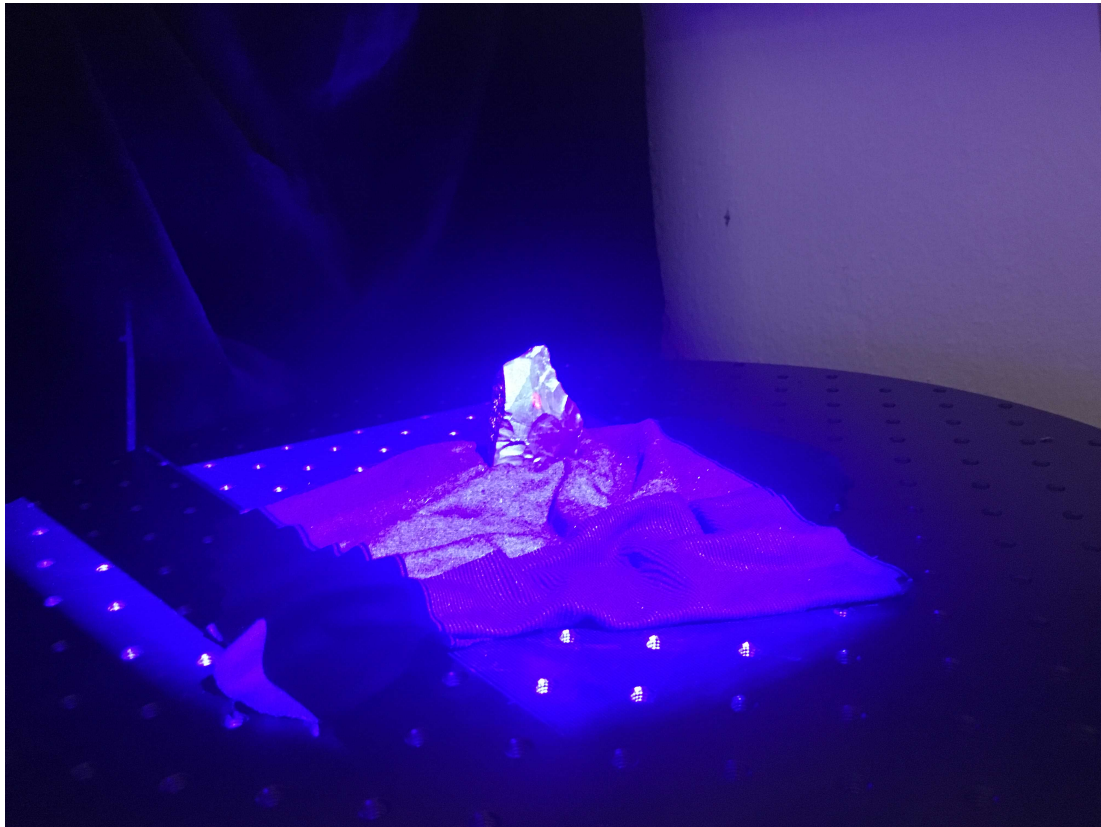

The sample (here Keilmesser) is placed upright with a sample holder, which is covered with a black cloth.

Note that for the calibrated standard angle some extra target marks on the sample were needed to make the aligning between the single scans possible. Thus, adhesive paper tags with symbols were placed on the calibrated standard angle. Despite ungreasing the metal surface of the sample, the adhesive paper tags showed a tendency to peel off after a period of time.

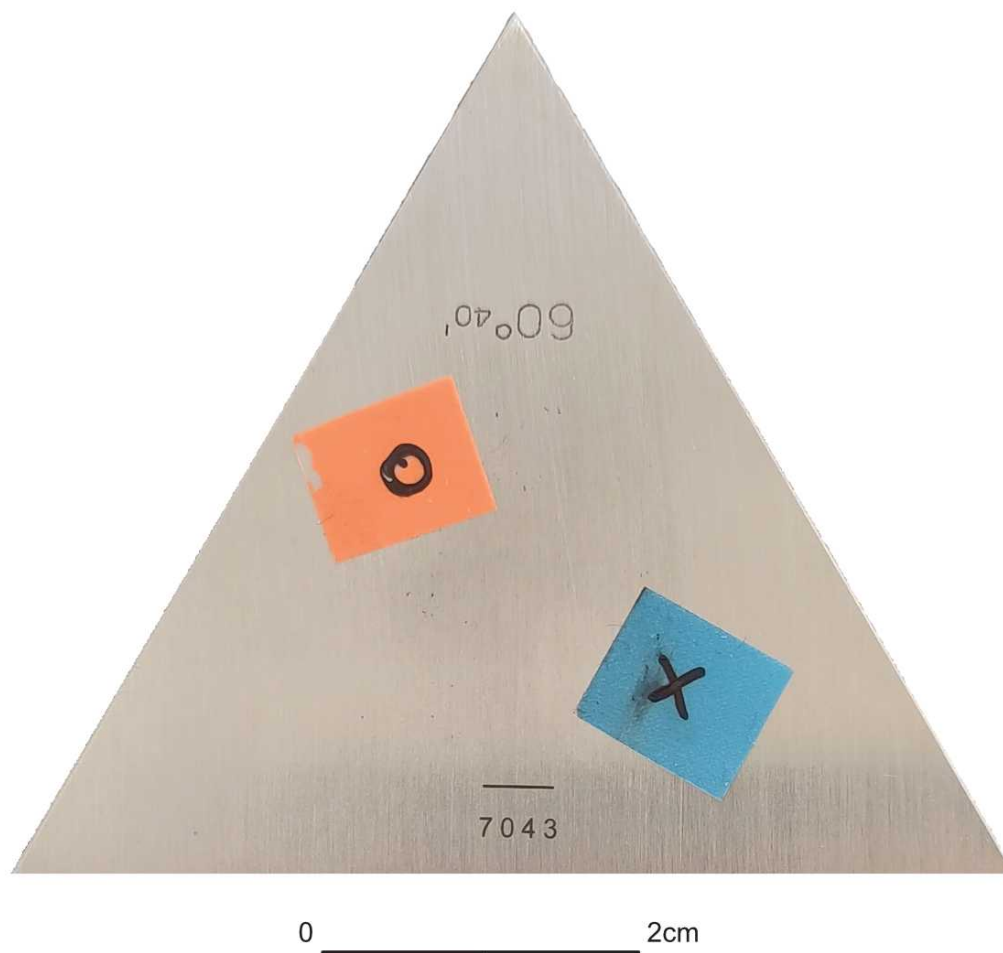

Calibrated standard angle with target marks (WEM).

IV. Acquisition settings

"Automatic " was chosen as a scanning mode. This means, based on the rotation of the turntable, the scans get aligned.


A new project was created named with the following structure

"objectnumber_FOV_yyyymmdd".


The software offers the possibility to use clip planes to define parts of the FOV that will not be considered during the processing of the 3D model (i.e. the turntable and the sample holder). Thus, clip planes were defined, eliminating everything visible below the sample (Z-axis). For doing so, the software needs to be taught where the rotation axis is relative to the object and scanner.

For initialising the scanning process, the number of grabs per scan (= different shutter times for an HDR scan) was selected with 4. The HDR-factor was set to manual so that the shutter time per grabs needed to be defined (light exposure). Different exposure times were selected between the grabs, guaranteeing that every part of the sample was sufficiently exposed for at

least one of the grabs. The gain was 1x.

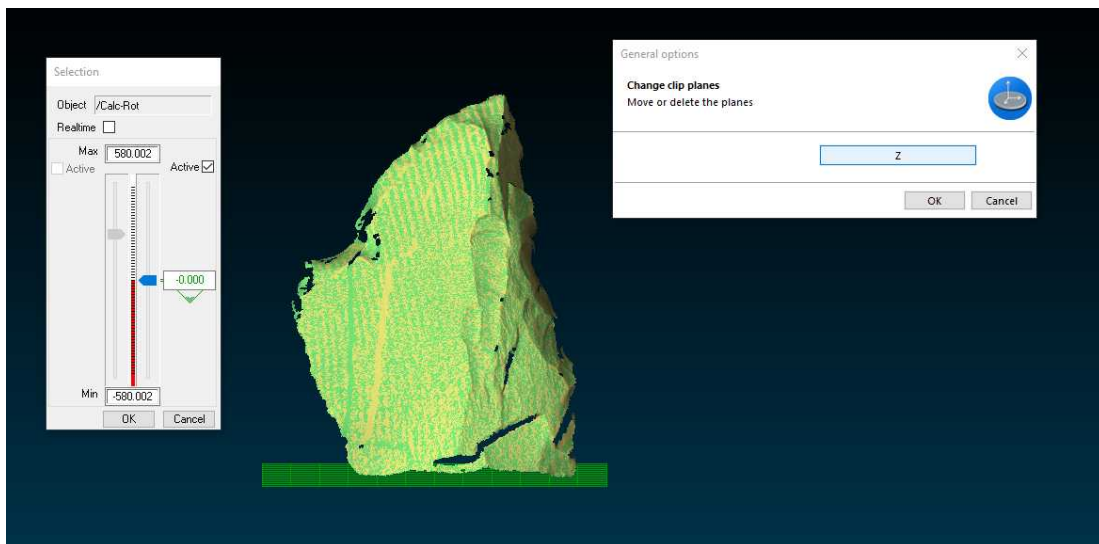

Definition of the clip planes (here Z-axis) in OptoCat 2018R1.

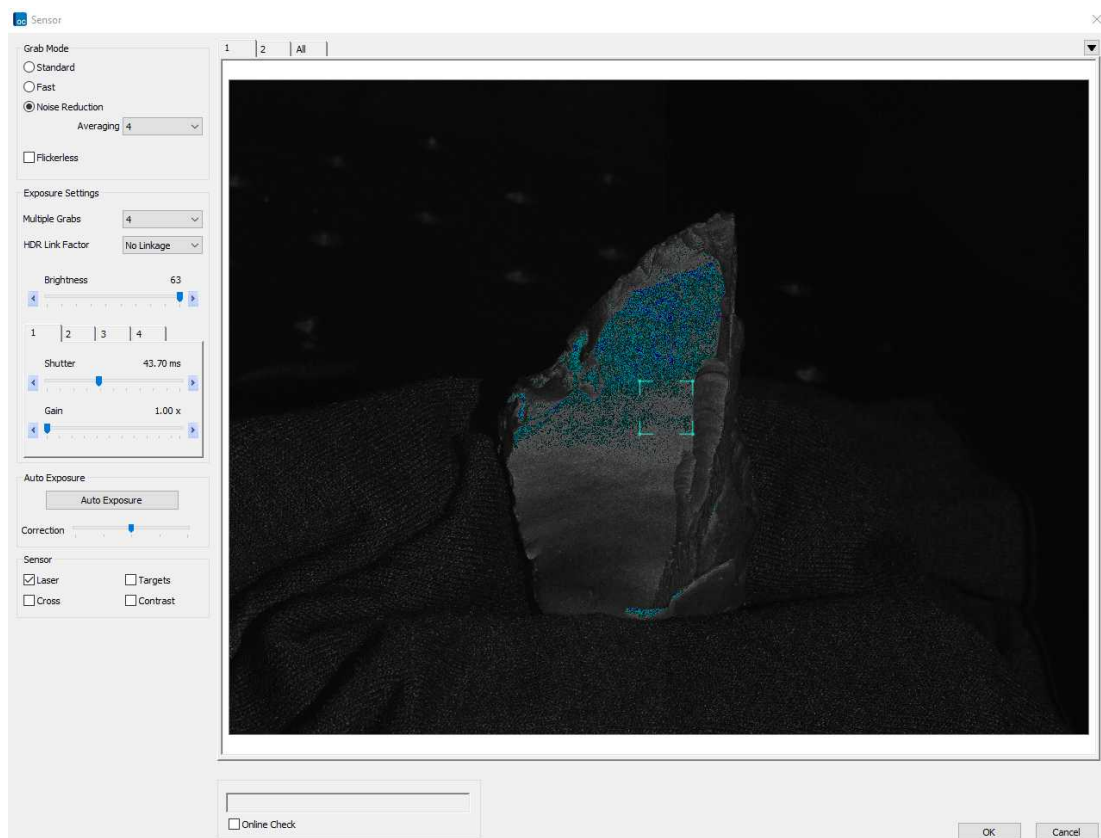

Selected sensor settings in OptoCat 2018R1.


V. Data acquisition

The number of scans per rotation (i.e. 360°) was defined with ten. For all scans, the chosen light settings stayed the same. After these ten scans, the sample was turned so that the previously covered parts of the samples were visible. A new 360° rotation with ten scans was performed. Sometimes, individual scans needed to be done before the turn of the sample could take place (when the alignment was not possible otherwise). After these two full rotations, the samples were scanned with as many single scans as needed to cover the entire surface of the samples adequately by moving the sample.

Approximately 20 to 27 scans per tool were needed to create a closed model.

Automatic moving unit

Please define the number of positions



Use taught positions

☐

Scans per rotation

Use taught presets

☐

Grab presets

Calculate 3D data

☒

Sensor to table changed

☐

Wait for data after each scan

☐

Wait to edit data

☐

Merge

☐


OK

Cancel

Selected acquisition settings in automatic mode in OptoCat 2018R1.

VI. Alignment

After acquiring the scans, the alignment was started. After starting the process, the software calculates a RMS (root mean square) value, which was $\leq 7.5 \mu\text{m} \times \text{deviation}$.

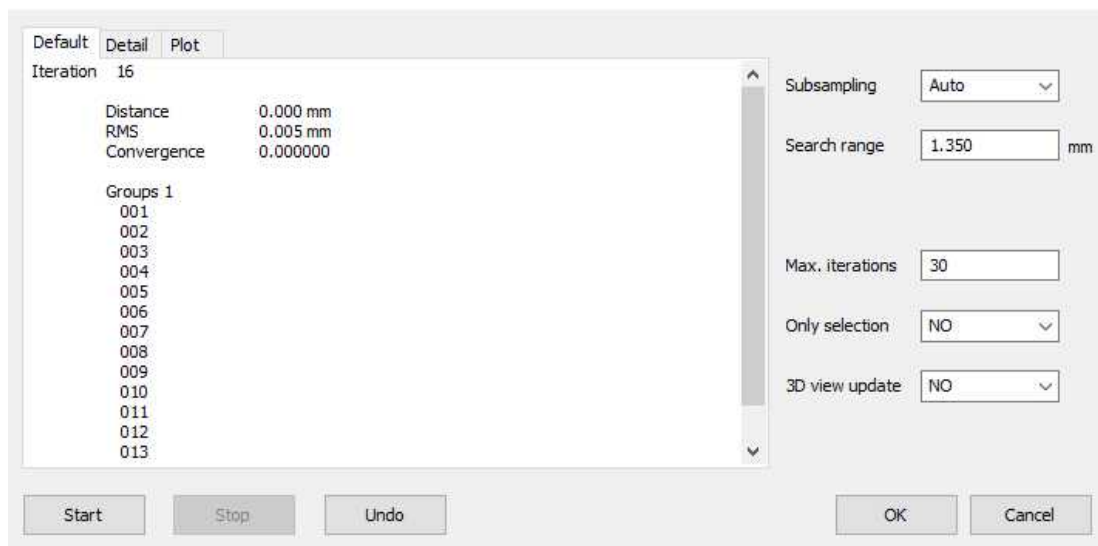


Completed alignment in OptoCat 2018R1, showing the RMS value.

VII. Data processing

The single 3D data (raw data) was opened in full resolution. The 3D data was cleaned by using the editing tools available in the software. Doing so, scanned data not belonging to the actual sample could be removed (e.g. background). After this, the data was aligned again (see step VI.). In a final step, the cleaned 3D data was merged in full resolution. The raw data was saved and the 3D model was exported in an STL-format.

Merge
Please define the parameters



Resolution
Full

Surface
Standard

Clean model
☒

Save in *.ply-format
☐

Name
BU-075_S-50_20190905

Remap color information
☐

Spotlight correction
☐

Screenshot of 3D result
☐

OK
Cancel

Selected merge settings in OptoCat 2018R1.

2.2 Editing of the 3D model

One prerequisite for the application of 3D-EdgeAngle method is a closed 3D model. Thus, the 3D model in the STL-format was edited in the free software GOM Inspect 2018 (2018 Hotfix 2, Rev. 111729).

Software		
GOM Inspect		NAME
Hotfix 2, Rev. 111729, build 2018-08-22		OS
GOM		DEVELOPER

The following steps were completed for each sample:


I. Import 3D model

The file was imported and opened as a STL-format, as new part. As "target element" type "mesh" was chosen.

In the right docking area, all information about the mesh could be seen.


II. Preparation 3D model

Within the editing menu, individual points, which are not entirely connected to the 3D mesh were selected and deleted (edit -> selection in 3D -> other selections -> select patches by number of points -> max. percentage 5.0%; edit -> delete selected points).

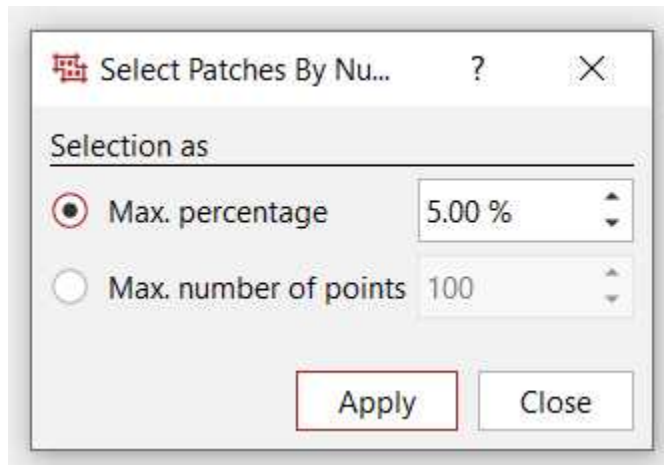

Maximal percentage chosen in the "selection in 3D" menu, GOM Inspect 2018.

Then, existing holes in the mesh were filled by selecting all points of elements of the mesh (mesh turns red and number of holes can be checked in information; operations -> mesh-> close holes -> automatically/interactively -> chose the max. hole size + max. number of edges).

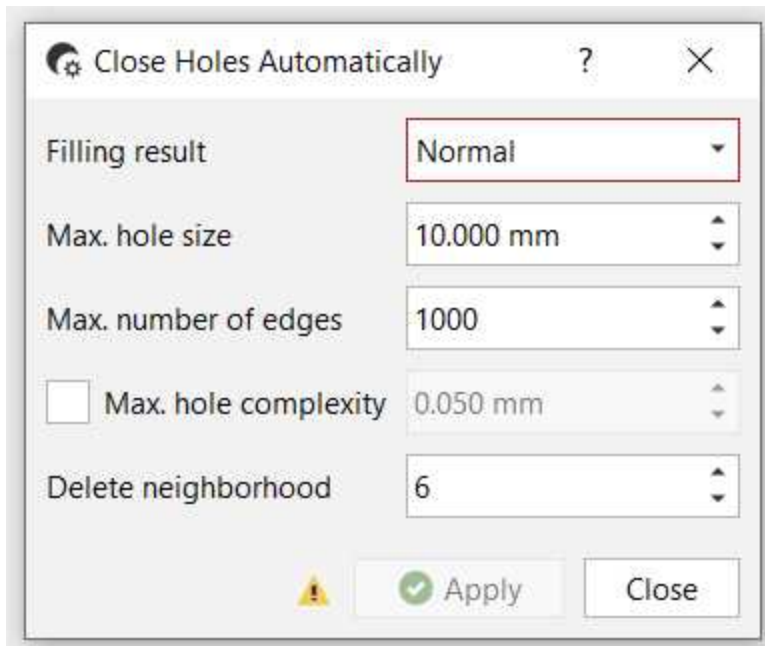

Settings for closing the holes in GOM Inspect 2018.

After doing this step once, mesh errors were additionally eliminated (operations -> mesh -> other -> eliminate mesh errors).

Depending on the number of holes these two steps were repeated until zero holes were left.

III. Alignment

In a final step, the samples were orientated (e.g., in case of the artefact the sample was orientated according to its technological features) following these steps: operations -> alignment -> manual alignment -> set matrix (first XYZ-rotation, then XYZ-translation)).

IV. Export

Each project was saved individually as a GOM Inspect file (*.ginspect). Moreover, each 3D model was again exported as STL.

Note

A dataset with a STL file for each 3D model can be found in open access on [Zenodo](https://doi.org/10.5281/zenodo.7326241): <https://doi.org/10.5281/zenodo.7326241>.

Definition of the edge

3 Surface curve

With a closed 3D model as a basis, some preparatory steps are necessary before calculating the edge angle: the digitalisation of the edge. The following step is manually done in GOM Inspect (GOM Software 2018, Hotfix 1, Rev. 111729, build 2018-08-22).

Software

GOM Inspect

NAME


Hotfix 2, Rev. 111729, build 2018-08-22

OS

GOM

DEVELOPER

The edge of interest for calculating the edge angle needs to be defined. This step has to be done manually to reach the intended accuracy. Therefore, a "surface curve" has been added to each 3D model. This means, a polyline was defined by tracking the edge of the object as precisely as possible and later exported in an IGES-format.

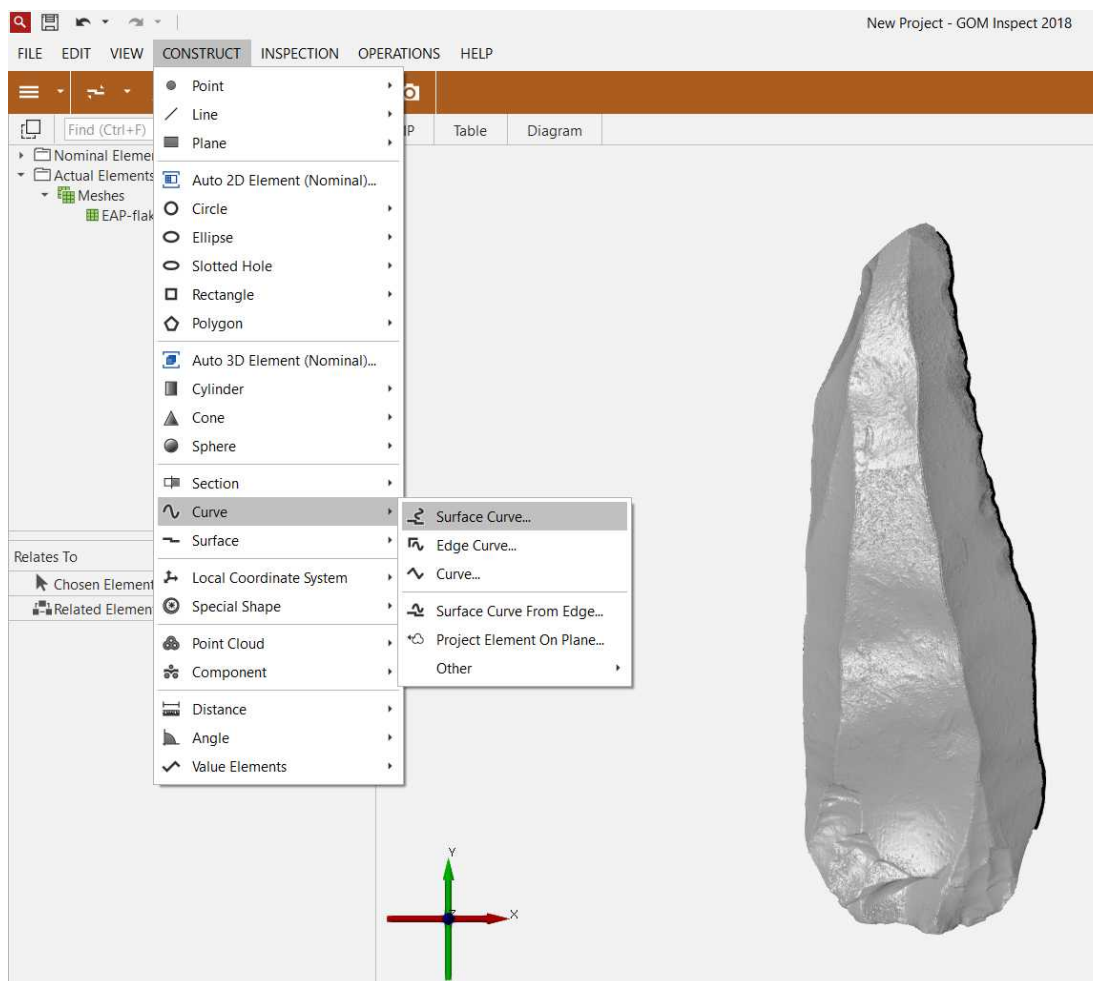

Construction of the surface curve in GOM Inspect 2018. Here the line is already defined for the experimental flake (see black line on the right side of the sample).

Step 3 includes a Step case.

Manual version

Script-based version

step case

Manual version

The steps described from #4 to #6 can be either done with scripts or manually.

4

Safety information


The following steps are done manually with the free version of GOM Inspect.

Software


GOM Inspect	NAME
Hotfix 2, Rev. 111729, build 2018-08-22	OS
GOM	DEVELOPER

Reference curve

Since the surface curve, which is following the tool's edge, can be complex with several changes in direction, a second, so-called "reference curve" has to be defined. The reference curve is a smoothed version of the surface curve. As a sample distance, 1 mm has been selected. The filter radius for smoothening the curve is defined as high.

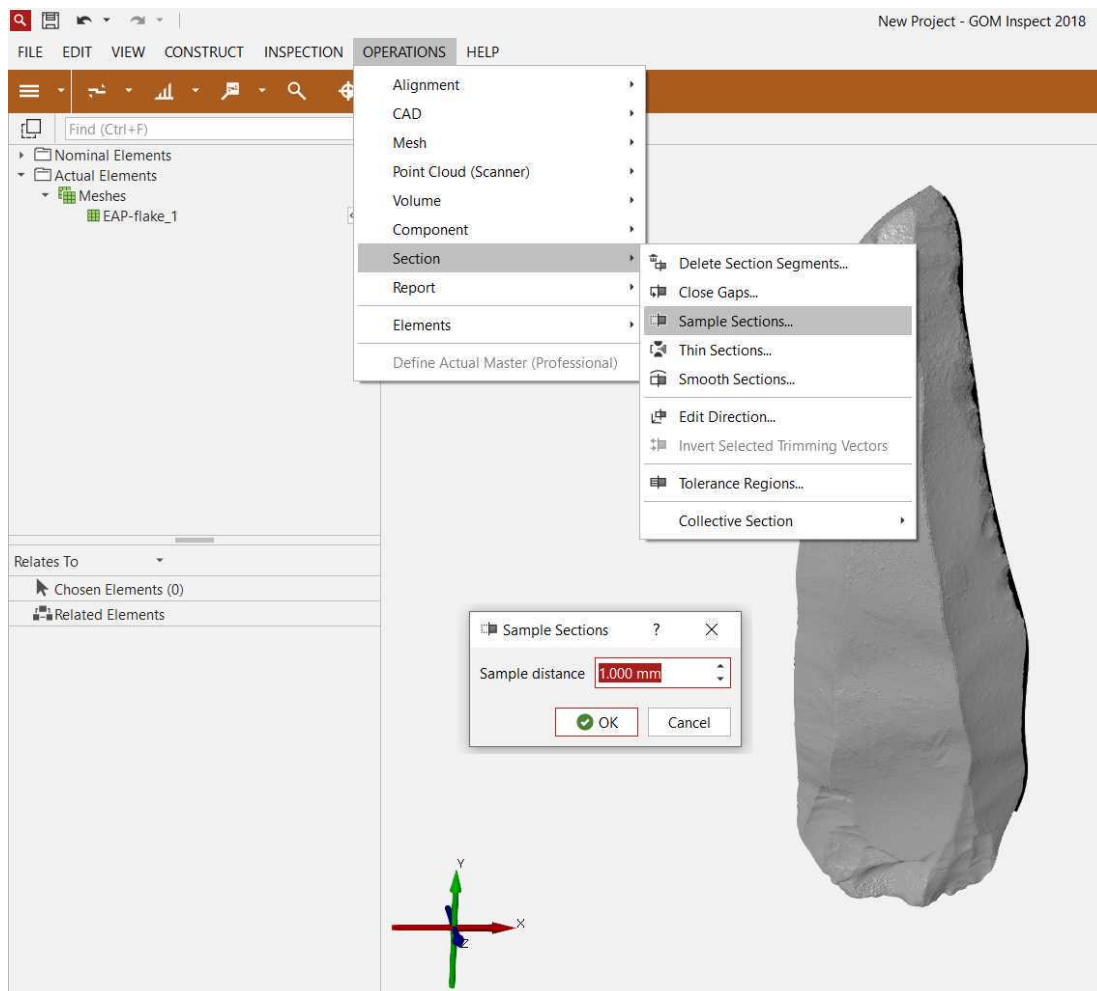

Sample section settings in GOM inspect 2018 to create a "reference curve".

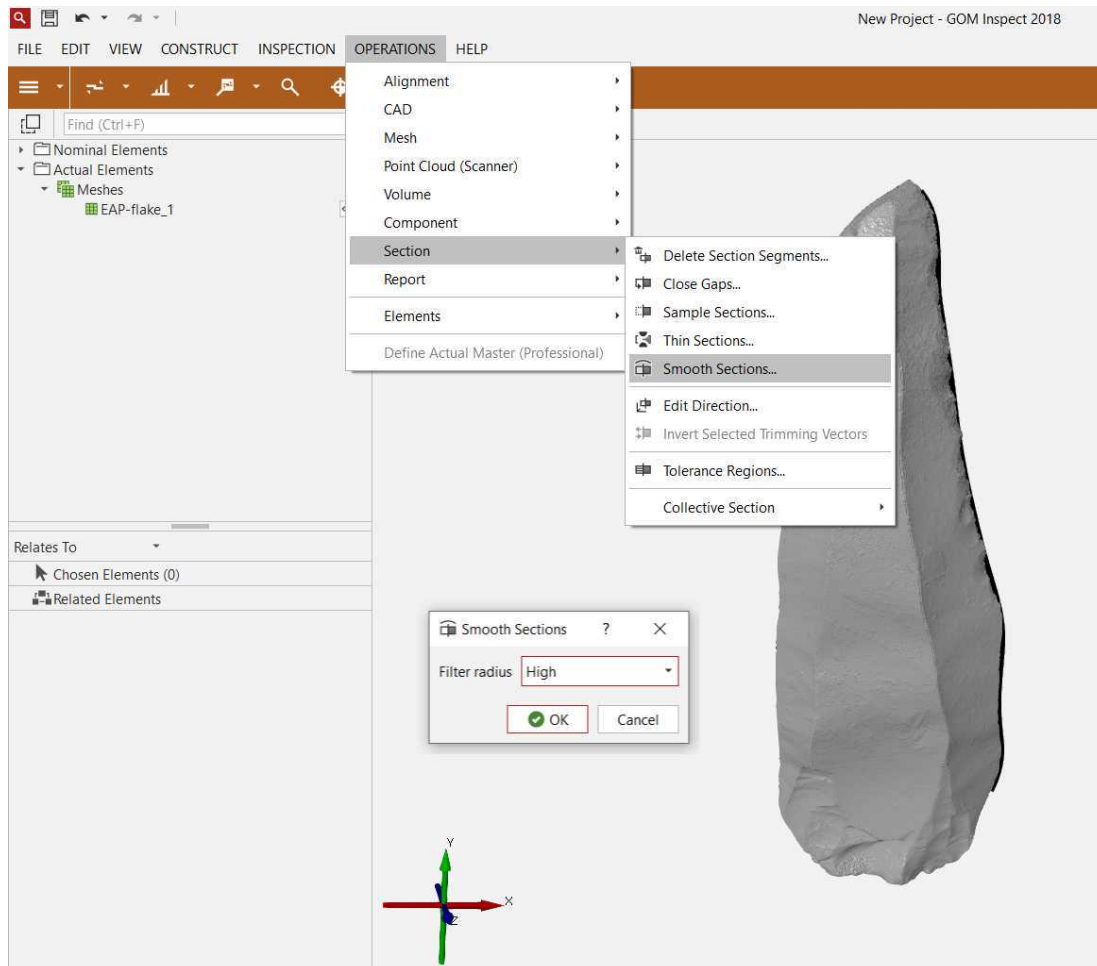

Settings for smoothening the curve in GOM Inspect 2018.

Note that these settings have been chosen based on the resolution of the 3D model. Depending on the resolution/size/accuracy other parameters have to be selected. The "reference curve" has been exported in an IGES-format.


Sections

5 Construction of the sections

The "reference curve" is used for only one purpose, which is the definition of the sections. The sections are the cross profiles originating perpendicular from the "reference curve" and following the artefact surface for a defined length. These sections are also polylines as the vertical polylines along the edge and are always orientated in a 90° angle towards the "reference curve".


To present 3D-EdgeAngle, **10 sections** per edge have been chosen as a parameter. In the manual version, however, the number of sections cannot be defined. Instead, the distance between the individual sections has to be determined. Thus, the distance has to be calculated.

To do so, the "reference curve" has to be converted into an "actual element".

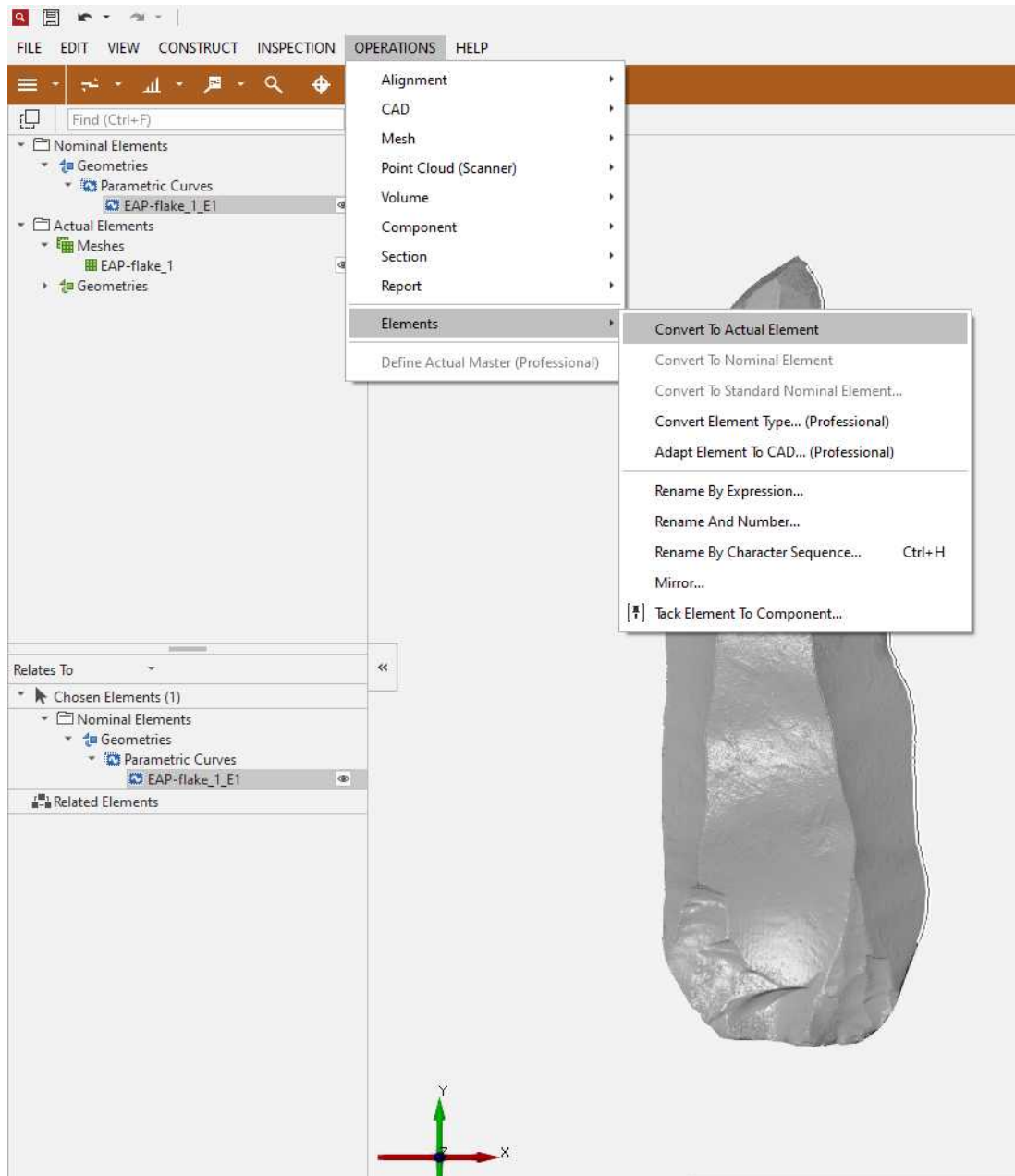


Conversion of the "reference curve" in GOM Inspect 2018.

As soon as the "reference curve" is converted into an "actual element", the length of the curve can be measured.

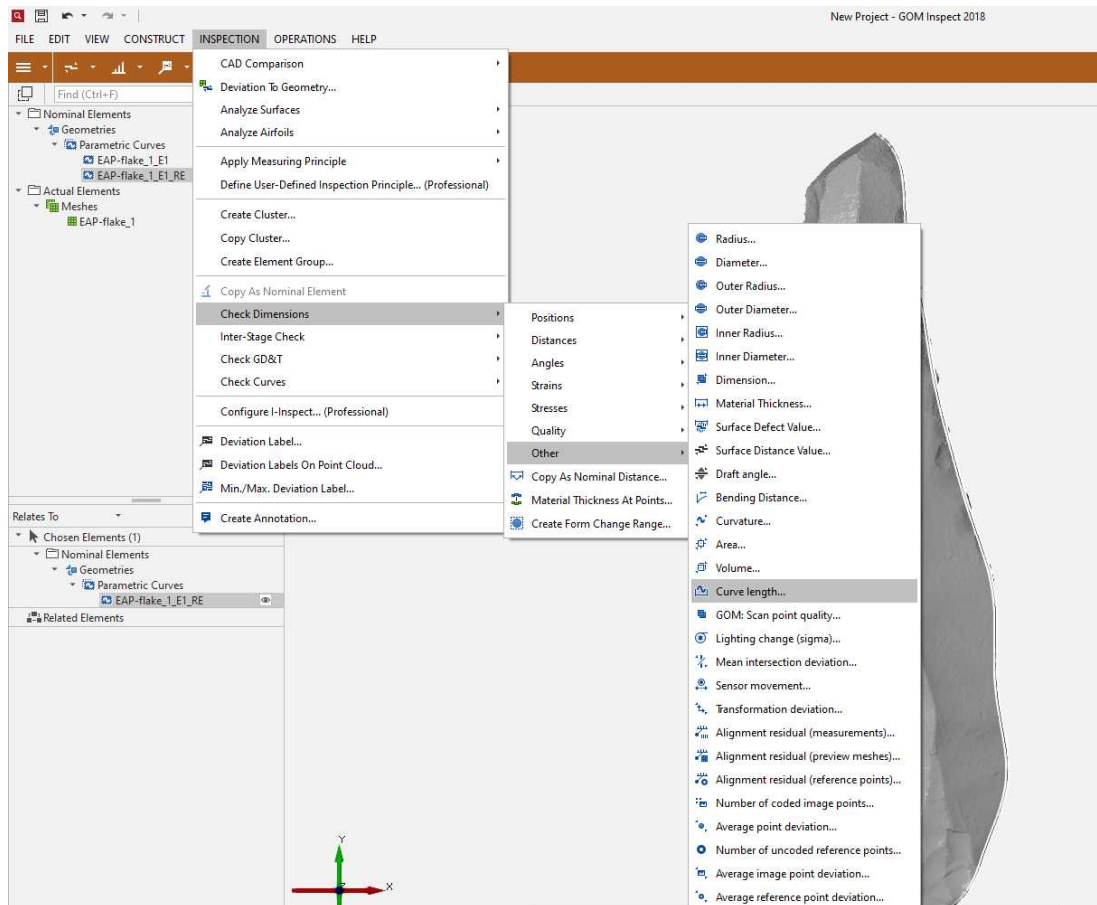


Measuring the length of the curve in GOM Inspect 2018.

The length in case of this "reference curve" is 82.29 mm. For ten sections in total, the length of the curve has to be divided by nine. Thus, the section distance is 9.14 mm.

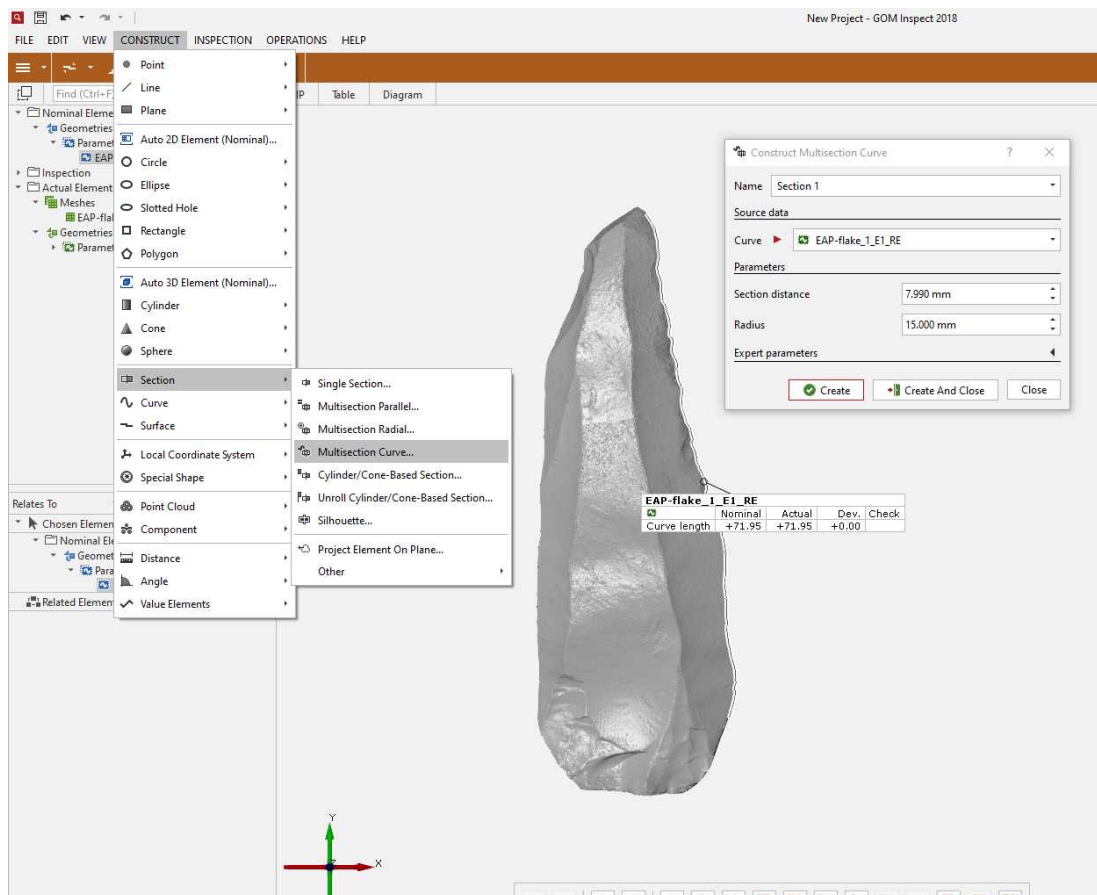

Construction of the sections in GOM Inspect 2018.

After the sections have been constructed, each section can be exported individually in an IGES format.

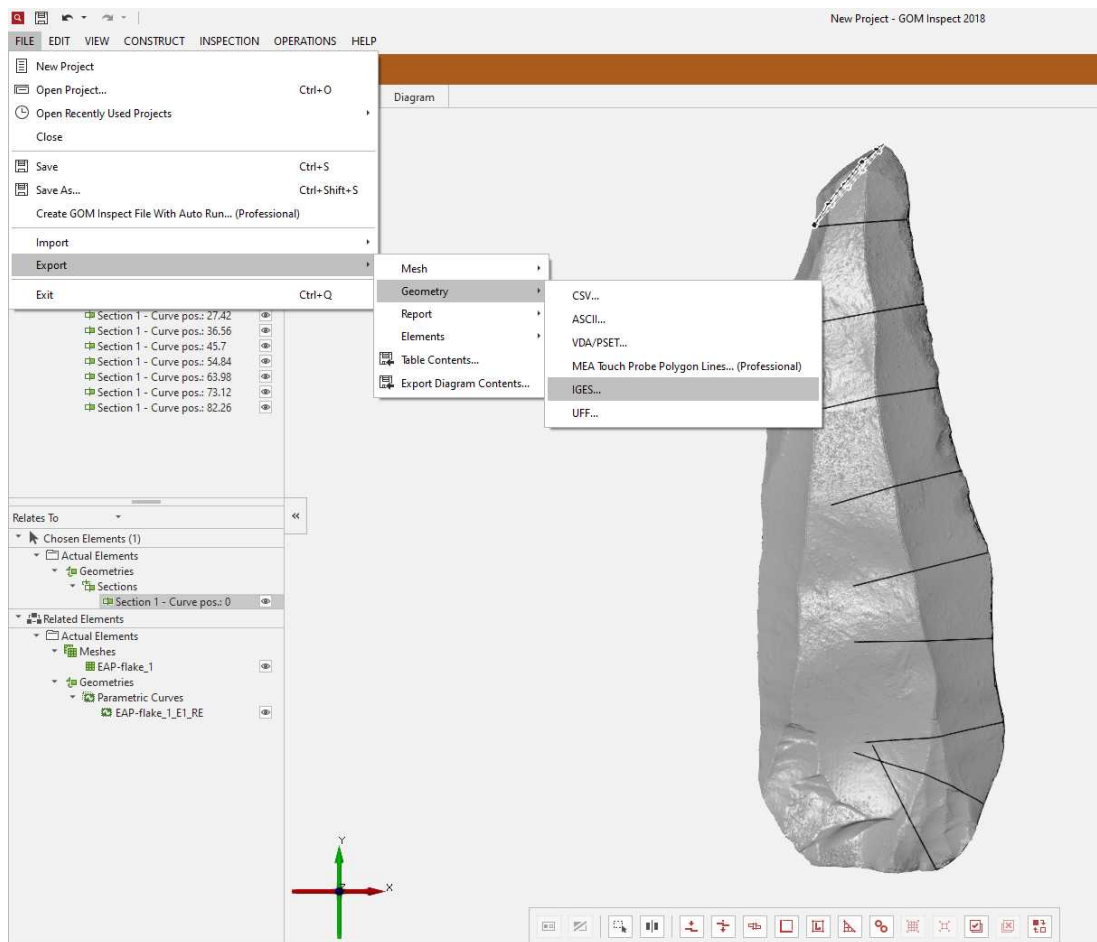

Export of the sections in GOM Inspect 2018.


5.1 Transformation of the sections

With the next steps, the individual sections will be projected onto a 2D plane. This is not necessary for the calculation of the edge angles, but helpful in case of possible further analyses (e.g. curve sketching).

I. Conversion

Since the "reference curve" was only created and used to define the sections, these steps are based on the original "surface curve" again.


First, this "surface curve" needs to be converted into an "actual element".

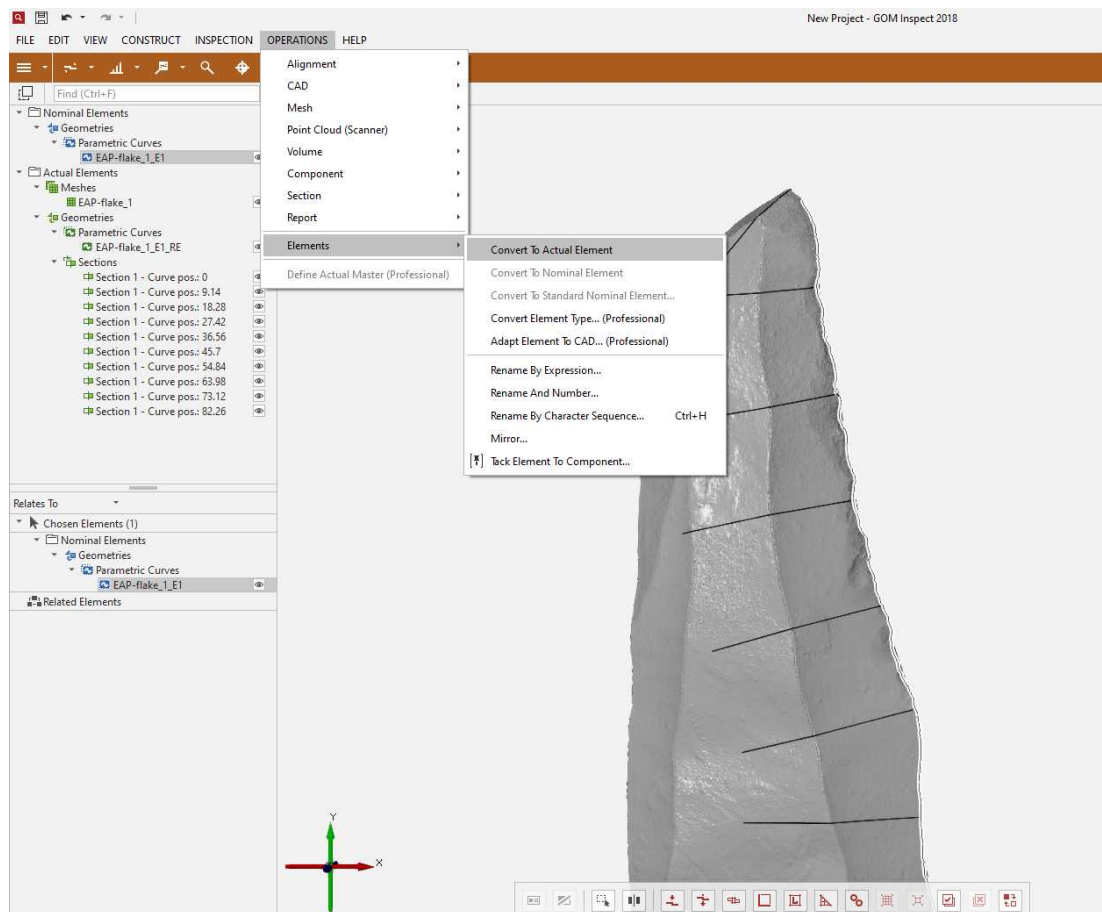

Conversion of the "surface curve" in GOM Inspect 2018.

II. Point definition


Three points need to be defined on each of the individual sections. The first point should be equal to the beginning of the section; the second equal to the end of the section. The third point should be somewhere between the first and the second point (~ middle of the section, but does not need to be precise).

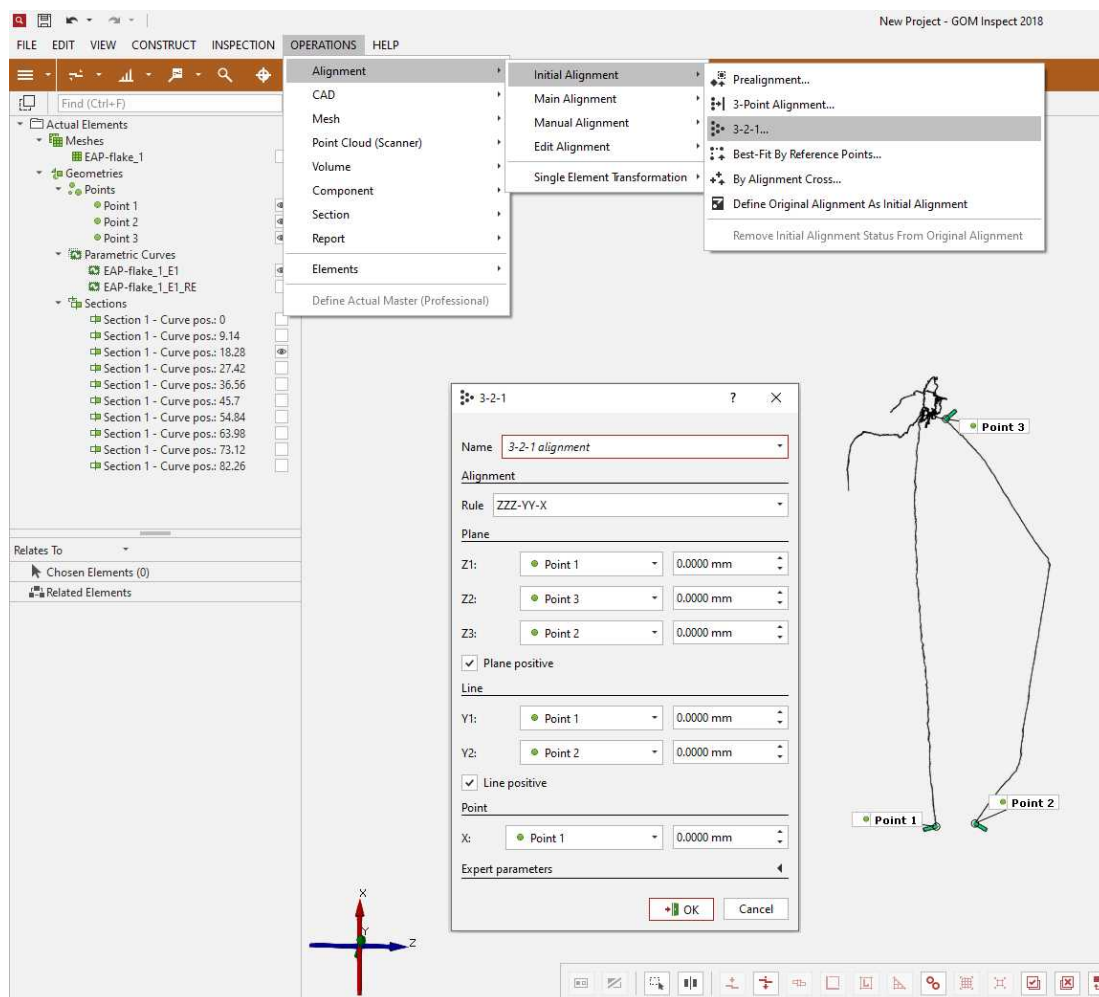

"3-2-1-alignment" in GOM Inspect 2018.


IV. Definition of the intersection

The intersection between the section and the "surface curve" needed to be defined. Thus, the distance between both elements was calculated in a first step.

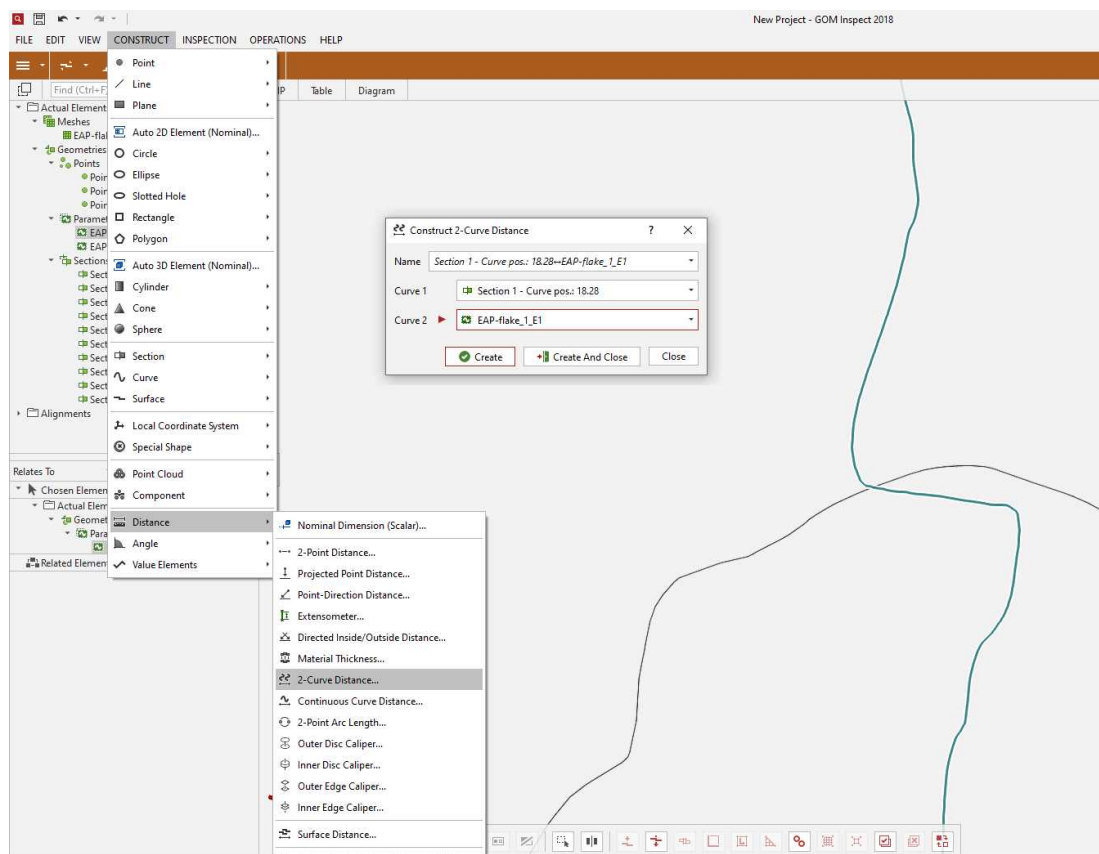


Calculation of the distance between the "surface curve" (here turquoise) and a section (here black) in GOM Inspect 2018.

At the intersection, a point has been constructed (zooming in needed to place the point correctly).

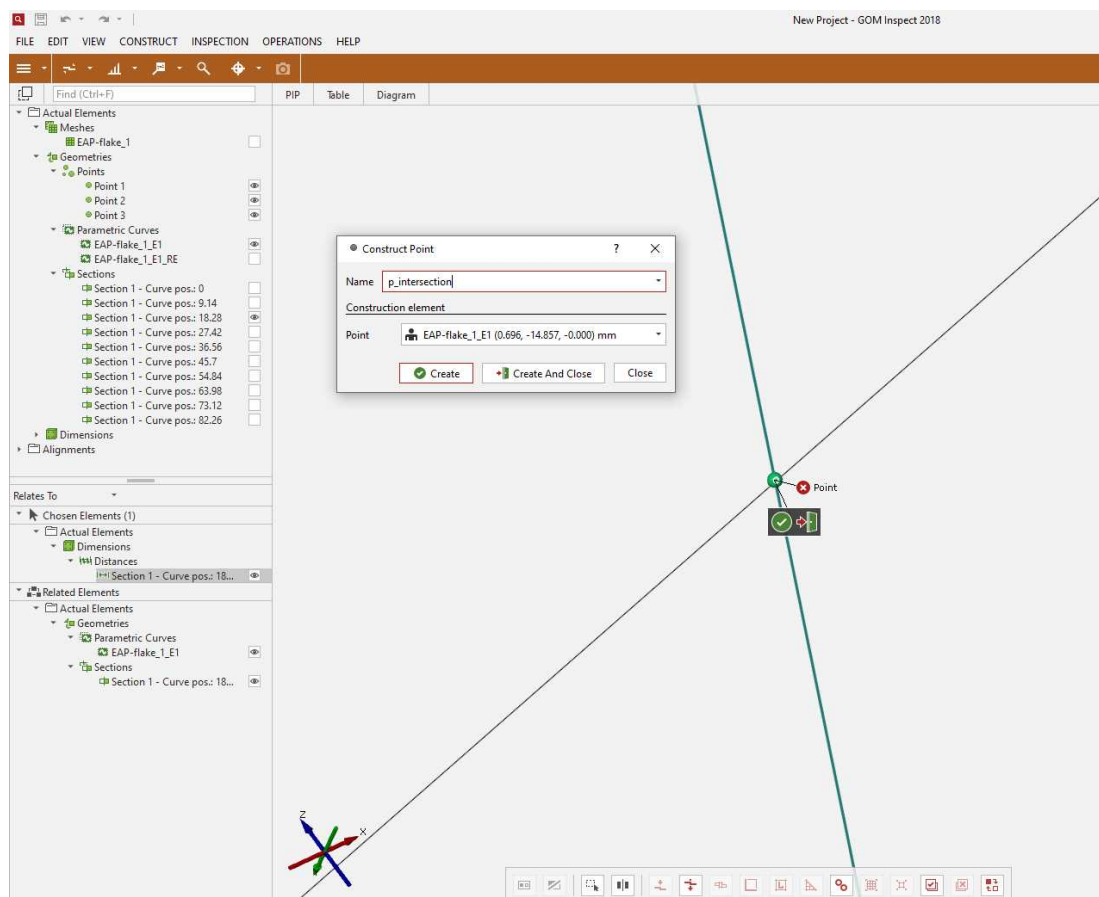

Construction of a point at the intersection in GOM Inspect 2018.

Finally, the constructed point at the intersection was set up as an origin of coordinates (thus, represents the digitalised surface curve on the section) .

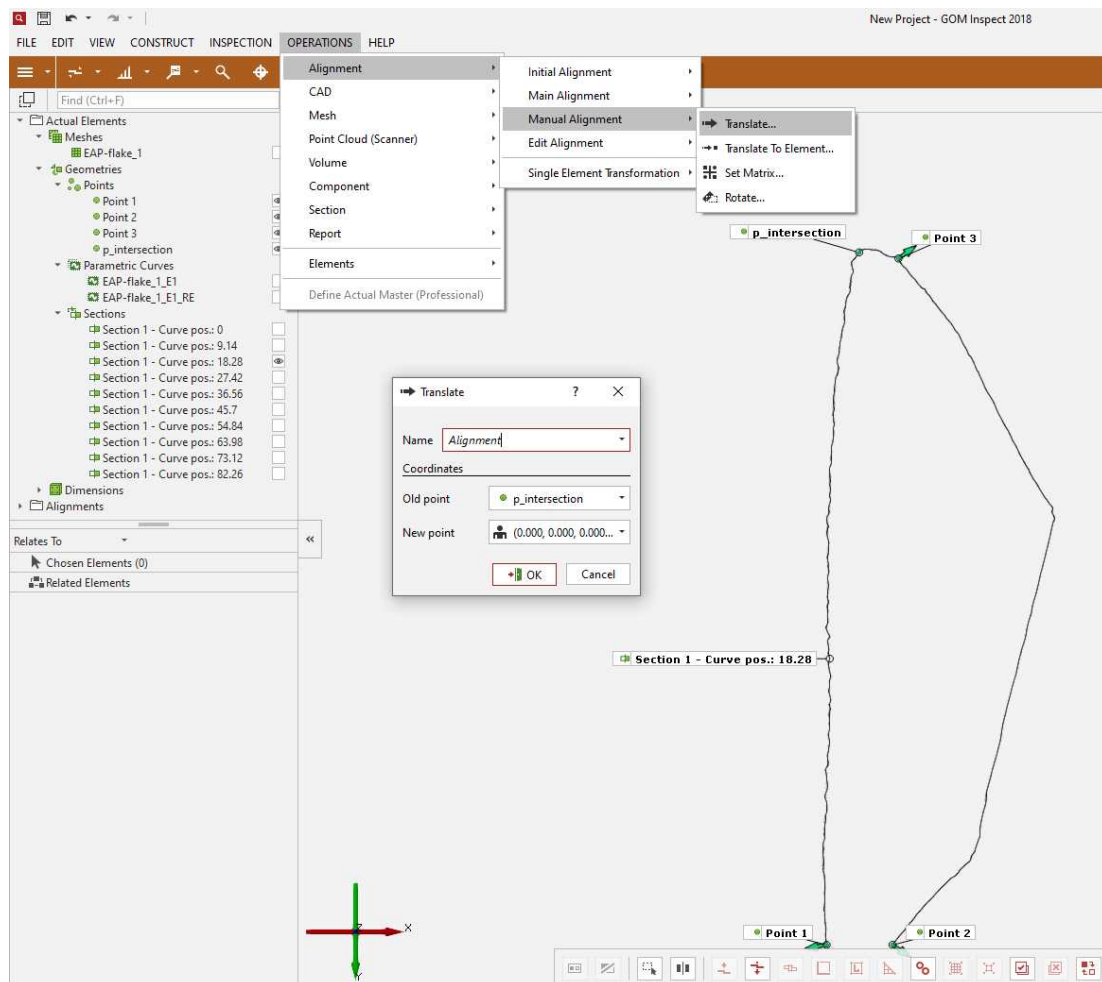

Definition of the origin of coordinates in GOM inspect 2018.

V. Export

In a last step, the new orientated curve was exported in different formats.

The formats are:

- PNG
- IGES
- ASCII
- gom inspect project

Note that these formats are not all needed to calculate the edge angle in a later step, but they can be useful for a quick visual control (PNG), additional analysis (e.g. curve sketching) etc.

Calculation of the edge angle

6 a) 3-points procedure

The "3-points" measurement procedure is modelled on the "*caliper method*" described by H. Dibble and M. Bernard in 1980.


CITATION

Harald Dibble, Mary C. Bernard (1980). A Comparative Study of Basic Edge Angle Measurement Techniques. *American Antiquity*, 45(4), 857-865.

LINK

<https://doi.org/10.2307/280156>


Whereas the result of the "caliper method" is based on a measurement taken by a special modified calliper at a known distance from the edge of the artefact by using a formula, the "3-points" digital measurement is using the same topographical indices, here represented by 3 points. The first point is the intersection between the vertical polyline (the "surface curve") and the horizontal polyline (the section). The two other points are placed on the section, in a defined distance away from the intersection (see #6.1). Following this triangulation, the edge angle can be calculated between these three points.

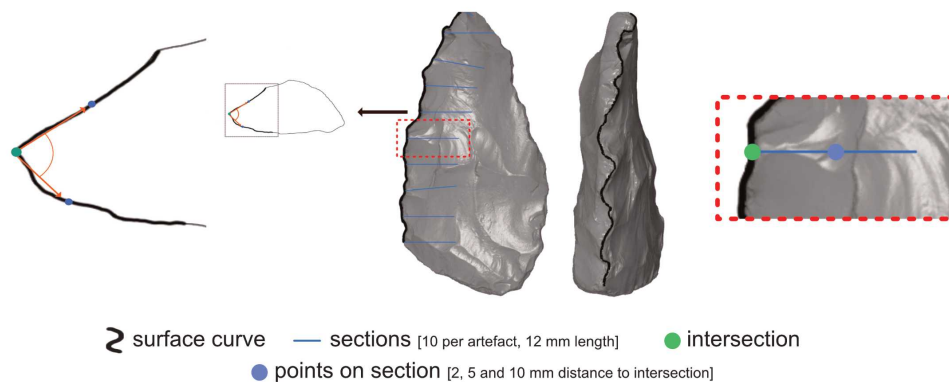

Schematic illustration of the "3-points" measurement procedure.

b) 2-lines procedure


The "2-lines" measurement procedure takes the intersection between the two polylines, as described in #6a, as one point. Beginning with this point, it follows the section on both surfaces for a given distance. Until this step, the two procedures are identical. Unlike the previous procedure, these points on the section (blue dots on the figure below) are used as starting points for the definition of two constructed lines. The constructed lines have a defined length (see #6.1) and take the points as a centre from where they expand in both directions. The extremities of these lines are placed on the section (blue lines on figure below) Thus, the calculation of the edge angle takes the two lines and the intersection into account.

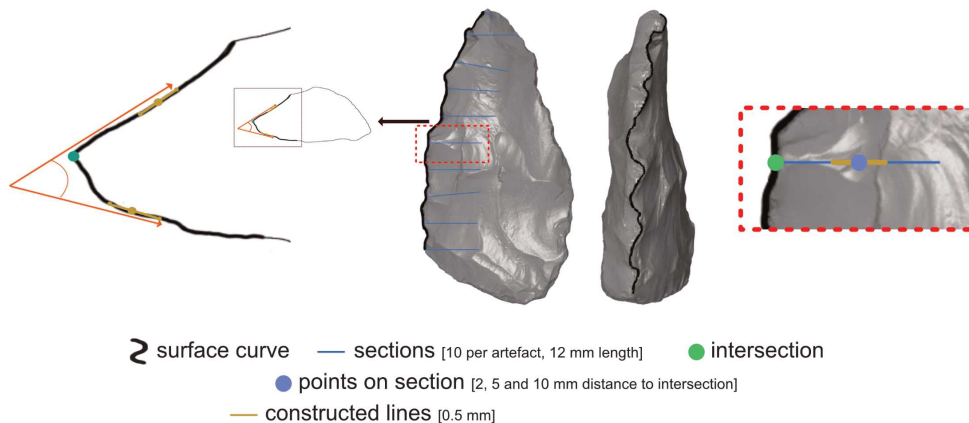

Schematic illustration of the "2-lines" measurement procedure.

c) best-fit procedure

The third way of calculating the edge angle, the "best-fit" procedure, is similar to the "2-lines" measurement. The difference is that the lines are constructed by finding the best fit along the defined length of the section. Therefore, the extremities of these lines are not necessarily placed on the section (blue lines on figure below). This has the advantage of reducing a possible error caused by small changes in the section direction.

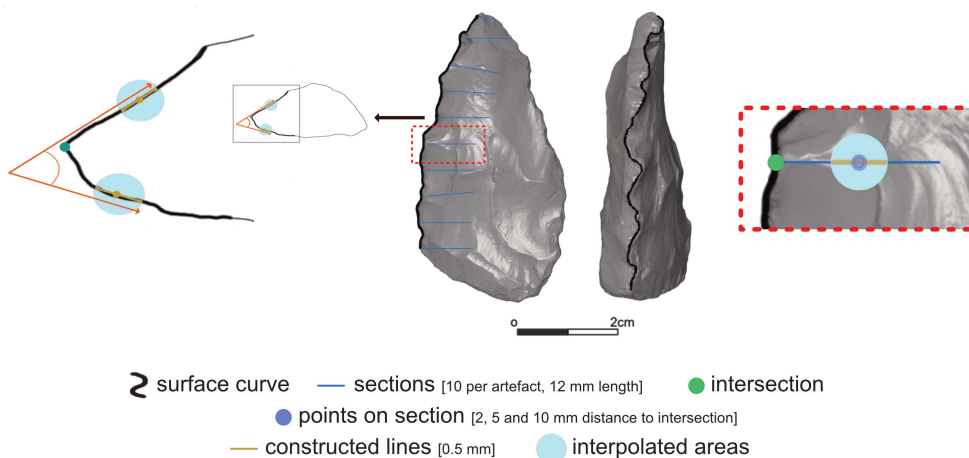

Schematic illustration of the "best-fit" measurement procedure.

6.1 Parameter

The parameters chosen for presenting the method are exemplary. In total, three different parameters are relevant for an edge angle calculation. One of these three is the number of sections defined along the tool's edge, chosen here to be ten sections per tool. It is possible to choose a fixed number of equidistant sections per tool, as done here, or to define steps, e.g., every few millimetres. Additionally, the length of the sections have to be

defined, which are 12 mm long in this case (see different options in the GitHub/Zenodo release: sections_SEC_NUM or sections_SEC_DIST).

Another parameter is the distance between the intersection and the points on the sections. In all three described measuring procedures, the intersection between the surface curve and the section determines one crucial point, but the other two points per measurement have to be defined numerically. Implemented as values in this study are 2 mm, 5 mm, and 10 mm.

The last parameter is the lengths of the two constructed lines used in the "2-lines" and "best-fit" measuring procedures, which is here always 0.5 mm.

Data analysis

7 Import and formatting of the CSV files

To format and sort the data acquired with the GOM Inspect pro script (see step #6), an R script was written.

This script reads all CSV files and imports them into single CSV files and single R files (one per sample).

Software

R Studio Desktop

NAME

The R Studio, Inc.

DEVELOPER

Available in open access on Zenodo: <https://doi.org/10.5281/zenodo.7961582>.

7.1 Statistical evaluation

A statistical analysis was then performed on the previously created three CSV files (see step #7). This step was done in the open-source software Python.

Software

python

NAME

3.6

OS

Guido van Rossum

DEVELOPER



Available in open access on Zenodo: <https://doi.org/10.5281/zenodo.7961582>.