

Aug 12, 2020

ISOLATING CELLS IN LOUPE BROWSER

Sumeet Gupta¹¹Whitehead Institute

1 Works for me

This document is published without a DOI.

Sumeet Gupta

ABSTRACT

ISOLATING CELLS IN LOUPE BROWSER

DOCUMENT CITATION

Sumeet Gupta 2020. ISOLATING CELLS IN LOUPE BROWSER. [protocols.io](https://protocols.io/view/isolating-cells-in-loupe-browser-bjmwkk7e)
<https://protocols.io/view/isolating-cells-in-loupe-browser-bjmwkk7e>

LICENSE

This is an open access document distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Aug 12, 2020

LAST MODIFIED

Aug 12, 2020

DOCUMENT INTEGER ID

40342

DISCLAIMER:


DISCLAIMER – FOR INFORMATIONAL PURPOSES ONLY; USE AT YOUR OWN RISK

The protocol content here is for informational purposes only and does not constitute legal, medical, clinical, or safety advice, or otherwise; content added to protocols.io is not peer reviewed and may not have undergone a formal approval of any kind. Information presented in this protocol should not substitute for independent professional judgment, advice, diagnosis, or treatment. Any action you take or refrain from taking using or relying upon the information presented here is strictly at your own risk. You agree that neither the Company nor any of the authors, contributors, administrators, or anyone else associated with protocols.io, can be held responsible for your use of the information contained in or linked to this protocol or any of our Sites/Apps and Services.


ISOLATING CELLS IN LOUPE BROWSER

Within Loupe Browser, you can follow these steps to isolate species specific cells. The example below isolates human cells.


1. Switch over to view the *Feature Plot*.

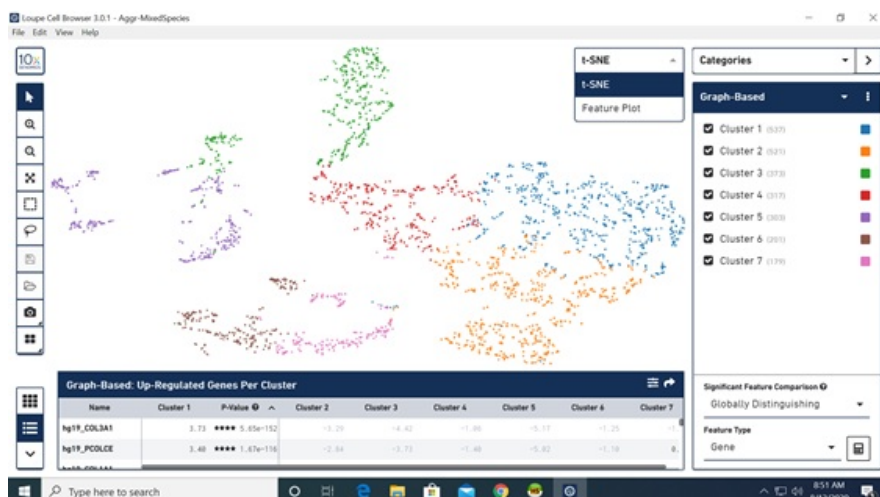


2. For the features, if you search for "sum", then you should see the human and mouse sums.

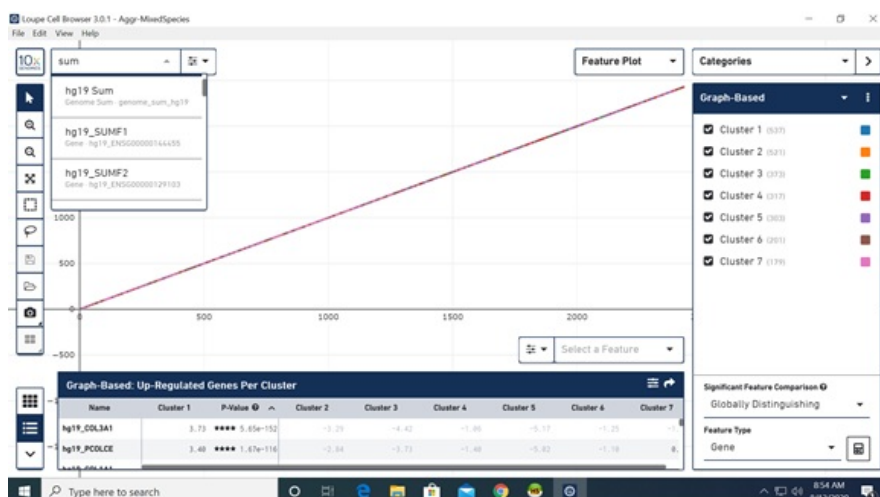


3. Select human sum for the y-axis, and mouse sum for the x-axis

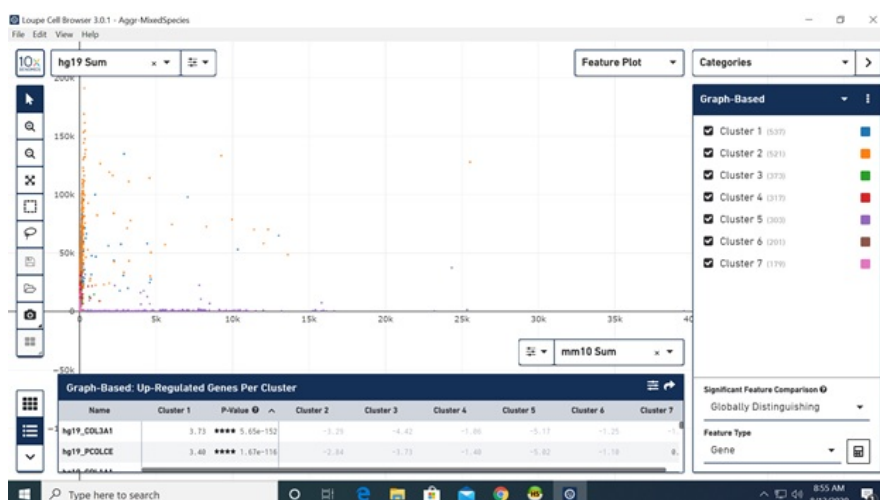

4. Now, make both axes logarithmic.

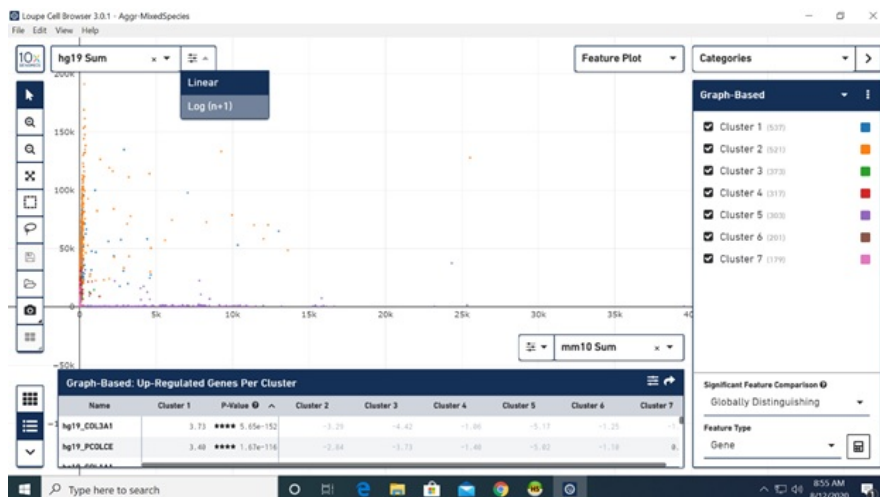

5. The human and mouse cells should be well separated, and you can select the desired cells using the lasso tool.

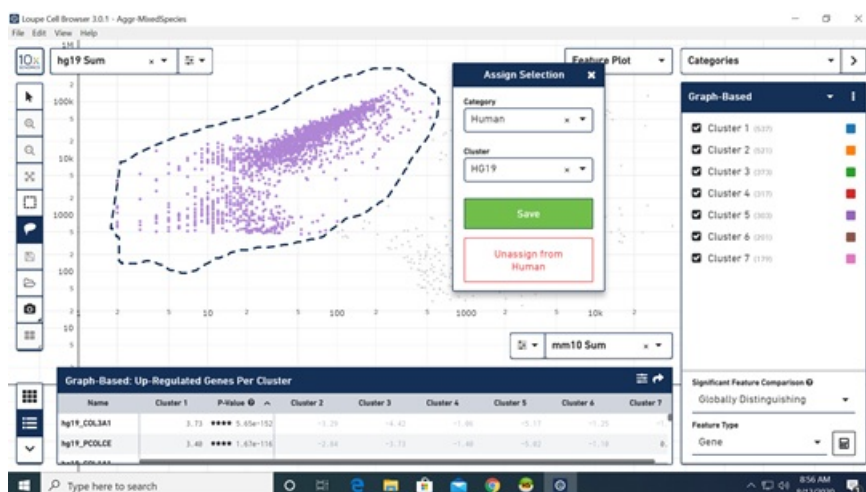

6. Once you have a custom cluster of human cells, you can switch back to t-SNE visualization if you wish.

