

JAN 22, 2024

OPEN ACCESS



Protocol Citation: Adrien Assie, Buck Samuel 2024. Mannitol

Agar. protocols.io

https://protocols.io/view/mannitolagar-zy8f7zw

License: This is an open access protocol distributed under the terms of the Creative Commons Attribution License, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working


Created: Apr 10, 2019

Last Modified: Jan 22, 2024

PROTOCOL integer ID: 22272

1 Start with 850mL distilled water

Mannitol Agar

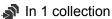

Adrien Assie¹, Buck Samuel¹

¹Baylor College of Medicine

Samuel Lab



Adrien Assie
Baylor College of Medicine

ABSTRACT

Mannitol agar recipe

2	3g Peptone
3	5g Yeast extract
4	25g Mannitol
5	Bring to 1000mL distilled water
6	(For LB Agar add 15 g agar)
7	Autoclave